

2023

*Product Catalog*

EST. 1923

**KNIT-RITE®**

**THUASNE**



**MEDROMA**

[www.medroma.co.za](http://www.medroma.co.za)

ORTHOTICS

PROSTHETICS

COMPRESSION GARMENTS

**SOLUTIONS THAT FIT**



Knit-Rite, LLC is a new subsidiary of Thuasne Group.


## LIFTING INDIVIDUALS TO NEW LEVELS SINCE 1923

Like many businesses, we have a mission—a mission that includes you. Whether you wear our products, prescribe them, or recommend them, Knit-Rite® is committed to products that help people live better lives. It might seem insignificant that a pair of compression socks could make life for someone a little easier and more mobile, but we know our products help lift individuals to a new level of living each day.

Knit-Rite® was founded in 1923 as an orthotic and prosthetic textile manufacturer with solutions that promote mobility, protection, and comfort. Therafirm®, part of the Knit-Rite® family of products, is a leading manufacturer of quality compression hosiery and socks. For more than 60 years, Therafirm® has made innovative compression garments, developing scientifically advanced products that are fashionable, comfortable and deliver the healthy benefits patients and consumers need. Our parent company, Thuasne®, shares in our mission of providing innovative products to help people master their health journey. Founded in 1847, Thuasne® designs, develops, and manufactures orthopedics, medical compression, home care and sport care products. Thuasne® is represented in more than 85 countries through its 15 brands in Europe, United States, and Asia.



### BRAND FAMILY

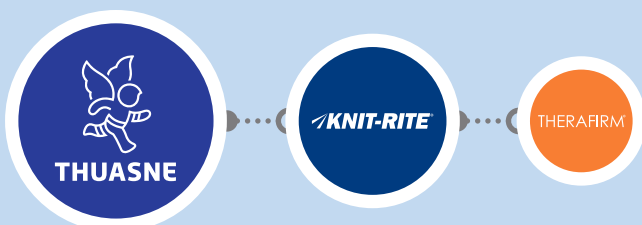

MEDROMA

Knit-Rite® is a Leading Designer, Manufacturer, and Marketer of Innovative Textiles for the Medical and Consumer Markets

Founded in 1923, Knit-Rite® remains committed to the continued research and development, manufacturing, and distributing of products that deliver healthy benefits for a better quality of life.

Therafirm®, a Division of Knit-Rite®, is a Leading Manufacturer of Quality Compression Hosiery and Socks

Since the 1950s, Therafirm® continues to innovate compression garments, developing scientifically advanced products that are fashionable and comfortable while delivering the healthy benefits patients and consumers need. Therafirm® offers a full line of true gradient compression hosiery in light, mild, moderate and firm support for a variety of markets.



**Interested in  
Therafirm compression?**

Scan to view and/or  
download our catalog.

### **ISO 13485 Certification**

is used by organizations involved in the design, production, installation and servicing of medical devices.



Manufactured under a BSI  
certified HACCP & GMP system.  
Certification #FM584357



Knit-Rite® owns manufacturing facilities in Kansas and North Carolina and is committed to producing products in the USA that allow a closer direct link to retailers and the ability to shorten our customers' supply chain management.

**Each plant has obtained ISO 13485 certification.**

**MEDROMA**

# TABLE OF CONTENTS

## PROSTHETIC TEXTILES | 6-41

### 6-25 PROSTHETIC SOCKS AND SHEATHS

- 7 Volume Management Solutions**
  - Shorty Soft-Sock®
  - Proximal Soft-Sock®
- 8-9 Soft-Sock® Prosthetic Sock**
  - Soft-Sock® Prosthetic Sock
  - Soft-Sock® Prosthetic Sock with Hole-In-Toe
  - Soft-Sock® Multi-Packs
- 10 X-Wool® Prosthetic Sock**
  - X-Wool® Prosthetic Sock
  - X-Wool® Prosthetic Sock with Hole-In-Toe
- 11 Stretch Spacer Fitting Sock**
  - Stretch Spacer Sock with Turn Down Cuff
  - Stretch Spacer Sock with Straight Top
  - Hole-In-Toe Stretch Spacer Sock with Turn Down Cuff
  - Hole-In-Toe Stretch Spacer Sock with Straight Top
- 12 Liner-Liner® Prosthetic Sock**
- 12 A/K Brim Sheath**
- 13 Prosthetic Sheaths**
  - Silver Sheath with Ionic+™
  - Super Stretch Sheath
  - Prosthetic Sheath
  - Prosthetic Sheath with Hole-In-Toe
- 14 Super Sock® Prosthetic Sock**
  - Super Sock® Prosthetic Sock
  - Super Sock® Prosthetic Sock with Hole-In-Toe
- 15 Sterling Stretch Prosthetic Sock**
  - Sterling Stretch Prosthetic Sock
  - Sterling Stretch Prosthetic Sock with Hole-In-Toe
- 16 Standard Symes Soft-Sock®**
  - Standard Sock Sizing Chart
- 17 Hip Disarticulation Prosthetic Sock**
- 18 Partial Foot Prosthetic Sock**
  - SmartKnit® Partial Foot Sock
  - SmartKnit® Partial Foot Sock with Heel
  - Partial Foot Sock with Heel
- 19 Post-Op Prosthetic Socks and Supplies**
  - Silver A/K and B/K Post-Op Socks
  - Original Unsterilized A/K and B/K Post-Op Socks
  - Symes and Chopart Unsterilized Post-Op Socks
- 20 Post-Op Supplies**
  - Distal End Pads
  - Lambs Wool

- 21 Materials and Care Instructions**
  - Knit-Rite Laundering Instructions
- 22 Prosthetic Sock Ordering Guidelines**
  - Knit-Rite Standard Prosthetic Sock Sizing
  - HPCPS Billing Codes
  - General Measurement Instructions
- 23-24 Special and Irregular Sock Ordering Guidelines**
  - How To Measure and Size
  - Special Order Fulfillment
  - Tailoring For Prosthetic and Orthotic Socks
- 25 Special Orders and Pricing**
  - Pricing On Irregular and Special Prosthetic Socks
  - Special Order Form

### 26-35 PROSTHETIC SHRINKERS

- 27 Compressogrip® B/K Elastic Shrinker**
- 28 Tapered Compressogrip® B/K Prosthetic Shrinker**
- 28 Bone Rings**
- 29 Compressogrip® A/K Prosthetic Shrinker**
- 30 Post-Op Dressing Retention Sock**
- 31-33 4-Way Stretch Shrinker**
  - B/K 4-Way Stretch Shrinker
  - A/K 4-Way Stretch Shrinker with Hip Attachment
  - A/K 4-Way Stretch Shrinker
- 34 Custom Prosthetic Shrinkers**
  - Made To Order A/K Shrinker Order Form
  - Made To Order Shrinker Order Form
- 35 Shrinker Private Logo Program**

### 36-41 PROSTHETIC SUSPENSION AND DONNING

- 37 Power Belt Auxiliary Suspension Belt®**
  - Original Power Belt® with Hand Loop
  - Power Belt® with Gunslinger Waist
  - Power Belt® with Adjustable Waist and Leg
- 38 Cool-Flex™ Auxiliary Suspension Belt**
- 38 B/K Suspension Belts**
  - Adjustable B/K Belt With Quick Release Fastener
  - B/K Suspension Belt
- 39 Power Pull® Donning Sock**
- 39 Shrinker Donning Tube**
- 40 Sleeve Saver™**
- 40 Foot Shell Sock**
- 41 Cosmetic Stockings**
  - A/K Cosmetic Stockings
  - B/K Cosmetic Stockings

## ORTHOTIC TEXTILES | 42–59

### 42–51 ORTHOTIC INTERFACES

- 43 Core-Spun by Therafirm® Patterned AFO Socks
- 44 SmartKnit® AFO and KAFO Interface Socks, Child
- 45 SmartKnit® AFO and KAFO Interface Socks, Adult
- 46 SmartKnit® Seamless SMO Interface Socks
- 47 SmartKnit® Big Toe Socks
- 48–49 Torso Interface®
  - Original V-Neck Tank
  - Original V-Neck Tank with Double Axilla Flaps
  - Original V-Neck Tank with Left Axilla Flap
  - Original Crew Neck Tank
  - Original Crew Neck Tank with Double Axilla Flaps
  - Original Crew Neck Tank with Sleeves
  - Original Strapless Shirt
  - Original Turtleneck with Sleeves
  - Lightweight V-Neck Tank
  - Lightweight V-Neck Tank with Left Axilla Flap
  - Lightweight V-Neck Tank with Double Axilla Flaps
- 50 Cranial Protection
  - Mullet Cranial Interface with Velcro Chin Closure
  - Balaclava Cranial Interface
  - Balaclava Scanning Interface
  - Cranial Scanning Sock
- 51 SmartKnit® Walker Boot Sock
- 51 Hip Orthosis Liner

### 52–53 FRACTURE INTERFACES

- 52 Cast-Rite™ Fracture Cast Sock
- 52 B/K Cast-Rite™ Fracture Brace Sock
- 53 Cast-Rite™ Humeral Sleeve
- 53 Cast-Rite™ Arm Sleeve

### 53–55 COMPRESSION INTERFACES

- 53 1-Inch Compression Bandage
- 54 SmartKnit® Knee Interface
  - SmartKnit® Knee Interface with Tapered Sleeves
  - SmartKnit® Knee Interface with Straight Sleeves
- 54 Compressogrip® Sleeves
- 55 Compressogrip® Tubular Elastic Bandage

### 56–59 DIABETIC AND SENSITIVITY SOCKS

- 57 SmartKnit® Seamless Diabetic Socks Information
- 58 SmartKnit® Seamless Diabetic Socks
  - SmartKnit® Mini Crew Socks
  - SmartKnit® Crew Socks
  - SmartKnit® Over-The-Calf Socks
- 59 SmartKnit® Seamless Wide Crew Diabetic Socks

## FABRICATION | 60–63

### 60–63 FABRICATION AND CASTING

- 61 Vacuum Forming Tubing and Hosiery
  - Nylon Vacuum Forming Tubing
  - Nyserts Vacuum Forming Hosiery
  - Sheer Vacuum Forming Hosiery
- 62–63 Stockinettes and Casting
  - Standard Rib-Knit Stockinette
  - Value Line Rib-Knit Stockinette
  - Cast Socks
  - A/K Casting Pants
  - AccuCast® Stretch Stockinette
  - AccuCast® Casting Socks

## TERMS AND CONDITIONS | 64

## MISSION STATEMENT | 65

# PROSTHETIC SOCKS AND SHEATHS

## PROSTHETIC SOCKS AND SHEATHS | 6-25

- 7** Volume Management Solutions
- 8-9** Soft-Sock® Prosthetic Sock
- 10** X-Wool® Prosthetic Sock
- 11** Stretch Spacer Fitting Sock
- 12** Liner-Liner® Prosthetic Sock
- 12** A/K Brim Sheath
- 13** Prosthetic Sheaths
- 14** Super Sock® Prosthetic Sock
- 15** Sterling Stretch Prosthetic Sock
- 16** Standard Symes Soft-Sock®
- 17** Hip Disarticulation Prosthetic Sock
- 18** Partial Foot Prosthetic Sock
- 19** Post-Op Prosthetic Socks and Supplies
- 20** Post-Op Supplies
- 21** Materials and Care Instructions
- 22** Prosthetic Sock Ordering Guidelines
- 23-24** Special and Irregular Sock Ordering Guidelines
- 25** Special Orders and Pricing



# VOLUME MANAGEMENT SOLUTIONS



SHORTY  
SOFT-SOCK®

- ✓ Designed for use with elevated vacuum and total surface bearing sockets
- ✓ Short lengths to reduce volume mainly in the distal one-third of the socket and remain below the suction sleeve
- ✓ Features distal hole
- ✓ Color-coded to indicate length for ease of patient instruction
- ✓ Knit with easy care polyester
- ✓ Sold each

## SHORTY SOFT-SOCK®

| ITEM #     | TOP STRIPE | TOP    | TOE    | LENGTH |
|------------|------------|--------|--------|--------|
| 1SPS-S1-81 | Red        | 4-3/4" | 2-3/4" | 7"     |
| 1SPS-S1-91 | Red        | 5-3/4" | 3-1/8" | 5-5/8" |
| 1SPS-S1-82 | Black      | 3-1/8" | 2-5/8" | 3-5/8" |
| 1SPS-S1-92 | Black      | 4-1/2" | 4"     | 3-1/4" |




PROXIMAL  
SOFT-SOCK®

- ✓ Designed to manage proximal volume changes
- ✓ Tapered and open-ended
- ✓ Knit with polyester and LYCRA®
- ✓ Sold each

## PROXIMAL SOFT-SOCK®

| LIGHTWEIGHT | 3-PLY      | 5-PLY      | SIZE           | FITS SOCK SIZES | LENGTH        |
|-------------|------------|------------|----------------|-----------------|---------------|
| 1PS1-RG-SH  | 1PS3-RG-SH | 1PS5-RG-SH | Regular Short  | 1, 2            | 10.50" (27cm) |
| 1PS1-RG-MD  | 1PS3-RG-MD | 1PS5-RG-MD | Regular Medium | 1, 2            | 14" (35cm)    |
| 1PS1-WD-MD  | 1PS3-WD-MD | 1PS5-WD-MD | Wide Medium    | 3+              | 14" (35cm)    |
| 1PS1-XW-MD  | 1PS3-XW-MD | 1PS5-XW-MD | X-Wide Medium  | 4+              | 14" (35cm)    |

## THICKNESS COLOR CODE (TOP SEAM)

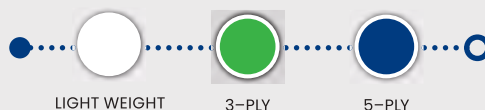

# SOFT-SOCK® PROSTHETIC SOCK



- ✓ Non-roll top looks great and eliminates bunching
- ✓ Knit-in 'Eye-Hole' reinforced distal hole option has no added thickness and minimal stretch
- ✓ Seamless three-dimensional, rounded toe
- ✓ Patent No. 7,363,778
- ✓ Soft-Sock® feels soft and cuddly worn next to the skin
- ✓ Wicks moisture away
- ✓ Helps inhibit odor in the sock
- ✓ Stretches and recovers to conform to limb
- ✓ Fine knit gives comfort and longevity
- ✓ Easy-to-fit sizes
- ✓ Easy care and shrinkage resistant

## THICKNESS COLOR CODE (TOP SEAM)



### SOFT-SOCK WITH WHITE COOLMAX®/LYCRA®

| LIGHTWEIGHT | 3-PLY      | 5-PLY      | 6-PLY      | SIZE            | FITS SOCK SIZE | LENGTH            |
|-------------|------------|------------|------------|-----------------|----------------|-------------------|
| 1SP1-CH-XS  | —          | —          | —          | Child X-Short   | A              | 6" (15cm)         |
| 1SP1-CH-SH  | 1SP3-CH-SH | 1SP5-CH-SH | 1SP6-CH-SH | Child Short     | A              | 8" (20cm)         |
| 1SP1-CH-12  | 1SP3-CH-12 | 1SP5-CH-12 | —          | Child Medium    | A              | 12" (30cm)        |
| 1SP1-NA-SH  | 1SP3-NA-SH | 1SP5-NA-SH | 1SP6-NA-SH | Narrow Short    | B, 0           | 10"-14" (25-35cm) |
| 1SP1-NA-MD  | 1SP3-NA-MD | 1SP5-NA-MD | 1SP6-NA-MD | Narrow Medium   | B, 0           | 16"-20" (40-50cm) |
| 1SP1-NA-LG  | 1SP3-NA-LG | 1SP5-NA-LG | 1SP6-NA-LG | Narrow Long     | B, 0           | 22"-28" (55-70cm) |
| 1SP1-RG-XS  | 1SP3-RG-XS | 1SP5-RG-XS | 1SP6-RG-XS | Regular X-Short | 1, 2           | 8"-12" (20-30cm)  |
| 1SP1-RG-SH  | 1SP3-RG-SH | 1SP5-RG-SH | 1SP6-RG-SH | Regular Short   | 1, 2           | 10"-14" (25-35cm) |
| 1SP1-RG-MD  | 1SP3-RG-MD | 1SP5-RG-MD | 1SP6-RG-MD | Regular Medium  | 1, 2           | 16"-20" (40-50cm) |
| 1SP1-RG-LG  | 1SP3-RG-LG | 1SP5-RG-LG | 1SP6-RG-LG | Regular Long    | 1, 2           | 22"-28" (55-70cm) |
| 1SP1-WD-SH  | 1SP3-WD-SH | 1SP5-WD-SH | 1SP6-WD-SH | Wide Short      | 3 +            | 10"-14" (25-35cm) |
| 1SP1-WD-MD  | 1SP3-WD-MD | 1SP5-WD-MD | 1SP6-WD-MD | Wide Medium     | 3 +            | 16"-20" (40-50cm) |
| 1SP1-XW-SH  | 1SP3-XW-SH | 1SP5-XW-SH | —          | X-Wide Short    | 4, 5+          | 10"-14" (25-35cm) |
| 1SP1-XW-MD  | 1SP3-XW-MD | 1SP5-XW-MD | —          | X-Wide Medium   | 4, 5+          | 16"-20" (40-50cm) |

### SOFT-SOCK WITH BLACK POLYESTER/LYCRA®

| LIGHTWEIGHT | 3-PLY      | 5-PLY      | SIZE            | FITS SOCK SIZE | LENGTH            |
|-------------|------------|------------|-----------------|----------------|-------------------|
| 1SB1-NA-SH  | 1SB3-NA-SH | 1SB5-NA-SH | Narrow Short    | B, 0           | 10"-14" (25-35cm) |
| 1SB1-NA-MD  | 1SB3-NA-MD | 1SB5-NA-MD | Narrow Medium   | B, 0           | 16"-20" (40-50cm) |
| 1SB1-NA-LG  | 1SB3-NA-LG | 1SB5-NA-LG | Narrow Long     | B, 0           | 22"-28" (55-70cm) |
| 1SB1-RG-XS  | 1SB3-RG-XS | 1SB5-RG-XS | Regular X-Short | 1, 2           | 8"-12" (20-30cm)  |
| 1SB1-RG-SH  | 1SB3-RG-SH | 1SB5-RG-SH | Regular Short   | 1, 2           | 10"-14" (25-35cm) |
| 1SB1-RG-MD  | 1SB3-RG-MD | 1SB5-RG-MD | Regular Medium  | 1, 2           | 16"-20" (40-50cm) |
| 1SB1-RG-LG  | 1SB3-RG-LG | 1SB5-RG-LG | Regular Long    | 1, 2           | 22"-28" (55-70cm) |
| 1SB1-WD-SH  | 1SB3-WD-SH | 1SB5-WD-SH | Wide Short      | 3 +            | 10"-14" (25-35cm) |
| 1SB1-WD-MD  | 1SB3-WD-MD | 1SB5-WD-MD | Wide Medium     | 3 +            | 16"-20" (40-50cm) |
| 1SB1-XW-SH  | 1SB3-XW-SH | 1SB5-XW-SH | X-Wide Short    | 4, 5+          | 10"-14" (25-35cm) |
| 1SB1-XW-MD  | 1SB3-XW-MD | 1SB5-XW-MD | X-Wide Medium   | 4, 5+          | 16"-20" (40-50cm) |

### SOFT-SOCK WITH POLYESTER/IONIC+™ AND SILVER/LYCRA®

| LIGHTWEIGHT | 3-PLY      | 5-PLY      | SIZE            | FITS SOCK SIZE | LENGTH            |
|-------------|------------|------------|-----------------|----------------|-------------------|
| 1SX1-NA-SH  | 1SX3-NA-SH | 1SX5-NA-SH | Narrow Short    | B, 0           | 10"-14" (25-35cm) |
| 1SX1-RG-XS  | 1SX3-RG-XS | 1SX5-RG-XS | Regular X-Short | 1, 2           | 8"-12" (20-30cm)  |
| 1SX1-RG-SH  | 1SX3-RG-SH | 1SX5-RG-SH | Regular Short   | 1, 2           | 10"-14" (25-35cm) |
| 1SX1-RG-MD  | 1SX3-RG-MD | 1SX5-RG-MD | Regular Medium  | 1, 2           | 16"-20" (40-50cm) |
| 1SX1-WD-MD  | 1SX3-WD-MD | 1SX5-WD-MD | Wide Medium     | 3 +            | 16"-20" (40-50cm) |



# SOFT-SOCK® PROSTHETIC SOCK

## HOLE-IN-TOE SOFT-SOCK® WITH WHITE COOLMAX®/LYCRA®

| LIGHTWEIGHT | 3-PLY      | 5-PLY      | 6-PLY      | SIZE            | FITS SOCK SIZE | LENGTH            |
|-------------|------------|------------|------------|-----------------|----------------|-------------------|
| 1SPS-CH-XS  | —          | —          | —          | Child X-Short   | A              | 6" (15cm)         |
| 1SPS-CH-SH  | 1SP3-SC-SH | 1SP5-SC-SH | 1SP6-SC-SH | Child Short     | A              | 8" (20cm)         |
| 1SPS-NA-SH  | 1SP3-SN-SH | 1SP5-SN-SH | 1SP6-SN-SH | Narrow Short    | B, 0           | 10"-14" (25-35cm) |
| 1SPS-NA-MD  | 1SP3-SN-MD | 1SP5-SN-MD | 1SP6-SN-MD | Narrow Medium   | B, 0           | 16"-20" (40-50cm) |
| 1SPS-RG-XS  | 1SP3-SR-XS | 1SP5-SR-XS | 1SP6-SR-XS | Regular X-Short | 1, 2           | 8"-12" (20-30cm)  |
| 1SPS-RG-SH  | 1SP3-SR-SH | 1SP5-SR-SH | 1SP6-SR-SH | Regular Short   | 1, 2           | 10"-14" (25-35cm) |
| 1SPS-RG-MD  | 1SP3-SR-MD | 1SP5-SR-MD | 1SP6-SR-MD | Regular Medium  | 1, 2           | 16"-20" (40-50cm) |
| 1SPS-RG-LG  | 1SP3-SR-LG | 1SP5-SR-LG | 1SP6-SR-LG | Regular Long    | 1, 2           | 22"-28" (55-70cm) |
| 1SPS-WD-SH  | 1SP3-SW-SH | 1SP5-SW-SH | 1SP6-SW-SH | Wide Short      | 3 +            | 10"-14" (25-35cm) |
| 1SPS-WD-MD  | 1SP3-SW-MD | 1SP5-SW-MD | 1SP6-SW-MD | Wide Medium     | 3 +            | 16"-20" (40-50cm) |

### CARE TIPS

Wash Soft-Socks  
fleece side in.  
No bleach. No  
liquid softeners.  
Tumble dry, cool  
temperatures.

## HOLE-IN-TOE SOFT-SOCK® WITH BLACK POLYESTER/LYCRA®

| LIGHTWEIGHT | 3-PLY      | 5-PLY      | SIZE            | FITS SOCK SIZE | LENGTH            |
|-------------|------------|------------|-----------------|----------------|-------------------|
| 1SBS-NA-SH  | 1SB3-SN-SH | 1SB5-SN-SH | Narrow Short    | B, 0           | 10"-14" (25-35cm) |
| 1SBS-NA-MD  | 1SB3-SN-MD | 1SB5-SN-MD | Narrow Medium   | B, 0           | 16"-20" (40-50cm) |
| 1SBS-RG-XS  | 1SB3-SR-XS | 1SB5-SR-XS | Regular X-Short | 1, 2           | 8"-12" (20-30cm)  |
| 1SBS-RG-SH  | 1SB3-SR-SH | 1SB5-SR-SH | Regular Short   | 1, 2           | 10"-14" (25-35cm) |
| 1SBS-RG-MD  | 1SB3-SR-MD | 1SB5-SR-MD | Regular Medium  | 1, 2           | 16"-20" (40-50cm) |
| 1SBS-RG-LG  | 1SB3-SR-LG | 1SB5-SR-LG | Regular Long    | 1, 2           | 22"-28" (55-70cm) |
| 1SBS-WD-SH  | 1SB3-SW-SH | 1SB5-SW-SH | Wide Short      | 3 +            | 10"-14" (25-35cm) |
| 1SBS-WD-MD  | 1SB3-SW-MD | 1SB5-SW-MD | Wide Medium     | 3 + 1          | 16"-20" (40-50cm) |

See page 22 for sizing  
information, and  
pages 23-25 for ordering  
information about special  
and irregular socks.

## MULTI-PACK SOFT-SOCK® PROSTHETIC SOCK

CONTAINS: (6) LIGHTWEIGHT SOCKS, (3)  
3-PLY SOCKS, (3) 5-PLY SOCKS

- ✓ Complete starter pack for new patients
- ✓ Multi-Packs available in original white COOLMAX® and black moisture-wicking Polyester (with or without 'Eye Hole')
- ✓ One part number to easily order a variety of plies
- ✓ Fewer SKUs to simplify inventory

## SOFT-SOCK® MULTI-PACK WITH WHITE COOLMAX®/LYCRA®

| WITH HOLE   | WITHOUT HOLE | SIZE            | FITS SOCK SIZE | LENGTH            |
|-------------|--------------|-----------------|----------------|-------------------|
| KRSP-RXS-WH | KRSP-RXS-NH  | Regular X-Short | #1, #2         | 8"-12" (20-30cm)  |
| KRSP-RS-WH  | KRSP-RS-NH   | Regular Short   | #1, #2         | 10"-14" (25-35cm) |
| KRSP-RM-WH  | KRSP-RM-NH   | Regular Medium  | #1, #2         | 16"-20" (40-50cm) |
| KRSP-WS-WH  | KRSP-WS-NH   | Wide Short      | #3+            | 10"-14" (25-35cm) |
| KRSP-WM-WH  | KRSP-WM-NH   | Wide Medium     | #3+            | 16"-20" (40-50cm) |

## SOFT-SOCK® MULTI-PACK WITH BLACK POLYESTER/LYCRA®

| WITH HOLE   | WITHOUT HOLE | SIZE            | FITS SOCK SIZE | LENGTH            |
|-------------|--------------|-----------------|----------------|-------------------|
| KRBP-RXS-WH | KRBP-RXS-NH  | Regular X-Short | #1, #2         | 8"-12" (20-30cm)  |
| KRBP-RS-WH  | KRBP-RS-NH   | Regular Short   | #1, #2         | 10"-14" (25-35cm) |
| KRBP-RM-WH  | KRBP-RM-NH   | Regular Medium  | #1, #2         | 16"-20" (40-50cm) |



## X-WOOL® PROSTHETIC SOCK




X-WOOL® PROSTHETIC SOCK WITH HOLE-IN-TOE



X-Wool® Prosthetic Sock

### THICKNESS COLOR CODE (TOP SEAM)

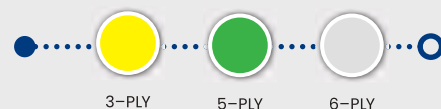

- ✓ Thick, absorbant virgin wool provides cushioning
- ✓ Contours to the limb without stretching out
- ✓ Fleeced lining adds softness
- ✓ Ionic+™ Mineral Antimicrobial assists in preventing the growth of odor-causing bacteria and fungi
- ✓ Ionic+™ Mineral Antimicrobial is thermodynamic, keeping you cooler in the summer and warmer in the winter
- ✓ Machine washable and resists shrinkage

### CARE TIPS

Machine washable.  
Resists shrinkage.

### X-WOOL® PROSTHETIC SOCK WITH HOLE-IN-TOE

| 3-PLY      | 5-PLY      | SIZE           | FITS SOCK SIZE | LENGTH            |
|------------|------------|----------------|----------------|-------------------|
| 1WX3-SN-XS | 1WX5-SN-XS | Narrow X-Short | 00, 0          | 8"-10" (20-25cm)  |
| 1WX3-SN-SH | 1WX5-SN-SH | Narrow Short   | 00, 0          | 10"-12" (25-30cm) |
| 1WX3-SN-RG | 1WX5-SN-RG | Narrow Regular | 00, 0          | 14"-16" (35-40cm) |
| 1WX3-SN-LG | 1WX5-SN-LG | Narrow Long    | 00, 0          | 18"-20" (45-50cm) |
| 1WX3-SM-XS | 1WX5-SM-XS | Medium X-Short | 1, 2           | 8"-10" (20-25cm)  |
| 1WX3-SM-SH | 1WX5-SM-SH | Medium Short   | 1, 2           | 10"-12" (25-30cm) |
| 1WX3-SM-RG | 1WX5-SM-RG | Medium Regular | 1, 2           | 14"-16" (35-40cm) |
| 1WX3-SM-LG | 1WX5-SM-LG | Medium Long    | 1, 2           | 18"-20" (45-50cm) |
| 1WX3-SW-SH | 1WX5-SW-SH | Wide Short     | 3, 4           | 10"-12" (25-30cm) |
| 1WX3-SW-MD | 1WX5-SW-MD | Wide Medium    | 3, 4           | 14"-16" (35-40cm) |

### X-WOOL® PROSTHETIC SOCK

| 3-PLY      | 5-PLY      | 6-PLY      | SIZE           | FITS SOCK SIZE | LENGTH            |
|------------|------------|------------|----------------|----------------|-------------------|
| 1WX3-NA-XS | 1WX5-NA-XS | —          | Narrow X-Short | 00, 0          | 8"-10" (20-25cm)  |
| 1WX3-NA-SH | 1WX5-NA-SH | 1WX6-NA-SH | Narrow Short   | 00, 0          | 10"-12" (25-30cm) |
| 1WX3-NA-RG | 1WX5-NA-RG | 1WX6-NA-RG | Narrow Regular | 00, 0          | 14"-16" (35-40cm) |
| 1WX3-NA-LG | 1WX5-NA-LG | 1WX6-NA-LG | Narrow Long    | 00, 0          | 18"-20" (45-50cm) |
| 1WX3-MD-XS | 1WX5-MD-XS | —          | Medium X-Short | 1, 2           | 8"-10" (20-25cm)  |
| 1WX3-MD-SH | 1WX5-MD-SH | 1WX6-MD-SH | Medium Short   | 1, 2           | 10"-12" (25-30cm) |
| 1WX3-MD-RG | 1WX5-MD-RG | 1WX6-MD-RG | Medium Regular | 1, 2           | 14"-16" (35-40cm) |
| 1WX3-MD-LG | 1WX5-MD-LG | 1WX6-MD-LG | Medium Long    | 1, 2           | 18"-20" (45-50cm) |
| 1WX3-MD-XL | 1WX5-MD-XL | —          | Medium X-Long  | 1, 2           | 22"-24" (55-60cm) |
| 1WX3-MD-XX | 1WX5-MD-XX | —          | Medium XX-Long | 1, 2           | 26"-28" (65-70cm) |
| 1WX3-WD-SH | 1WX5-WD-SH | 1WX6-WD-SH | Wide Short     | 3, 4           | 10"-12" (25-30cm) |
| 1WX3-WD-MD | 1WX5-WD-MD | 1WX6-WD-MD | Wide Medium    | 3, 4           | 14"-16" (35-40cm) |
| 1WX3-WD-LG | 1WX5-WD-LG | 1WX6-WD-LG | Wide Long      | 3, 4           | 18"-20" (45-50cm) |
| 1WX3-WD-XL | 1WX5-WD-XL | —          | Wide X-Long    | 3, 4           | 22"-24" (55-60cm) |
| 1WX3-WD-XX | 1WX5-WD-XX | —          | Wide XX-Long   | 3, 4           | 26"-28" (65-70cm) |
| 1WX3-EW-SH | 1WX5-EW-SH | —          | X-Wide Short   | 5, 6           | 12"-14" (30-35cm) |
| 1WX3-EW-MD | 1WX5-EW-MD | —          | X-Wide Medium  | 5, 6           | 16"-18" (40-45cm) |

**GREY IONIC+™ STYLES  
CONTAIN IONIC+™  
MINERAL ANTIMICROBIAL  
FABRIC KNOWN FOR  
ITS MOISTURE WICKING  
PROPERTIES**

# STRETCH SPACER FITTING SOCKS

- ✓ Available in 1-ply thickness
- ✓ Rounded toe is reinforced for good fit and longer wear
- ✓ Toe seam is one inch above the distal end, like a sheath
- ✓ Blended cotton polyester fabric is soft, comfortable, and absorbent
- ✓ Provides improved fitting qualities, range of motion, and reduced migration
- ✓ Grips the residuum for an intimate fit
- ✓ Smooth surface makes socks ideal for casting and computer modeling
- ✓ Pre-washed so size remains after washing
- ✓ Sold in pack of 6 socks



Stretch Spacer Sock with Nylon Toe (left) and Stretch Spacer Sock Hole-in-Toe with Reinforced Toe (right)

## COTTON/LYCRA® STRETCH SPACER SOCKS WITH TURN DOWN CUFF

| ITEM #     | SIZE         | FITS SIZES | UNSTRETCHED LENGTH | TOE MATERIAL      |
|------------|--------------|------------|--------------------|-------------------|
| 6S-SS35-SH | Narrow Short | #A, #B     | 12" (30cm)         | Cotton/Lycra® Toe |
| 6S-SS40-SH | Medium Short | #0, #1     | 12" (30cm)         | Cotton/Lycra® Toe |
| 6S-SS45-WS | Wide Short   | #2, #3, #4 | 12" (30cm)         | Cotton/Lycra® Toe |

## COTTON/LYCRA® STRETCH SPACER SOCKS WITH STRAIGHT TOP

| ITEM #     | SIZE       | FITS SIZES | UNSTRETCHED LENGTH | TOE MATERIAL |
|------------|------------|------------|--------------------|--------------|
| 6S-8S35-MD | Narrow     | #A, #B,    | 12" (30cm)         | Nylon Toe    |
| 6S-8S40-MD | Medium     | #0, #1     | 15" (38cm)         | Nylon Toe    |
| 6S-8S45-SH | Wide Short | #2, #3, #4 | 12" (30cm)         | Nylon Toe    |
| 6S-8S45-MD | Wide       | #2, #3, #4 | 18" (45cm)         | Nylon Toe    |
| 6S-8C45-MD | Wide       | #2, #3, #4 | 18" (45cm)         | Cotton Toe   |
| 6S-8C50-18 | X-Wide     | #4, #5, #6 | 18" (45cm)         | Cotton Toe   |

## COTTON/LYCRA® STRETCH SPACER SOCKS HOLE-IN-TOE WITH TURN DOWN CUFF


| ITEM #     | SIZE         | FITS SIZES | UNSTRETCHED LENGTH |
|------------|--------------|------------|--------------------|
| 6S-8S35-SH | Narrow Short | #A, #B     | 12" (30cm)         |
| 6S-8S40-SH | Medium Short | #0, #1     | 12" (30cm)         |
| 6S-8S45-WS | Wide Short   | #2, #3, #4 | 12" (30cm)         |

## COTTON/LYCRA® STRETCH SPACER SOCKS HOLE-IN-TOE WITH STRAIGHT TOP


| ITEM #     | SIZE       | FITS SIZES | UNSTRETCHED LENGTH | TOE MATERIAL |
|------------|------------|------------|--------------------|--------------|
| 6S-SS35-MD | Narrow     | #A, #B,    | 12" (30cm)         | Nylon Toe    |
| 6S-SS40-MD | Medium     | #0, #1     | 15" (38cm)         | Nylon Toe    |
| 6S-SS45-SH | Wide Short | #2, #3, #4 | 12" (30cm)         | Nylon Toe    |
| 6S-SS45-MD | Wide       | #2, #3, #4 | 18" (45cm)         | Nylon Toe    |
| 6S-SC45-MD | Wide       | #2, #3, #4 | 18" (45cm)         | Cotton Toe   |
| 6S-SC50-18 | X-Wide     | #4, #5, #6 | 18" (45cm)         | Cotton Toe   |



ORIGINAL

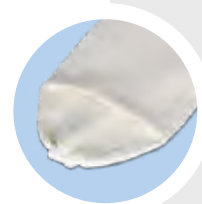

HOLE-IN-TOE

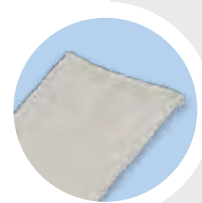

STRAIGHT TOP



TURN DOWN CUFF




## LINER-LINER® PROSTHETIC SOCK



**A PROSTHETIC SOCK DESIGNED TO BE WORN UNDER A SUSPENSION LINER NEXT TO THE SKIN.**

- ✓ Relieves skin shear irritations and improves comfort with liners
- ✓ Tested and proven beneficial for amputees who wear liners as the primary interface
- ✓ Ultra-thin material minimizes compromise to liner's suction grip on skin
- ✓ Moisture-wicking fibers for a cool and comfortable garment
- ✓ Ionic+™ Mineral Antimicrobial assists in preventing odor in the garment
- ✓ Washable interface keeps liner cleaner and helps control skin irritations
- ✓ High stretch for excellent fitting qualities

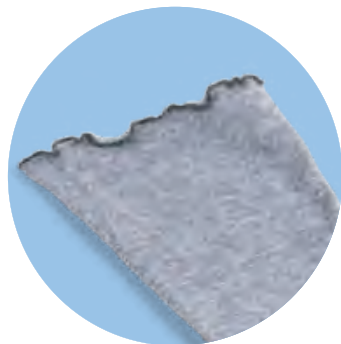

**NOTE TO PROSTHETISTS:** Due to the sheerness of sock material, liner suspension should not be affected; however, suspension should be examined and verified for each wearer. Advise and select candidates accordingly.

### CARE TIPS

Wear clean socks each day. Pre-spot in cool water. Machine wash and dry, warm, all-purpose detergent. No bleach. Rinse thoroughly.

### LINER-LINER® PROSTHETIC SOCK

| ITEM #     | SIZE             | TOP WIDTH | TOE WIDTH | LENGTH             |
|------------|------------------|-----------|-----------|--------------------|
| 1S-XLNA-XS | Narrow X-Short   | 3.5"      | 2"        | 8"-12" (20-31cm)   |
| 1S-XLNA-SH | Narrow Short     | 3.5"      | 2"        | 10"-14" (25-35cm)  |
| 1S-XLRX-XS | Regular XX-Short | 4.5"      | 2"        | 2"6"-10" (15-25cm) |
| 1S-XLRG-XS | Regular X-Short  | 4.75"     | 2"        | 8"-12" (20-31cm)   |
| 1S-XLRG-SH | Regular Short    | 5"        | 2"        | 10"-14" (25-35cm)  |
| 1S-XLRG-MD | Regular Medium   | 5"        | 2"        | 16"-20" (40-50cm)  |
| 1S-XLWD-SH | Wide Short       | 7"        | 3.5"      | 10"-14" (25-35cm)  |
| 1S-XLWD-MD | Wide Medium      | 7"        | 3.5"      | 16"-20" (40-50cm)  |

## A/K BRIM SHEATH

- ✓ For above knee patients, for all socket styles
- ✓ Designed to separate and protect the soft tissue in the top brim of the socket
- ✓ Used to bridge the gap between the top edge of any roll on liner and the groin
- ✓ Works on pin suspension and seal in liners.
- ✓ Uses stretchy, slick nylon material, constructed in two angled sections to allow the fabric to fit the high lateral and lower medial shape appropriately
- ✓ Includes gripper elastic to secure to outer socket



A/K Brim Sheath



### A/K BRIM SHEATH

| ITEM #     | SIZE    | SOCKET CIRCUMFERENCE        |
|------------|---------|-----------------------------|
| 2H-BSAK-SM | Small   | 14" - 18" (35.5 - 45.7 cm)  |
| 2H-BSAK-MD | Medium  | 18" - 22" (45.7 - 55.8 cm)  |
| 2H-BSAK-LG | Large   | 22" - 30" (55.8 - 76.2 cm)  |
| 2H-BSAK-XL | X-Large | 30" - 40" (76.2 - 101.6 cm) |

# PROSTHETIC SHEATHS

- ✓ Sheaths control friction, adjust fit, and aid in skin care
- ✓ Silver Sheath with Ionic+™ Mineral Antimicrobial inhibits odor in the garment, wicks moisture away from skin, and aids in thermal consistency
- ✓ Super Stretch Sheaths add LYCRA® to accommodate an even greater range and more comfortable fit, enhance flexion, movement and improve suspension



## CARE TIPS

Sheaths may be hand or machine washed. Use soap or mild detergent in warm water. No bleach. Rinse thoroughly. Air dry overnight or tumble dry low.

## SILVER SHEATH WITH IONIC+™

| SILVER     | SOCK SIZE                 | LENGTH  |
|------------|---------------------------|---------|
| 2H-PXNA-XS | Narrow X-Short #B, #O, #1 | 8"-10"  |
| 2H-PXNA-SH | Narrow Short #B, #O, #1   | 10"-14" |
| 2H-PXNA-MD | Narrow Medium #B, #O, #1  | 16"-20" |
| 2H-PXRG-XS | Regular X-Short #1, #2    | 8"-10"  |
| 2H-PXRG-SH | Regular Short #1, #2      | 10"-14" |
| 2H-PXRG-MD | Regular Medium #1, #2     | 16"-20" |
| 2H-PXWD-SH | Wide Short #3 and above   | 10"-14" |

- ✓ 91% Nylon
- ✓ 9% pure Silver Ionic+™
- ✓ Sold per box (3 pair)

## SUPER STRETCH SHEATH

| DUSK       | WHITE      | SOCK SIZE             | LENGTH   |
|------------|------------|-----------------------|----------|
| 2H-PSBR-SH | 2H-PSWR-SH | Regular Short #1, #2  | 10"-14"+ |
| 2H-PSBR-MD | 2H-PSWR-MD | Regular Medium #1, #2 | 16"-20"+ |

- ✓ 98% Nylon
- ✓ 2% LYCRA®
- ✓ Sold per box (3 pair)

## PROSTHETIC SHEATH

| DUSK       | WHITE      | SOCK SIZE                          | LENGTH  |
|------------|------------|------------------------------------|---------|
| 2H-PBCH-SH | 2H-PWCH-SH | Child Short #A, LYCRA® Spandex top | 6"-8"   |
| 2H-PBCH-12 | 2H-PWCH-12 | Child 12" #A, non-elastic top      | 12"     |
| 2H-PBNA-SH | 2H-PWNA-SH | Narrow Short #B, #O, #1            | 10"-14" |
| 2H-PBNA-MD | 2H-PWNA-MD | Narrow Medium #B, #O, #1           | 16"-20" |
| 2H-PBNA-LG | 2H-PWNA-LG | Narrow Long #B, #O, #1             | 22"-28" |
| 2H-PBRG-XS | 2H-PWRG-XS | Regular X-Short #1, #2             | 8"-10"  |
| 2H-PBRG-SH | 2H-PWRG-SH | Regular Short #1, #2               | 10"-14" |
| 2H-PBRG-MD | 2H-PWRG-MD | Regular Medium #1, #2              | 16"-20" |
| 2H-PBRG-LG | 2H-PWRG-LG | Regular Long #1, #2                | 22"-28" |
| 2H-PBWD-SH | 2H-PWWD-SH | Wide Short #3 and above            | 10"-14" |
| 2H-PBWD-MD | 2H-PWWD-MD | Wide Medium #3 and above           | 16"-20" |
| 2H-PBWD-LG | 2H-PWWD-LG | Wide Long #3 and above             | 22"-28" |

- ✓ 100% Stretch Nylon
- ✓ Sized 0 socks
- ✓ Sold per box (3 pair)

## PROSTHETIC SHEATH WITH HOLE-IN-TOE

| DUSK       | WHITE      | SOCK SIZE               | LENGTH  |
|------------|------------|-------------------------|---------|
| 2H-PBRH-SH | 2H-PWRH-SH | Regular Short #1, #2    | 10"-14" |
| 2H-PBWH-SH | 2H-PWWH-SH | Wide Short #3 and above | 10"-14" |

- ✓ 100% Stretch Nylon
- ✓ Sold per pair



# SUPER SOCK® PROSTHETIC SOCK

See pages 23-25 for information needed for irregular and special socks.



Super Sock  
Prosthetic Sock

- ✓ 100% Virgin Wool – uniquely soft
- ✓ Machine wash and dry
- ✓ Sold each

## STANDARD PROSTHETIC SOCK ORDERING

### INFORMATION NEEDED

1. Quantity
2. Length
3. Top and Toe Width (Flat)
4. Material Content

## SUPER SOCK® PROSTHETIC SOCK

| SIZE | 3-PLY      | 5-PLY      | TOP X TOE X LENGTH              |
|------|------------|------------|---------------------------------|
| 0    | 1K13-00-10 | —          | 6 x 3.50 x 10" (15 x 9 x 25cm)  |
| 0    | 1K13-00-12 | —          | 6 x 3.50 x 12" (15 x 9 x 30cm)  |
| 0    | 1K13-00-14 | 1K15-00-14 | 6 x 3.50 x 14" (15 x 9 x 35cm)  |
| 0    | 1K13-00-16 | 1K15-00-16 | 6 x 3.50 x 16" (15 x 9 x 40cm)  |
| 1    | 1K13-01-12 | 1K15-01-12 | 7 x 4 x 12" (17.50 x 10 x 30cm) |
| 1    | 1K13-01-14 | 1K15-01-14 | 7 x 4 x 14" (17.50 x 10 x 35cm) |
| 1    | 1K13-01-16 | 1K15-01-16 | 7 x 4 x 16" (17.50 x 10 x 40cm) |
| 1    | 1K13-01-18 | 1K15-01-18 | 7 x 4 x 18" (17.50 x 10 x 45cm) |
| 1    | 1K13-01-20 | 1K15-01-20 | 7 x 4 x 20" (17.50 x 10 x 50cm) |
| 1    | 1K13-01-22 | 1K15-01-22 | 7 x 4 x 22" (17.50 x 10 x 50cm) |
| 1    | 1K13-01-24 | —          | 7 x 4 x 24" (17.50 x 10 x 60cm) |
| 2    | 1K13-02-10 | 1K15-02-10 | 8 x 5 x 10" (20 x 12.5 x 25cm)  |
| 2    | 1K13-02-12 | 1K15-02-12 | 8 x 5 x 12" (20 x 12.5 x 30cm)  |
| 2    | 1K13-02-14 | 1K15-02-14 | 8 x 5 x 14" (20 x 12.5 x 35cm)  |
| 2    | 1K13-02-16 | 1K15-02-16 | 8 x 5 x 16" (20 x 12.5 x 40cm)  |
| 2    | —          | 1K15-02-18 | 8 x 5 x 18" (20 x 12.5 x 45cm)  |
| 2    | —          | 1K15-02-20 | 8 x 5 x 20" (20 x 12.5 x 50cm)  |
| 3    | —          | 1K15-03-10 | 9 x 6 x 10" (22.5 x 15 x 25cm)  |
| 3    | 1K13-03-12 | 1K15-03-12 | 9 x 6 x 12" (22.5 x 15 x 30cm)  |
| 3    | 1K13-03-14 | 1K15-03-14 | 9 x 6 x 14" (22.5 x 15 x 35cm)  |
| 3    | 1K13-03-16 | 1K15-03-16 | 9 x 6 x 16" (22.5 x 15 x 40cm)  |
| 3    | —          | 1K15-03-18 | 9 x 6 x 18" (22.5 x 15 x 45cm)  |
| 3    | —          | 1K15-03-20 | 9 x 6 x 20" (22.5 x 15 x 50cm)  |

## SUPER SOCK® PROSTHETIC SOCK WITH HOLE-IN-TOE

| SIZE | 3-PLY      | 5-PLY      | TOP X TOE X LENGTH             |
|------|------------|------------|--------------------------------|
| 1    | 1K13-S1-12 | 1K15-S1-12 | 7 x 4 x 12" (17.5 x 10 x 30cm) |
| 2    | 1K13-S2-12 | 1K15-S2-12 | 8 x 5 x 12" (20 x 12.5 x 30cm) |
| 2    | 1K13-S2-14 | —          | 8 x 5 x 14" (20 x 12.5 x 35cm) |
| 3    | 1K13-S3-14 | 1K15-S3-14 | 9 x 6 x 14" (23 x 15 x 35cm)   |

# STERLING STRETCH PROSTHETIC SOCK

- ✓ COOLMAX®/LYCRA® plaited to the inside of the high-bulk acrylic sock
- ✓ Conforms to the shape of the limb
- ✓ Wicks moisture away from the skin
- ✓ Fleeing on the inside adds softness next to the skin
- ✓ Contour and cushion the limb
- ✓ Available in 3-ply and 5-ply

## CARE TIPS

Machine wash warm. Do not use bleach. Tumble dry on medium.

3-Ply and 5-Ply Sterling Stretch Prosthetic Socks



## STERLING STRETCH PROSTHETIC SOCK

| 3-PLY      | 5-PLY      | SIZE           | FITS SOCK SIZE | TOP X TOE X LENGTH                                     |
|------------|------------|----------------|----------------|--|
| 1SS3-NA-SH | 1SS5-NA-SH | Narrow Short   | 00, 0, AS, AL  | 4-6" x 3-3.50" x 10-12" (10-15 x 7.6-9 x 25-30cm)      |
| 1SS3-NA-RG | 1SS5-NA-RG | Narrow Regular | 00, 0, AS, AL  | 4-6" x 3-3.50" x 14-16" (10-15 x 7.6-9 x 35-40cm)      |
| 1SS3-NA-LG | 1SS5-NA-LG | Narrow Long    | 00, 0, AS, AL  | 4-6" x 3-3.50" x 18-20" (10-15 x 7.6-9 x 45-50cm)      |
| 1SS3-MD-XS | 1SS5-MD-XS | Medium X-Short | 1, 2           | 7-8" x 4-5" x 8-10" (17.5-20 x 10-13 x 20-25cm)        |
| 1SS3-MD-SH | 1SS5-MD-SH | Medium Short   | 1, 2           | 7-8" x 4-5" x 10-12" (17.5-20 x 10-13 x 25-30cm)       |
| 1SS3-MD-RG | 1SS5-MD-RG | Medium Regular | 1, 2           | 7-8" x 4-5" x 14-16" (17.5-20 x 10-13 x 35-40cm)       |
| 1SS3-MD-LG | 1SS5-MD-LG | Medium Long    | 1, 2           | 7-8" x 4-5" x 18-20" (17.5-20 x 10-13 x 45-50cm)       |
| 1SS3-WD-SH | 1SS5-WD-SH | Wide Short     | 3, 4           | 9-10" x 6-7" x 10-12" (23-25.5 x 15-17.5 x 25-30cm)    |
| 1SS3-WD-MD | 1SS5-WD-MD | Wide Medium    | 3, 4           | 9-10" x 6-7" x 14-16" (23-25.5 x 15-17.5 x 35-40cm)    |
| —          | 1SS5-WD-LG | Wide Long      | 3, 4           | 9-10" x 4.50-5" x 18-20" (23-25.5 x 11.5-13 x 45-50cm) |
| —          | 1SS5-WD-XL | Wide X-Long    | 3, 4           | 9-10" x 4.50-5" x 22-24" (23-25.5 x 11.5-13 x 55-60cm) |



5-Ply Sterling Stretch Prosthetic Sock with Hole-In-Toe

## THICKNESS COLOR CODE (TOP SEAM)



## STERLING STRETCH PROSTHETIC SOCK WITH HOLE-IN-TOE

### STERLING STRETCH PROSTHETIC SOCK WITH HOLE-IN-TOE

| 3-PLY      | 5-PLY      | SIZE              | FITS SOCK SIZE | TOP X TOE X LENGTH                                     |
|------------|------------|-------------------|----------------|--|
| 1HS3-NA-SH | 1HS5-NA-SH | Narrow Short      | 00, 0, AS, AL  | 4-6" x 3-3.50" x 10-12" (10-15 x 7.6-9 x 25-30cm)      |
| 1HS3-MD-SH | 1HS5-MD-SH | Medium Short      | 1, 2           | 7-8" x 4-5" x 10-12" (17.5-20 x 10-13 x 25-30cm)       |
| 1HS3-MD-RG | 1HS5-MD-RG | Medium Regular    | 1, 2           | 7-8" x 4-5" x 14-16" (17.5-20 x 10-13 x 35-40cm)       |
| 1HS3-MD-LG | 1HS5-MD-LG | Medium Long       | 1, 2           | 7-8" x 4-5" x 18-20" (17.5-20 x 10-13 x 45-50cm)       |
| 1HS3-WD-SH | 1HS5-WD-SH | Wide Short        | 3, 4           | 9-10" x 6-7" x 10-12" (23-25.5 x 15-17.5 x 25-30cm)    |
| 1HS3-WD-MD | 1HS5-WD-MD | Wide Medium       | 3, 4           | 9-10" x 6-7" x 14-16" (23-25.5 x 15-17.5 x 35-40cm)    |
| 1HS3-WD-LG | 1HS5-WD-LG | Wide Long         | 3, 4           | 9-10" x 4.50-5" x 18-20" (23-25.5 x 11.5-13 x 45-50cm) |
| 1HS3-EW-MD | 1HS5-EW-MD | Extra Wide Medium | 5, 6           | 11-12" x 7-8" x 16-18" (28-30.5 x 17.5-20 x 40-45cm)   |



# STANDARD SYMES PROSTHETIC SOFT-SOCK



STANDARD SYMES  
SOFT-SOCK


The standard size Symes Soft-Sock will provide a well fitting sock for most amputees, but if special sizes are needed, please refer to pages 23- 25.

- ✓ Provides a well-fitting sock for most amputees
- ✓ Special sizes available
- ✓ Sold each.

## STANDARD SYMES PROSTHETIC SOFT-SOCK

| SIZE | 3-PLY      | 5-PLY      | TOP x ANKLE x BALL               |
|------|------------|------------|----------------------------------|
| 0    | 1SP3-A0-18 | 1SP5-A0-18 | 6" x 3" x 3.50" (15 x 7.5 x 9cm) |
| 0    | 1SP3-A0-20 | 1SP5-A0-20 | 6" x 3" x 3.50" (15 x 7.5 x 9cm) |
| 0    | 1SP3-A0-22 | 1SP5-A0-22 | 6" x 3" x 3.50" (15 x 7.5 x 9cm) |
| 0    | 1SP3-A0-24 | 1SP5-A0-24 | 6" x 3" x 3.50" (15 x 7.5 x 9cm) |

## THICKNESS COLOR CODE (TOP SEAM + LOGO)

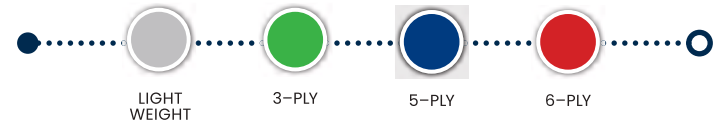


## STANDARD SOCK SIZING – APPROXIMATE RANGE

| SIZE | ANKLE WIDTH | BALL WIDTH   | ANKLE CIRCUMFERENCE      | BALL CIRCUMFERENCE       |
|------|-------------|--------------|--------------------------|--------------------------|
| A    | 3" (7.5cm)  | 3.50" (9cm)  | 5"-6.50" (12.5-16.5cm)   | 6.50"-8" (16.5-20cm)     |
| B    | 3" (7.5cm)  | 4" (10cm)    | 5"-6.50" (12.5-16.5cm)   | 7.50"-9.50" (19-22.5cm)  |
| C    | 3.50" (9cm) | 4" (10cm)    | 6"-8" (15-20cm)          | 8.50"-9.50" (20-24cm)    |
| D    | 3.50" (9cm) | 4.50" (11cm) | 6"-8" (15-20cm)          | 8.50"-10" (21-25cm)      |
| E    | 4" (10cm)   | 4.50" (11cm) | 7.50"-10.50" (19-26.5cm) | 9-10.50" (22.5-26.5cm)   |
| H    | 5" (12.5cm) | 5.50" (14cm) | 9.50"-11" (24-28cm)      | 10.50"-12" (26.5-30.5cm) |

The toe width at the very tip end of the Symes Sock is wider (on both standard Symes sizes and special Symes sizes) than the width at the very tip end of regular B/K size prosthetic socks.

Use of regular B/K size socks for Symes amputees may result in wrinkles at the bottom of the limb and/or break through the sock near the toe.

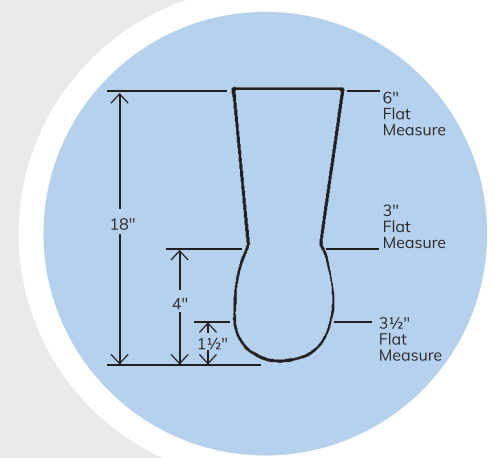

Example of sizing 1S15-A0-18:  
18" x 6" x 3" x 3.50"  
Graduated for larger sizes.

# HIP DISARTICULATION PROSTHETIC SOCK

- ✓ Seamless
- ✓ Made in Poly/Ionic+™ and Soft-Sock fabric.
- ✓ Waist has a turned over hem, plus a 3" leg extension
- ✓ Special features available for hip disarticulation socks
- ✓ Sold each


Special features available for hip disarticulation socks are non-standard sizing, hemi-pelvectomy, fly openings (state whether right or left amputation) and elastic for leg or waist.

Specify if bilateral hip disarticulation or hemi-pelvectomy sock is desired.

The flat measurement width of standard sizes is listed. Extra length and/or width is available for an additional charge.

## HIP DISARTICULATION PROSTHETIC SOCK STANDARD CONVENTIONAL UNILATERAL

| FLAT WIDTH | LENGTH     | 3-PLY<br>SOFT-SOCK | 5-PLY<br>SOFT-SOCK | 3-PLY<br>IONIC+™ |
|------------|------------|--------------------|--------------------|------------------|
| 8" (20cm)  | 12" (30cm) | —                  | —                  | 1HX3-08-12       |
| 10" (25cm) | 12" (30cm) | —                  | —                  | 1HX3-10-12       |
| 12" (30cm) | 14" (36cm) | —                  | —                  | 1HX3-12-14       |
| 12" (30cm) | 16" (40cm) | 1HP3-12-16         | 1HP5-12-16         | —                |
| 13" (33cm) | 14" (36cm) | —                  | —                  | 1HX3-13-14       |
| 13" (33cm) | 16" (40cm) | 1HP3-13-16         | 1HP5-13-16         | —                |
| 14" (36cm) | 14" (36cm) | —                  | —                  | 1HX3-14-14       |
| 14" (36cm) | 16" (40cm) | 1HP3-14-16         | 1HP5-14-16         | —                |
| 15" (38cm) | 16" (40cm) | 1HP3-15-16         | 1HP5-15-16         | 1HX3-15-16       |
| 16" (40cm) | 16" (40cm) | 1HP3-16-16         | 1HP5-16-16         | 1HX3-16-16       |
| 17" (44cm) | 16" (40cm) | 1HP3-17-16         | 1HP5-17-16         | 1HX3-17-16       |
| 18" (46cm) | 16" (40cm) | 1HP3-18-16         | 1HP5-18-16         | —                |
| 18" (46cm) | 18" (46cm) | —                  | —                  | 1HX3-18-18       |
| 19" (48cm) | 16" (40cm) | 1HP3-19-16         | 1HP5-19-16         | —                |
| 19" (48cm) | 18" (46cm) | —                  | —                  | 1HX3-19-18       |

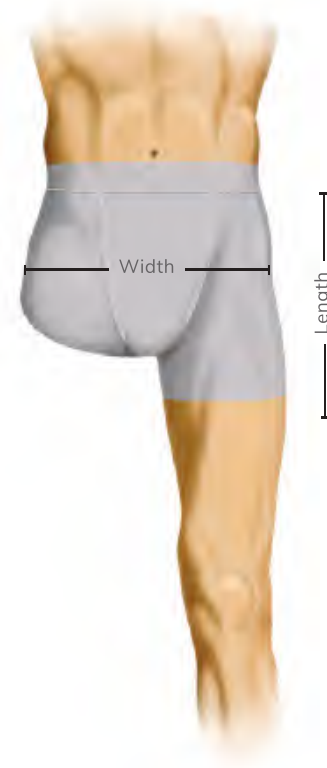

## SIZING TIPS

To determine which flat measurement of hip socks to order, take 10–20% less than half the actual circumference hip measurement in order to accommodate both the hip and the waist, and make the best use of the stretch of the fabric.

Each size will fit circumferences smaller than those stated if stretched in the length, or larger circumferences if stretched in the width.

For instance, a 14" width flat measurement, 16" length sock would fit a tall person with a 26" hip as well as a short person with up to a 32" hip measurement.

Stretching in the width shortens the length.

Stretching in the length reduces the width.


Allowances should be made in sizing.

Because of the increased elasticity of these fabrics, flat widths of the socks will be up to 10% less than measurement stated.

# PARTIAL FOOT PROSTHETIC SOCK



US Patent  
#6,158,253

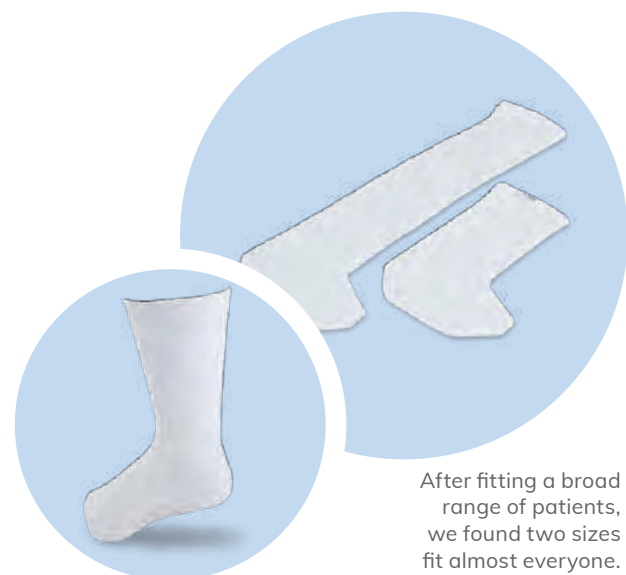

After fitting a broad  
range of patients,  
we found two sizes  
fit almost everyone.

## SMARTKNIT® PARTIAL FOOT SOCK

| WHITE | BLACK/IONIC+™ | WIDTH / LENGTH            |
|-------|---------------|---------------------------|
| 84107 | 84707         | 3.50" x 14.50" (9 x 37cm) |

- ✓ No heel assures custom fit
- ✓ For all levels of partial foot amputation
- ✓ Completely seamless
- ✓ Knit with core-spun yarns
- ✓ Fits like a second skin
- ✓ Multi-directional stretch for smooth, comfortable fit
- ✓ Breathable, washable COOLMAX®/LYCRA® in white and Ionic+™ in black
- ✓ Wicks moisture from the skin
- ✓ Antimicrobial to inhibit odor in the sock
- ✓ Non-binding Halo™ top
- ✓ One size fits all
- ✓ Sold each

## SMARTKNIT® PARTIAL FOOT SOCK WITH HEEL

| WHITE | BLACK/IONIC+™ | STYLE         | WIDTH / LENGTH           |
|-------|---------------|---------------|--------------------------|
| 74103 | 74703         | Crew          | 3.50" x 10" (9 x 25cm)   |
| 74105 | 74705         | Over-The-Calf | 3.75" x 18" (9.5 x 46cm) |

## PARTIAL FOOT SOCK WITH HEEL

| ITEM # | SIZE | FOOT WIDTH / TOP / LENGTH |
|--------|------|---------------------------|
|--------|------|---------------------------|

### LIGHTWEIGHT SOFT-SOCK – COOLMAX/LYCRA®

|            |       |                                |
|------------|-------|--------------------------------|
| 1SP1-4S-10 | Short | 6" x 4" x 10" (15 x 10 x 25cm) |
| 1SP1-4S-20 | Long  | 6" x 4" x 20" (15 x 10 x 50cm) |

### 2 PLY A-PLUS – ACRYLIC/LYCRA®

|            |       |                                   |
|------------|-------|-----------------------------------|
| 1KL2-4S-10 | Short | 6" x 4.50" x 10" (15 x 11 x 25cm) |
| 1KL2-4S-20 | Long  | 6" x 4.50" x 20" (15 x 11 x 50cm) |

- ✓ Completely seamless sock
- ✓ Smooth and soft
- ✓ Lightweight Soft-Sock or 2-ply A-Plus available
- ✓ Contours without wrinkles
- ✓ Wicks moisture to keep the skin dry
- ✓ Machine washable
- ✓ Sold each.

## CARE TIPS

Machine wash  
and tumble dry  
warm. No bleach.

# POST-OP PROSTHETIC SOCKS AND SUPPLIES

## SILVER A/K AND B/K POST-OP SOCKS

PROSTHETIC INTERFACE FOR  
IMMEDIATE POST-OPERATIVE CARE.\*

| ITEM #     | LENGTH     | TOP WIDTH   | TOE WIDTH   |
|------------|------------|-------------|-------------|
| 1PL3-X5-18 | 18" (45cm) | 6" (15cm)   | 5" (12.5cm) |
| 1PL3-X6-18 | 18" (45cm) | 8" (20cm)   | 6" (15cm)   |
| 1PL3-X7-18 | 18" (45cm) | 9" (22.5cm) | 7" (17.5cm) |

- ✓ Made with 10% Ionic+™ Antimicrobial Technology.
- ✓ Packed in individual polyethylene bags, four (4) to a box
- ✓ Rolled for application
- ✓ 87.5% Polyester, 10% Ionic+™ Silver, 2.5% LYCRA®
- ✓ Patent Pending



Third party testing showed 99.98% reduction of bacteria directly on the surface of the sock (versus 93% for control sample)\*\*

## ORIGINAL UNSTERILIZED A/K AND B/K POST-OP SOCKS

PROSTHETIC INTERFACE FOR  
IMMEDIATE POST-OPERATIVE CARE.\*

| ITEM #     | LENGTH     | TOP WIDTH   | TOE WIDTH   |
|------------|------------|-------------|-------------|
| 1PL3-BS-18 | 18" (45cm) | 5" (12.5cm) | 4" (10cm)   |
| 1PL3-05-18 | 18" (45cm) | 6" (15cm)   | 5" (12.5cm) |
| 1PL3-26-18 | 18" (45cm) | 8" (20cm)   | 6" (15cm)   |
| 1PL3-37-18 | 18" (45cm) | 9" (22.5cm) | 7" (17.5cm) |

- ✓ Packed in individual polyethylene bags, four (4) to a box
- ✓ 97% Acrylic, 3% LYCRA®
- ✓ Modified box toe
- ✓ Rolled for application



Original Unsterilized A/K and B/K Post-Op Socks

- ✓ 3-Ply weight
- ✓ Fleeced inside
- ✓ Latex free materials

**NOTE:** Socks knitted from these materials are sized to fit per sizing of other materials such as wool. Due to the stretch and contraction of the material the flat measurements will be approximately 10% less in both length and width.



## SYMES AND CHOPART UNSTERILIZED POST-OP SOCKS

**WARNING — MUST BE GAS AUTOCLAVED TO BE STERILE.**

| ITEM #         | LENGTH     | TOP WIDTH    | ANKLE WIDTH  | BALL        |
|----------------|------------|--------------|--------------|-------------|
| <b>SYMES</b>   |            |              |              |             |
| 1SL3-E1-24     | 24" (60cm) | 7.50" (19cm) | 4.50" (11cm) | 5" (12.5cm) |
| <b>CHOPART</b> |            |              |              |             |
| 1CL3-E1-24     | 24" (60cm) | 7" (18cm)    | 4" (10cm)    | 6" (15cm)   |

- ✓ For immediate post operative use
- ✓ Seamless for optimum fitting qualities
- ✓ 97% acrylic, 3% LYCRA®
- ✓ Sold each

\*Not a primary sterile dressing

\*\* Data on file

Symes and Chopart  
Post-Op Socks



## POST-OP SUPPLIES

- ✓ Reticulated, Polyester Urethane Foam, with 20 pores per lineal inch, (PPI)
- ✓ 97% void
- ✓ Weighs only 2 pounds per cubic foot
- ✓ Inert, and in itself will not support fungus growth
- ✓ Most suitable for the intended purpose as an interface agent in the rigid dressing
- ✓ Greatly simplifies interface material application since it gives uniformity to all users
- ✓ Latex free materials

**NOTE — MAY BE STERILIZED IN THE SAME MANNER AS THE POST-OP SOCKS OR BY STEAM AUTOCLAVE**



### DISTAL END PADS

| ITEM #     | SIZE                |
|------------|---------------------|
| 3M-2RSP-ST | Felt Pads, Set of 3 |



### LAMBS WOOL

| ITEM #     | SIZE  |
|------------|-------|
| 2M-S0LW-10 | 1 oz. |

**GREY IONIC+™  
STYLES CONTAIN  
IONIC+™ MINERAL  
ANTIMICROBIAL  
FABRIC KNOWN FOR  
ITS MOISTURE WICKING  
PROPERTIES**



Ionic+™ is made by permanently bonding 99.9% pure silver to the surface of fibers. Socks powered by Ionic+™ inhibit the growth of bacteria and fungi directly on the surface of the fabric. Because it is permanent, it will not wash out and performs for the life of the product.

Ionic+™ technology bonds a layer of pure silver to the surface of a textile fiber, permanently retaining the best characteristics of both.

Garments made with Ionic+™ silver fibers receive rave performance reviews and enjoy widespread usage in medicine, the US Military, NASA astronauts and over 60 Olympic teams.

See page 38 for Post-Op  
Suspension Belts



# MATERIALS AND CARE INSTRUCTIONS FOR PROSTHETIC AND ORTHOTIC SOCKS AND INTERFACES

This section lists the entire Knit-Rite family of orthotic and prosthetic sock materials. Because the textile interface is a vital component of the orthotic/prosthetic device system, no one sock can meet all the needs of wearers of prosthetic and orthotic devices.

Our carefully designed interfaces vary in knit construction, sizing, weight, fiber content, and special features to provide the ultimate selection of textile interfaces for your patients. Selection of the proper socks contributes significantly to

the comfort, function and mobility of the individual, and the performance of the device. You may find that different members of the Knit-Rite sock family fulfill different needs as individuals change jobs, activities, age, gain or lose weight.

Initially, a planned supply for a new amputee may include several kinds of socks. As patients are re-evaluated, the type of socks that would best meet the current needs of the individual should be reconsidered.

## 100% WOOL SUPER SOCKS

The most universally accepted prosthetic socks. Knit-Rite uses a blend of the finest wool fibers available to make the socks soft, resilient, absorbent, thick and protective. Easy care Super Socks are available in 3, 5 or 6-ply.

## ACRYLIC/LYCRA® (A/L) SOCKS

These soft, stretchable socks provide a more intimate fit with fewer wrinkles. 2-ply approximates 3-ply wool in thickness; 3-ply approximates 5-ply wool; 4-ply approximates 6-ply wool.

Stretch accommodates irregular shapes and congenital deformities. Flat measurements of A/L leg/arm orthotic interfaces will be approximately 10% less in width than the stated size to adjust for garment stretch.

## SOFT-SOCK® WITH COOLMAX® AND SOFT-SOCK® WITH IONIC+™

These lightweight, 3-ply, 5-ply, and 6-ply stretchable socks add comfort by wicking perspiration away from the skin and inhibiting odor-causing bacteria.

## X-WOOL® SOCKS

The resiliency of wool, the stretch of LYCRA® spandex, and the high performance advantages of Ionic+™ – Mineral Antimicrobial fit a greater size range; accommodate flexion; resist shrinking and felting. 3-, 5- and 6-ply.

## KNIT-RITE LAUNDERING INSTRUCTIONS

| FABRIC                      | LOAD  | WASH         |                            |                                     | RINSE       |            | DRY           |                                      |
|-----------------------------|---|--------------|----------------------------|-------------------------------------|-------------|------------|---------------|--------------------------------------|
|                             |   | TIME         | TEMP                       | PRODUCT                             | TIME        | TEMP       | TIME          | TEMP                                 |
| WOOL SUPER SOCK®<br>X-WOOL® | Turn transfer side in; wash with white laundry or as wool | 5–10 Minutes | Warm 100°F<br>Excess water | 1/2 Cup detergent<br>NO BLEACH      | 3–5 Minutes | Warm 100°F | 30–40 Minutes | Tumble Dry Low or Hang Dry           |
| STERLING STRETCH            | White laundry. Machine wash sock right side out           | 8–10 Minutes | Warm on gentle cycle       | 1/2 Cup mild detergent<br>NO BLEACH | 3–5 Minutes | Warm 100°F | 20–30 Minutes | Tumble Dry Medium or permanent press |
| COTTON                      | White laundry   | 8–10 Minutes | Hot 130°F                  | 1/2 Cup detergent<br>Bleach O.K.    | 3–5 Minutes | Warm 100°F | 40–50 Minutes | Tumble Dry High                      |

## CARE TIPS

1. Wear a fresh sock and sheath each day.
2. Pre-spot socks in cool water.
3. Turn transfer label to inside for wool; leave sewn-in label inside for Acrylic/LYCRA®, cotton and Soft-Sock.
4. Wash according to chart above.
5. After washing, take tip of toe in one hand, top in other hand and pull taut to smooth sock before drying.

# PROSTHETIC SOCK ORDERING GUIDELINES

## KNIT-RITE STANDARD PROSTHETIC SOCK SIZING

| SIZE  | TOP WIDTH    | TOE WIDTH   |
|-------|--------------|---|
| NO. 0 | 6" (15cm)    | 3.50" (9cm)                                       |
| NO. 1 | 7" (17.5cm)  | 4" (10cm)   |
| NO. 2 | 8" (20cm)    | 5" (12.5cm) [over 24" (60cm) length 4" (10cm)]    |
| NO. 3 | 9" (23cm)    | 6" (15cm) [over 20" (50cm) length 5.50" (14cm)]   |
| NO. 4 | 10" (25cm)   | 7" (17.5cm) [over 20" (50cm) length 6" (15cm)]    |
| NO. 5 | 11" (27.5cm) | 7" (17.5cm) [over 18" (45cm) length 5.50" (14cm)] |
| NO. 6 | 12" (30cm)   | 8" (20cm) [over 18" (45cm) length 6" (15cm)]      |

## KNIT-RITE STANDARD ARM SOCK SIZING

| SIZE | TOP WIDTH   | TOE WIDTH  |
|------|-------------|------------|
| A    | 4" (10cm)   | 3" (7.5cm) |
| B    | 5.5" (14cm) | 3.5" (9cm) |

\*These sizes shipped unless otherwise specified.

## HCPCS BILLING CODES

| CODES | DESCRIPTIONS   |
|-------|--|
| L0984 | Protective Body Sock, each   |
| L2820 | Addition to lower extremity orthosis soft interface for molded plastic, below the knee section (Knit-Rite AFO Liner Sock).   |
| L2840 | Addition to lower extremity orthosis, tibial length sock, fracture or equal, each  |
| L2850 | Addition to lower extremity orthosis, sock, femoral length sock, fracture or equal, each   |
| L3995 | Addition to upper extremity orthosis, sock, fracture or equal, each  |
| L8400 | Prosthetic Sheath, below knee, each  |
| L8410 | Prosthetic Sheath, above knee, each  |
| L8415 | Prosthetic Sheath, upper limb, each  |
| L8420 | Prosthetic Sock, multiple ply, below knee, each*   |
| L8430 | Prosthetic Sock, multiple ply, above knee, each*   |
| L8435 | Prosthetic Sock, multiple ply, upper limb, each*<br>*Traditional Knit-Rite Prosthetic Socks with knit-in taper, including Knit-Rite Liner Liner. Multiple (1-8) plies of natural or synthetic yarns. Knit on flat bed knitting machines. |
| L8440 | Prosthetic Shrinker, below knee, each  |
| L8460 | Prosthetic Shrinker, above knee, each  |
| L8465 | Prosthetic Shrinker, upper limb, each  |
| L8470 | Fitting Sock, single ply, below knee, each   |
| L8480 | Fitting Sock, single ply, above knee, each   |
| L8485 | Fitting Sock, single ply, upper limb, each<br>Fitting Socks, single ply, tubular construction. Knit on circular knitting machines.   |
| L7600 | Prosthetic Donning Sleeve, any material, each  |

We have included the above codes to assist the provider in determining and obtaining proper reimbursement for our products. It is the provider's responsibility to ascertain the appropriateness of individual applications of these billing codes.



## GENERAL MEASUREMENT INSTRUCTIONS

To provide the proper fitting prosthetic sock for your patient, it is important that a correct measurement be made. There are three basic measurements that are necessary for a proper fitting prosthetic sock: length, top width and toe width.

### A-B MEASUREMENT (TOP WIDTH)

Should be taken from top left corner of sock lying flat to top right corner and correspond to 1/2 the circumference measurement at the point where the top of the sock should be on the leg or arm.

### C-D MEASUREMENT (LENGTH)

Should be taken from top of sock lying flat to bottom of toe and should correspond to the length from distal to where top of sock should come.

### E-F MEASUREMENT (TOE WIDTH)

Should be taken from left side of sock lying flat, perpendicular to toe fashionings, to right side of sock and should correspond to 1/2 the circumference 2" up from the distal end of the limb.

## REMEMBER

A proper measurement means a correct fitting prosthetic sock for your patient.

**NOTE:** AFO and KAFO socks found on pages 44–45. "Fits Circumferences" guidelines for Torso-Interface, found on pages 48–49. "Approximate Calf Circumference" fitting guidelines for fracture brace socks found on pages 52–53.

Refer to the Index for pages on Symes, Hip Socks, X-Wool, and Soft-Sock

# SPECIAL AND IRREGULAR SOCK ORDERING GUIDELINES

## PROSTHETIC SOCK SIZES

Nearly 90% of the amputees can be correctly fitted with a stock size sock because the knitted construction has a stretch which allows a fitting range. Most prosthetists can judge the needed size by working with the limb and the cast. Nevertheless, some measurement guidelines can be helpful. First, measure the limb circumference at the top of the socket; second, measure the toe limb circumference 2 inches from the distal end; third, measure from the distal end to the top of the socket. The length of the sock should be 2 to 4 inches more than this measurement depending on the type of fitting and suspension preference. (Wool and A-Plus will shrink more than Soft Sock.) The patient should not have to pull the sock too taut.

## SOCKS ARE SIZED IN FLAT MEASUREMENTS

i.e. half the circumference measurement (less 5% to give a snug fit) for top and toe—then matching those figures to the closest comparable size measurements in all the charts—see page 22. Custom socks need to be ordered only when patient's measurements do not correspond with any of the coded stock sizes or when a sock needs an addition or adjustment, i.e. split for hip joint, suspension system holes, elastic at waist or leg, fly, etc. (see page 24.)

## HIP JOINT/PELVIC BELT PROSTHETIC SOCKS

An order for socks with a split for hip joint should state the size of the sock and whether or not a tie or just the edged split is desired. Slits are 3 inches deep unless otherwise specified.

## ISCHIAL CONTAINMENT / HIGH LATERAL BRIM

To eliminate excessive rollover, medial cuts for sockets with high lateral wall can be made when requested. Specify the depth of the cut desired.

## IRREGULAR SHAPES

For prosthetic socks which are not evenly graduated from the top to the toe (as those pictured on next page), take circumference measurements at the irregular points, lengths from the distal end to each of those points and determine total length. Sending a drawing is very helpful. State fabric and ply choice: Super Sock, A-Plus, Soft-Sock, X-Wool, Sterling Wool and Sterling Stretch.

## FITTING SOCKS AND RE-ORDERS

Please state patient's I.D. on all orders for custom socks. A "fitting" sock will normally be sent. It is suggested that only two socks be ordered until the right size and material can be established. Then the quantity can be specified. Records for each customer are kept and reorders can duplicate the original or alterations can be made, as needed, usually without another pattern.

## ORTHOTIC SOCK SIZES

Although socks to interface orthotic appliances can be made of any material, most are made of yarn with a LYCRA® spandex core and an acrylic, COOLMAX® or cotton covering. The extra stretch of these fabrics allows wide fitting ranges.

## TAILORING

When socks need to be very irregular, an outline pattern with measurements every two inches can be sent or an exact pattern can be duplicated. A single flat seam will be used in construction with all edges being finished. Tailoring of the knitted material can produce socks that are impossible to fashion on the knitting machines, especially when rounding is desired. If knitted socks or tailored socks are preferred, please specify.

### 1. & 2. SYMES AND CHOPART

Symes and Chopart socks have a bulbous end. Standard sizes listed on pages 16. If a custom size is needed, prepare a profile drawing of the limb, giving circumference measurements at the top, ankle and ball. Note on drawing if limb is symmetrical or asymmetrical. If asymmetrical, record a circumference measurement of the partial foot extension.

### 3. KNEE DISARTICULATION (BOX TOE)

State the length desired, the top circumference at the perineum and the toe circumference 1 inch from the end of the limb. Specify box toe.

### 4. LARGE KNEE AND THIGH

Send a drawing of the limb, front and back view. Specify length of sock with allowance for turn down. State length from the end to the beginning of the bulge below the knee and from the end to the largest part of the knee. State circumferences at top, knee, below condyles and 2 inches from end of the limb.

### 5. FLEXED KNEE

Prepare a profile drawing. Specify the length of the sock. Specify circumferences at top, above knee, around knee, and at the end of the limb. If flexion angle is desired, please specify.

### 6. ABOVE ELBOW ARM SOCKS

Prepare a drawing of the front view. Note the total length and the length from the axilla to the distal end of the limb. Take circumference measurements 1" from the distal end and at the axilla. Specify a flat measurement from the upper arm to the point on the chest where the sock should come.

### 7. HIP DISARTICULATION AND HEMI-PELVECTOMY


Prepare a drawing as illustrated. Specify length from groin to waist. Hip disarticulation socks are made with turnover hem at waist. Give circumference measurements for waist, hips, and sound leg (take at angle). If fly or colostomy hole is desired, specify which leg has been amputated. Hip socks are made with 3" leg cuffs, unless otherwise specified. Standard sizes listed on page 16.




# SPECIAL AND IRREGULAR SOCK ORDERING GUIDELINES

## SPECIFY LENGTHS BETWEEN CIRCUMFERENCE MEASUREMENTS ON ALL CUSTOM SOCKS.


### 1. CHOPART

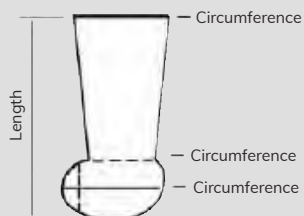


### 2. SYMES

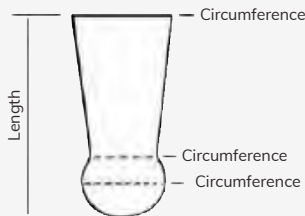


### 3. KNEE BEARING

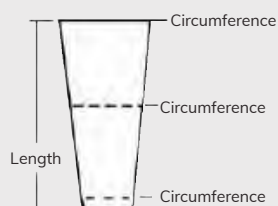


### 4. LARGE KNEE

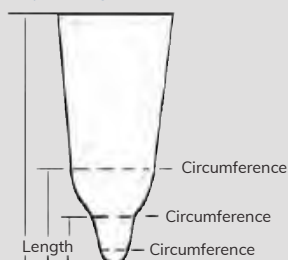


### 5. FLEXED KNEE

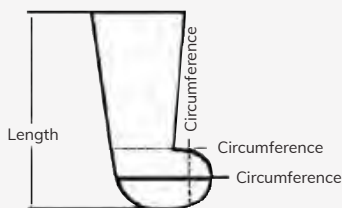

### 6. ABOVE ELBOW ARM SOCK

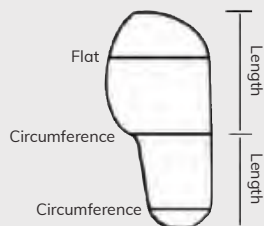

### 7. HIP DISARTICULATION

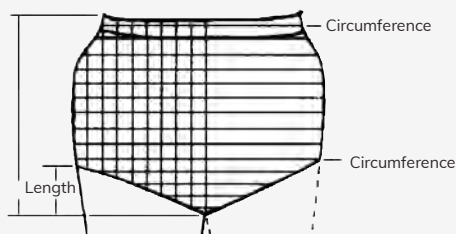

## SPECIAL ORDER FULFILLMENT

All orders will be filled in their entirety if at all possible. Custom orders take an estimated delivery time of 3 to 5 working days to make before shipping. If a trial fitting sock is desired to determine the correct size and material of a custom sock, only two socks should be ordered until the correct size and fabric are established. Then additional quantities can be ordered. Custom socks are not returnable.

## TAILORING FOR PROSTHETIC AND ORTHOTIC SOCKS

Tailoring socks is a technique developed to overcome the limitations imposed by the knitting machines. The knitted fabrics used for tailoring are the same as those of the machine fashioned socks; however, the shaping is a custom process to provide a smooth contoured fit not possible in some knitted socks. Tailored socks are cut and seamed to the exact shapes stated in measurements and/or drawings. Providing a front and side view pattern is helpful. Seams open flat and are less distinguishable in feel than the fashionings in knitting. Tailored socks may be made in Super Sock, A-Plus and Soft-Sock. Socks with special contours achieved by knitting are available in most instances if requested. Tailoring offers the prosthetist and orthotist a way to obtain a smoother, more even fit over irregular contours.

Tailored socks are cut and seamed to the exact shapes stated in measurements and/or drawings.



# SPECIAL ORDERS AND PRICING

## PRICING ON IRREGULAR AND SPECIAL PROSTHETIC SOCKS

The base price (of either the knitted or tailored sock) is the price of the coded standard prosthetic sock nearest to the irregular or special sock, plus all the additional modifications requested. Special prosthetic socks cannot be returned. Special prosthetic socks must be ordered in multiples of 2. The base price is determined from:

**A. TOP FLAT WIDTH AND LENGTH DIMENSIONS** as matched to standard sizes (page 22). On fractional top widths (i.e. when circumference doesn't divide evenly by 2) use the next lowest flat width dimension to allow for stretch. When top circumference is 20" or more, take 10% from the circumference measurement and divide by 2 for the flat measurement.

**B. MATERIAL AND PLY (WEIGHT)** Using (A) and (B), the computer selects the nearest regular sock code and price, then adds specified features.

BASE PRICE EACH \_\_\_\_\_  
 ADDITIONAL FEATURES \_\_\_\_\_  
 \_\_\_\_\_  
 \_\_\_\_\_  
 TOTAL \_\_\_\_\_

### DESCRIPTION

#### EACH ADDED TO BASE PRICE

Knitted or tailored to dimensions other than standard sizes, but having a top width of 12" or less (flat measurement) and a length of less than standard length. This includes box toe and non-standard brace socks

Knitted or tailored with pelvic split or high lateral wall

Distal hole, non-standard size

Distal hole, standard size if glue reinforced

Folded and sewn top edge

Any alteration to a knitted or tailored sock

Knitted or tailored with pelvic split and ties

Knitted or tailored and given extra processing or fleecing

Knitted or tailored sock over 12" in width (flat measurement) and greater than standard length

Knitted sock reinforced at toe or complete sock

Hip sock over 19" in width (flat measurement) and/or over 16" in length. Tailored hip sock with leg extension; other than coded material. (This charge may occur more than once.)

Hip sock with colostomy—specify right or left amputation

Hip sock with fly—specify right or left amputation, and whether or not velcro closing is desired:

Right ☐ Left ☐ Velcro ☐

Sock with elastic binding (if hip sock requests elastic binding at both waist and leg, this charge will occur twice)

Unhemmed (Sterling Socks)

Silicone-dotted band

### PLEASE SEND THE FOLLOWING INFORMATION FOR ALL IRREGULAR AND SPECIAL PROSTHETIC SOCKS

Facility \_\_\_\_\_

Address \_\_\_\_\_

City \_\_\_\_\_

State \_\_\_\_\_

Zip \_\_\_\_\_

Telephone \_\_\_\_\_

Contact Person At Facility \_\_\_\_\_

1. Patient's 1SPEC # or I.D. \_\_\_\_\_

2. Site Of Amputation \_\_\_\_\_

3. ☐ New patient

☐ Reorder

4. Quantity \_\_\_\_\_

5. Material \_\_\_\_\_

6. Ply (Weight) \_\_\_\_\_

7. Specify one

☐ 3, 5 or 6-ply Wool Super Sock

☐ 2, 3 or 4-ply acrylic/LYCRA®

☐ Lightweight, 3 or 5-ply COOLMAX®/LYCRA®

☐ Ionic+™/LYCRA® Soft-Sock

8. Measurement — (See guidelines on page 21.)

Record measurements in "inches" unless otherwise noted. Remakes due to wrong unit of measurement will be at customer's expense.

If reorder, state flat measurements as on old package (noting any modifications).

Length \_\_\_\_\_ Top \_\_\_\_\_ Toe \_\_\_\_\_

SPEC. No. \_\_\_\_\_

For new patient drawing (outline) with measurements provided:

☐ Yes ☐ No

Length \_\_\_\_\_ Top Circ. \_\_\_\_\_ Toe Circ. \_\_\_\_\_

For irregular shapes, also provide circumferences from toe:

At: 2" \_\_\_\_\_ 4" \_\_\_\_\_ 6" \_\_\_\_\_

8" \_\_\_\_\_ 10" \_\_\_\_\_ 12" \_\_\_\_\_

9. Additional features, special instructions, shipping:



# PROSTHETIC SHRINKERS

## PROSTHETIC SHRINKERS | 26-35

- 27** Compressogrip® B/K Prosthetic Shrinker
- 28** Tapered Compressogrip® B/K Prosthetic Shrinker
- 28** Bone Rings
- 29** Compressogrip® A/K Prosthetic Shrinker
- 30** Post-Op Dressing Retention Sock
- 31-33** 4-Way Stretch Shrinker
- 34** Custom Prosthetic Shrinkers
- 35** Private Logo Program



For custom logo shrinkers, see page 35, or ask customer service or your distributor for program details.

# COMPRESSOGRIP® B/K PROSTHETIC SHRINKER

THE ORIGINAL TUBULAR COMPRESSION PROSTHETIC SHRINKER USED FOR CONTROL OR PREVENTION OF EDEMA FOLLOWING AMPUTATION SURGERY, FOR LIMB VOLUMETRIC MAINTENANCE, TO MINIMIZE HYPERTROPHIC SCARRING AND TO GIVE COMPRESSIVE SUPPORT FOR PAIN RELIEF.

- ✓ Assure controlled compression, ease of application and wearer acceptance
- ✓ Efficacy clinically tested (reprints available)
- ✓ Ionic+™ adds permanent antimicrobial benefits to inhibit odor in the garment and effectively wicks moisture away from the skin
- ✓ Made from latex-free materials
- ✓ Washable

## MATERIALS

White, Beige — 92% polyester, 8% spandex

Polyester/Ionic+™ Grey — 87.4% polyester, 4.6% Ionic+™, 8% spandex

Compressogrip®  
B/K Elastic  
Shrinker

PROSTHETIC SHRINKERS



## FITTING INSTRUCTIONS



1. Apply the first layer.



2. Position the ring snugly.



3. Apply the second layer to a point of desired maximum compression.

## COMPRESSOGRIP® B/K ELASTIC SHRINKER MEASUREMENT OF B/K CIRCUMFERENCE 2" BELOW MEDIAL TIBIA TUBERCLE

| WHITE      | BEIGE      | GREY/<br>GREY IONIC+™ | SIZE | LENGTH     | FITS CIRCUMFERENCE | BONE RING SIZE |
|------------|------------|-----------------------|------|------------|--------------------|----------------|
| 2H0S-02-18 | 2H0B-02-18 | 2H-0X02-18            | #2   | 18" (45cm) | 5"-7" (13-18cm)    | 3/4"           |
| —          | —          | 2H-0X0S-18            | #2.5 | 18" (45cm) | 7"-10" (18-25cm)   | 3/4"           |
| 2H0S-0S-24 | —          | —                     | #2.5 | 24" (60cm) | 7"-10" (18-25cm)   | 3/4"           |
| 2H0S-03-18 | 2H0B-03-18 | 2H-0X03-18            | #3   | 18" (45cm) | 9"-12" (23-30cm)   | 7/8"           |
| 2H0S-03-24 | 2H0B-03-24 | 2H-0X03-24            | #3   | 24" (60cm) | 9"-12" (23-30cm)   | 7/8"           |
| 2H0S-03-30 | 2H0B-03-30 | —                     | #3   | 30" (75cm) | 9"-12" (23-30cm)   | 7/8"           |
| 2H0S-S3-18 | 2H0B-S3-18 | —                     | #3.5 | 18" (45cm) | 10"-13" (25-33cm)  | 7/8"           |
| 2H0S-S3-24 | 2H0B-S3-24 | —                     | #3.5 | 24" (60cm) | 10"-13" (25-33cm)  | 7/8"           |
| —          | —          | 2H-0XS3-30            | #3.5 | 30" (75cm) | 10"-13" (25-33cm)  | 7/8"           |
| 2H0S-04-18 | 2H0B-04-18 | 2H-0X04-18            | #4   | 18" (45cm) | 12"-16" (30-40cm)  | 7/8"           |
| 2H0S-04-24 | 2H0B-04-24 | 2H-0X04-24            | #4   | 24" (60cm) | 12"-16" (30-40cm)  | 7/8"           |
| 2H0S-04-30 | 2H0B-04-30 | 2H-0X04-30            | #4   | 30" (75cm) | 12"-16" (30-40cm)  | 7/8"           |
| 2H0S-05-18 | 2H0B-05-18 | 2H-0X05-18            | #5   | 18" (45cm) | 15"-20" (38-50cm)  | 1"             |
| 2H0S-05-24 | 2H0B-05-24 | 2H-0X05-24            | #5   | 24" (60cm) | 15"-20" (38-50cm)  | 1"             |
| 2H0S-05-30 | 2H0B-05-30 | 2H-0X05-30            | #5   | 30" (75cm) | 15"-20" (38-50cm)  | 1"             |
| 2H0S-06-24 | 2H0B-06-24 | 2H-0X06-24            | #6   | 24" (60cm) | 18"-24" (45-60cm)  | 1 1/8"         |
| 2H0S-06-30 | 2H0B-06-30 | 2H-0X06-30            | #6   | 30" (75cm) | 18"-24" (45-60cm)  | 1 1/8"         |


See page 28 for Bone Rings, and page 39 for Shrinker Donning Tubes.

NOTE: Compressogrip® shrinker fabric gives up to 20–25 mmHg compression at 50% stretch. At 100% stretch the compression is almost doubled. Double layers will almost double compression at both 50% and 100% stretch. Application of more than two layers may exceed the limits of desired therapeutic pressure.

When determining sizing, allow for swelling



# TAPERED COMPRESSOGRIP B/K PROSTHETIC SHRINKER

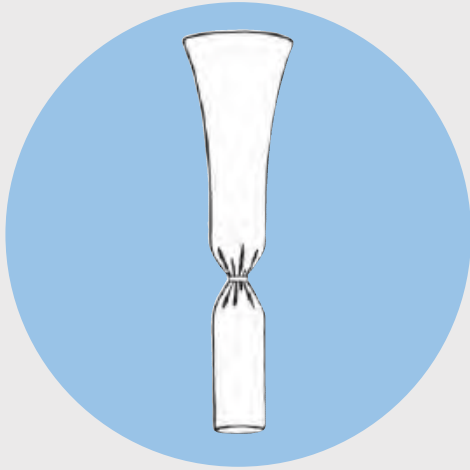

TAPERED  
COMPRESSOGRIP

PROVIDES ADDITIONAL WIDTH AT THE TOP TO ACCOMMODATE FULLER THIGHS AND REDUCE ROLLING AND BUNCHING BEHIND THE KNEE CAUSED BY TOO MUCH COMPRESSION ABOVE THE KNEE.

- ✓ Made from latex-free materials
- ✓ Washable
- ✓ Patent pending

## MATERIALS

White, Beige — 92% polyester, 8% spandex

Polyester/Ionic™ Grey — 87.4% polyester, 4.6% Ionic™, 8% spandex

## FITTING INSTRUCTIONS



1. Apply the first layer.



2. Position the ring snugly.



3. Apply the second layer to a point of desired maximum compression.

## TAPERED COMPRESSOGRIP® B/K ELASTIC SHRINKER

### MEASUREMENT OF B/K CIRCUMFERENCE 2" BELOW MEDIAL TIBIA TUBERCLE

| WHITE      | BEIGE      | GREY/<br>IONIC™ | SIZE | LENGTH     | FITS B/K CIRC.    | FLAT WIDTH  | TOP  |
|------------|------------|-----------------|------|------------|-------------------|-------------|--|
| 2H05-T4-24 | 2H0B-T4-24 | 2H-0XT4-24      | #4   | 24" (60cm) | 12"-16" (30-40cm) | 4" (10.2cm) | 5" Flat Width; Tapers Last 4" of Length    |
| 2H05-T4-30 | 2H0B-T4-30 | 2H-0XT4-30      | #4   | 30" (75cm) | 12"-16" (30-40cm) | 4" (10.2cm) | 5" Flat Width; Tapers Last 4" of Length    |
| —          | 2H0B-T5-24 | 2H-0XT5-24      | #5   | 24" (60cm) | 15"-20" (38-50cm) | 5" (12.7cm) | 6.50" Flat Width; Tapers Last 5" of Length |
| 2H05-T5-30 | 2H0B-T5-30 | 2H-0XT5-30      | #5   | 30" (75cm) | 15"-20" (38-50cm) | 5" (12.7cm) | 6.50" Flat Width; Tapers Last 5" of Length |
| —          | —          | 2H-0XT6-24      | #6   | 24" (60cm) | 18"-24" (45-60cm) | 6" (15.2cm) | 7.50" Flat Width; Tapers Last 6" of Length |
| —          | —          | 2H-0XT6-30      | #6   | 30" (75cm) | 18"-24" (45-60cm) | 6" (15.2cm) | 7.50" Flat Width; Tapers Last 6" of Length |

When determining size, allow for swelling.

## BONE RINGS

### BONE RINGS

| ITEM #     | SIZE   |
|------------|--------|
| 5B70-55-07 | 7/8"   |
| 5B70-55-10 | 1"     |
| 5B70-55-11 | 1-1/8" |
| 5B70-55-15 | 1-1/2" |
| 5B70-55-06 | 3/4"   |



See page 39 for  
Shrinker Donning Tubes.

# COMPRESSOGRIP A/K PROSTHETIC SHRINKER



Compressogrip  
A/K Shrinker  
U.S. Patent No. 7,575,561

- ✓ Clinically effective and easy to use and apply
- ✓ Lower cost than traditional shrinkers — fewer sizes fit more patients
- ✓ Easy and affordable to keep on hand for same day delivery
- ✓ The fabric used for these shrinkers is circular knit with increased power in the elastic, producing a heavier fabric than that used for B/K Compressogrip shrinkers
- ✓ Multi-directional stretch for improved fitting qualities and a universal waist belt
- ✓ Available in standard white Compressogrip and Ionic+™ styles
- ✓ Gusset is a non-compressive 1-ply fabric at medial aspect of transfemoral limb close to the groin
- ✓ 25-30mmHg Compression / 50% Stretch
- ✓ Garment features a 2" waist belt

## COMPRESSOGRIP® A/K SHRINKER

### MEASUREMENT OF A/K CIRCUMFERENCE 2" FROM DISTAL END

| WHITE                     | IONIC+™     | SIZE    | FITS CIRCUMFERENCE | FLAT WIDTH |
|---------------------------|-------------|---------|--------------------|------------|
| 2H-0S05-XS                | 2H-0X05-XS  | X-Small | 12"-17" (30-43cm)  | 5"         |
| 2H-0S06-SM                | 2H-0X06-SM  | Small   | 15"-20" (38-50cm)  | 6"         |
| 2H-0SS7-MD                | 2H-0XS7-MD  | Medium  | 18"-24" (45-60cm)  | 7.50"      |
| 2H-0S09-LG*               | 2H-0X09-LG* | Large   | 24"-32" (60-80cm)  | 9"         |
| 2H-GS05-XS (WITH GUSSET)  | —           | X-Small | 12"-17" (30-43cm)  | 5"         |
| 2H-GS06-SM (WITH GUSSET)  | —           | Small   | 15"-20" (38-50cm)  | 6"         |
| 2H-GSS7-MD (WITH GUSSET)  | —           | Medium  | 18"-24" (45-60cm)  | 7.50"      |
| 2H-GS09-LG* (WITH GUSSET) | —           | Large   | 24"-32" (60-80cm)  | 9"         |
| 2H-GS09-L2 (WITH GUSSET)  | —           | Large   | 24"-32" (60-80cm)  | 9"         |
| 2H-GS10-XL (WITH GUSSET)* | —           | X-Large | 30"-40" (76-102cm) | 10.5"      |

\*Garment features a 4" waist belt



For custom logo shrinkers, see page 35, or ask customer service or your distributor for program details.

# POST-OP DRESSING RETENTION SOCK



POST-OP DRESSING  
RETENTION SOCK

UTILIZES A VERY LIGHT COMPRESSION FIT, WITH A BELT SUSPENSION TO COVER AND RETAIN POSTOPERATIVE DRESSINGS FOR TRANS-FEMORAL AMPUTEES.

Surgeons often want to limit compression therapy, while still needing to retain post-operative dressings. In the past, many methods have been used with limited success. The Post-Op Dressing Retention Sock solves this problem.

- ✓ Very soft LYCRA® and COOLMAX® body for very low compression
- ✓ Universal right or left belt system to maintain secure suspension
- ✓ Aids in retaining surgical dressing in place—particularly with obese patients
- ✓ Can be used over dressing and then use an ace bandage over the top to apply “controlled” compression
- ✓ Does not irritate or stress sutures during donning
- ✓ Trains patients in sock management
- ✓ Can be used in combination with removable rigid dressings
- ✓ Helps assess patient compliance with future pre-prosthetic regimens
- ✓ Low cost
- ✓ Little staff training required at hospital
- ✓ Gets patient into the prosthetic services earlier
- ✓ Does not compromise peripheral circulation in marginal bypass vascular grafts; at times limited compression is desired for the first week or so

## POST-OP DRESSING RETENTION SOCK

| ITEM #     | SIZE            | FITS CIRCUMFERENCE | LENGTH |
|------------|-----------------|--------------------|--------|
| 1R-P104-28 | 4" Post Op Sock | 10"-12" (25-30cm)  | 28"    |
| 1R-P106-32 | 6" Post Op Sock | 14"-18" (35-45cm)  | 32"    |
| 1R-P108-32 | 8" Post Op Sock | 18"-24" (40-60cm)  | 32"    |

For custom logo shrinkers, see page 35, or ask customer service or your distributor for program details.

## 4-WAY STRETCH SHRINKER



- ✓ Features rounded toe with stretchy polyester/LYCRA® Core-Spun yarns and dot silicone suspension band
- ✓ Available in Medium (20–30 mmHg) and Heavy (30–40 mmHg) compression, Transtibial and Transfemoral sizes
- ✓ Multidirectional stretch provides improved fit qualities and comfort
- ✓ Softer, more comfortable and easier to don
- ✓ Moisture-wicking properties
- ✓ A/K version also available with waist belt
- ✓ Black shrinker includes Ionic+™ Mineral Antimicrobial fibers to help prevent odor in the garment, provide thermal properties, and is anti-static

### MATERIALS

#### 20–30mmHg

**White** — 82% polyester, 18% spandex


**Black** — 72% polyester, 4% Ionic+™, 18% spandex

#### 30–40 mmHg

**White** — 76% polyester, 24% spandex

**Black** — 72% polyester, 4% Ionic+™, 24% spandex

### ORDERING GUIDELINES

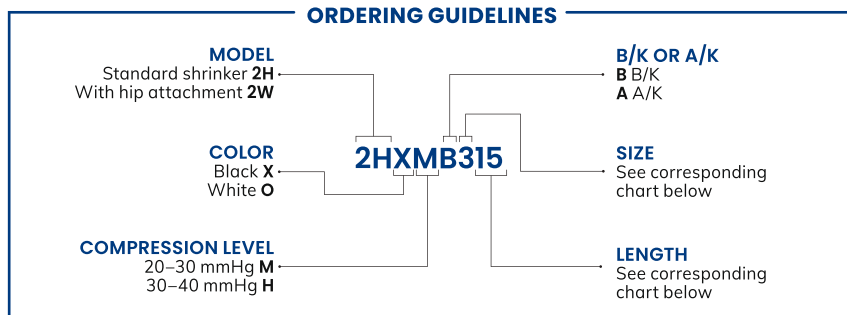

### B/K SHRINKER SIZE CHART

1. Measure the circumference of areas F, E, and C to find the correct size.

2. Measure from the end of the limb up to the F area to find the correct length.

|                             | F                                | E                               | C                                |
|-----------------------------|----------------------------------|---------------------------------|----------------------------------|
| SIZE                        | CIRCUMFERENCES                   |                                 |                                  |
| 0                           | 14" – 15.75"<br>(35.5 – 40cm)    | 10.75" – 11.75"<br>(27 – 30cm)  | 10" – 11"<br>(25.5 – 28cm)       |
| 1                           | 15.75" – 16.75"<br>(40 – 42.5cm) | 11.75" – 12.75"<br>(30 – 32cm)  | 11" – 12"<br>(28 – 30cm)         |
| 2                           | 16.75" – 17.75"<br>(42.5 – 45cm) | 12.75" – 14"<br>(32 – 35.5cm)   | 12" – 13"<br>(30 – 33cm)         |
| 3                           | 17.75" – 19"<br>(45 – 48 cm)     | 14" – 15"<br>(35.5 – 38cm)      | 13" – 14.5"<br>(33 – 37cm)       |
| 4                           | 19" – 20.25"<br>(48 – 51.5cm)    | 15" – 16.25"<br>(38 – 41cm)     | 14.5" – 15.5"<br>(37 – 39cm)     |
| 5                           | 20.25" – 21.25"<br>(51.5 – 54cm) | 16.25" – 17.5"<br>(41 – 44.5cm) | 15.5" – 16.75"<br>(39 – 42.5cm)  |
| 6                           | 21.25" – 22.75"<br>(54 – 58cm)   | 17.5" – 18.5"<br>(44.5 – 47cm)  | 16.75" – 17.75"<br>(42.5 – 45cm) |
| BELOW KNEE SHRINKER LENGTHS |                                  |                                 |                                  |
| SHORT                       | MEDIUM                           | LONG                            |                                  |
| 12" (30cm)                  | 15" (38cm)                       | 18" (46cm)                      |                                  |

\*Extra Long is not available in size 0.

### AK SHRINKER SIZE CHART

1. Measure the circumference of areas G, F, and E to find the correct size.

2. Measure from the end of the limb up to the F area to find the correct length.

|                             | G                                | F                                | E                              |
|-----------------------------|----------------------------------|----------------------------------|--------------------------------|
| SIZE                        | CIRCUMFERENCES                   |                                  |                                |
| 0                           | 15.5" – 16.75"<br>(39 – 42.5cm)  | 14.75" – 15.75"<br>(37.5 – 40cm) | 10.75" – 11.75"<br>(27 – 30cm) |
| 1                           | 16.75" – 18"<br>(42.5 – 46cm)    | 15.75" – 16.75"<br>(40 – 42.5cm) | 11.75" – 12.75"<br>(30 – 32cm) |
| 2                           | 18" – 19.75"<br>(46 – 50cm)      | 16.75" – 17.75"<br>(42.5 – 45cm) | 12.75" – 14"<br>(32 – 35.5cm)  |
| 3                           | 19.75" – 21.25"<br>(50 – 54 cm)  | 17.75" – 19"<br>(45 – 48cm)      | 14" – 15"<br>(35.5 – 38cm)     |
| 4                           | 21.25" – 23"<br>(54 – 58.5cm)    | 19" – 20.25"<br>(48 – 51.5cm)    | 15" – 16.5"<br>(38 – 42cm)     |
| 5                           | 23" – 24.5"<br>(58.5 – 62cm)     | 20.25" – 21.25"<br>(51.5 – 54cm) | 16.5" – 17.5"<br>(42 – 44.5cm) |
| 6                           | 24.25" – 26.25"<br>(62 – 66.5cm) | 21.25" – 22.75"<br>(54 – 58cm)   | 17.5" – 18.5"<br>(44.5 – 47cm) |
| ABOVE KNEE SHRINKER LENGTHS |                                  |                                  |                                |
| SHORT                       | MEDIUM                           | LONG                             | X-LONG*                        |
| 8" (20cm)                   | 10" (25cm)                       | 12" (30cm)                       | 14" (35cm)                     |



# 4-WAY STRETCH SHRINKER

See next page for more  
4-Way Stretch Shrinkers.

## B/K 4-WAY STRETCH SHRINKER

| WHITE<br>POLYESTER/<br>LYCRA® | BLACK IONIC+™<br>POLYESTER/<br>IONIC+™/LYCRA® | SIZE | FITS<br>DISTAL LENGTH | LENGTH |
|-------------------------------|---|------|-----------------------|--------|
| 20–30 mmHg                    |   |      |                       |        |
| 2H-OMB0-12                    | 2H-XMB0-12                                    | 0    | 12"; 30.5 cm          | Short  |
| 2H-OMB0-15                    | 2H-XMB0-15                                    | 0    | 15"; 38.1cm           | Medium |
| 2H-OMB1-12                    | 2H-XMB1-12                                    | I    | 12"; 30.5 cm          | Short  |
| 2H-OMB1-15                    | 2H-XMB1-15                                    | I    | 15"; 38.1 cm          | Medium |
| 2H-OMB2-12                    | 2H-XMB2-12                                    | II   | 12"; 30.5 cm          | Short  |
| 2H-OMB2-15                    | 2H-XMB2-15                                    | II   | 15"; 38.1 cm          | Medium |
| 2H-OMB3-12                    | 2H-XMB3-12                                    | III  | 12"; 30.5 cm          | Short  |
| 2H-OMB3-15                    | 2H-XMB3-15                                    | III  | 15"; 38.1 cm          | Medium |
| 2H-OMB4-12                    | 2H-XMB4-12                                    | IV   | 12"; 30.5 cm          | Short  |
| 2H-OMB4-15                    | 2H-XMB4-15                                    | IV   | 15"; 38.1 cm          | Medium |
| 2H-OMB5-12                    | 2H-XMB5-12                                    | V    | 12"; 30.5 cm          | Short  |
| 2H-OMB5-15                    | 2H-XMB5-15                                    | V    | 15"; 38.1 cm          | Medium |
| 2H-OMB6-12                    | 2H-XMB6-12                                    | VI   | 12"; 30.5 cm          | Short  |
| 2H-OMB6-15                    | 2H-XMB6-15                                    | VI   | 15"; 38.1 cm          | Medium |
| 30–40 mmHg                    |   |      |                       |        |
| 2H-OHB0-12                    | 2H-XHB0-12                                    | 0    | 12"; 30.5 cm          | Short  |
| 2H-OHB0-15                    | 2H-XHB0-15                                    | 0    | 15"; 38.1 cm          | Medium |
| —                             | 2H-XHB0-18                                    | 0    | 18"; 45.7cm           | Long   |
| 2H-OHB1-12                    | 2H-XHB1-12                                    | I    | 12"; 30.5 cm          | Short  |
| 2H-OHB1-15                    | 2H-XHB1-15                                    | I    | 15"; 38.1 cm          | Medium |
| —                             | 2H-XHB1-18                                    | I    | 18"; 45.7cm           | Long   |
| 2H-OHB2-12                    | 2H-XHB2-12                                    | II   | 12"; 30.5 cm          | Short  |
| 2H-OHB2-15                    | 2H-XHB2-15                                    | II   | 15"; 38.1 cm          | Medium |
| —                             | 2H-XHB2-18                                    | II   | 18"; 45.7cm           | Long   |
| 2H-OHB3-12                    | 2H-XHB3-12                                    | III  | 12"; 30.5 cm          | Short  |
| 2H-OHB3-15                    | 2H-XHB3-15                                    | III  | 15"; 38.1 cm          | Medium |
| —                             | 2H-XHB3-18                                    | III  | 18"; 45.7cm           | Long   |
| 2H-OHB4-12                    | 2H-XHB4-12                                    | IV   | 12"; 30.5 cm          | Short  |
| 2H-OHB4-15                    | 2H-XHB4-15                                    | IV   | 15"; 38.1 cm          | Medium |
| —                             | 2H-XHB4-18                                    | IV   | 18"; 45.7cm           | Long   |
| 2H-OHB5-12                    | 2H-XHB5-12                                    | V    | 12"; 30.5 cm          | Short  |
| 2H-OHB5-15                    | 2H-XHB5-15                                    | V    | 15"; 38.1 cm          | Medium |
| —                             | 2H-XHB5-18                                    | V    | 18"; 45.7cm           | Long   |
| 2H-OHB6-12                    | 2H-XHB6-12                                    | VI   | 12"; 30.5 cm          | Short  |
| 2H-OHB6-15                    | 2H-XHB6-15                                    | VI   | 15"; 38.1 cm          | Medium |
| —                             | 2H-XHB6-18                                    | VI   | 18"; 45.7cm           | Long   |

## A/K 4-WAY STRETCH SHRINKER W/ HIP ATTACHMENT\*

| WHITE<br>POLYESTER/<br>LYCRA® | BLACK IONIC+™<br>POLYESTER/<br>IONIC+™/LYCRA® | SIZE | FITS<br>DISTAL LENGTH | LENGTH |
|-------------------------------|---|------|-----------------------|--------|
| 20–30 mmHg                    |   |      |                       |        |
| 2W-OMA3-10                    | 2W-XMA3-10                                    | III  | 10"; 25 cm            | Medium |
| 2W-OMA3-12                    | 2W-XMA3-12                                    | III  | 12"; 30 cm            | Long   |
| 2W-OMA4-10                    | 2W-XMA4-10                                    | IV   | 10"; 25 cm            | Medium |
| 2W-OMA4-12                    | 2W-XMA4-12                                    | IV   | 12"; 30 cm            | Long   |
| 2W-OMA5-10                    | 2W-XMA5-10                                    | V    | 10"; 25 cm            | Medium |
| 2W-OMA5-12                    | 2W-XMA5-12                                    | V    | 12"; 30 cm            | Long   |
| 2W-OMA6-10                    | 2W-XMA6-10                                    | VI   | 10"; 25 cm            | Medium |
| —                             | 2W-XMA6-12                                    | VI   | 12"; 30 cm            | Long   |
| —                             | 2W-XMA6-14                                    | VI   | 14"; 35 cm            | X-Long |

\*A/K version with  
waist belt is not  
available in size 0.



# 4-WAY STRETCH SHRINKER

## A/K 4-WAY STRETCH SHRINKER


| WHITE<br>POLYESTER/<br>LYCRA® | BLACK IONIC+™<br>POLYESTER/<br>IONIC+™/LYCRA® | SIZE | FITS<br>DISTAL LENGTH | LENGTH |
|-------------------------------|---|------|-----------------------|--------|
| 20–30mmHg                     |   |      |                       |        |
| 2H-0MA0-08                    | 2HXMA008                                      | 0    | 8"; 20 cm             | Short  |
| 2H-0MA0-10                    | 2HXMA010                                      | 0    | 10"; 25 cm            | Medium |
| 2H-0MA0-12                    | 2HXMA012                                      | 0    | 12"; 30.5 cm          | Long   |
| 2H-0MA1-08                    | 2H-XMA1-08                                    | I    | 8"; 20 cm             | Short  |
| 2H-0MA1-10                    | 2H-XMA1-10                                    | I    | 10"; 25 cm            | Medium |
| 2H-0MA1-12                    | 2H-XMA1-12                                    | I    | 12"; 30 cm            | Long   |
| 2H-0MA1-14                    | 2H-XMA1-14                                    | I    | 14"; 35 cm            | X-Long |
| 2H-0MA2-08                    | 2H-XMA2-08                                    | II   | 8"; 20 cm             | Short  |
| 2H-0MA2-10                    | 2H-XMA2-10                                    | II   | 10"; 25 cm            | Medium |
| 2H-0MA2-12                    | 2H-XMA2-12                                    | II   | 12"; 30 cm            | Long   |
| 2H-0MA2-14                    | 2H-XMA2-14                                    | II   | 14"; 35 cm            | X-Long |
| 2H-0MA3-08                    | 2H-XMA3-08                                    | III  | 8"; 20 cm             | Short  |
| 2H-0MA3-10                    | 2H-XMA3-10                                    | III  | 10"; 25 cm            | Medium |
| 2H-0MA3-12                    | 2H-XMA3-12                                    | III  | 12"; 30 cm            | Long   |
| 2H-0MA3-14                    | 2H-XMA3-14                                    | III  | 14"; 35 cm            | X-Long |
| 2H-0MA4-08                    | 2H-XMA4-08                                    | IV   | 8"; 20 cm             | Short  |
| 2H-0MA4-10                    | 2H-XMA4-10                                    | IV   | 10"; 25 cm            | Medium |
| 2H-0MA4-12                    | 2H-XMA4-12                                    | IV   | 12"; 30 cm            | Long   |
| 2H-0MA4-14                    | 2H-XMA4-14                                    | IV   | 14"; 35 cm            | X-Long |
| 2H-0MA5-08                    | 2H-XMA5-08                                    | V    | 8"; 20 cm             | Short  |
| 2H-0MA5-10                    | 2H-XMA5-10                                    | V    | 10"; 25 cm            | Medium |
| 2H-0MA5-12                    | 2H-XMA5-12                                    | V    | 12"; 30 cm            | Long   |
| 2H-0MA5-14                    | 2H-XMA5-14                                    | V    | 14"; 35 cm            | X-Long |
| 2H-0MA6-08                    | 2H-XMA6-08                                    | VI   | 8"; 20 cm             | Short  |
| 2H-0MA6-10                    | 2H-XMA6-10                                    | VI   | 10"; 25 cm            | Medium |
| 2H-0MA6-12                    | 2H-XMA6-12                                    | VI   | 12"; 30 cm            | Long   |
| 2H-0MA6-14                    | 2H-XMA6-14                                    | VI   | 14"; 35 cm            | X-Long |
| 30–40mmHg*                    |   |      |                       |        |
| 2H-0HA1-08                    | 2H-XHA1-08                                    | I    | 8"; 20 cm             | Short  |
| 2H-0HA1-10                    | 2H-XHA1-10                                    | I    | 10"; 25 cm            | Medium |
| 2H-0HA1-12                    | 2H-XHA1-12                                    | I    | 12"; 30 cm            | Long   |
| 2H-0HA1-14                    | 2H-XHA1-14                                    | I    | 14"; 35 cm            | X-Long |
| 2H-0HA2-08                    | 2H-XHA2-08                                    | II   | 8"; 20 cm             | Short  |
| 2H-0HA2-10                    | 2H-XHA2-10                                    | II   | 10"; 25 cm            | Medium |
| 2H-0HA2-12                    | 2H-XHA2-12                                    | II   | 12"; 30 cm            | Long   |
| 2H-0HA2-14                    | 2H-XHA2-14                                    | II   | 14"; 35 cm            | X-Long |
| 2H-0HA3-08                    | 2H-XHA3-08                                    | III  | 8"; 20 cm             | Short  |
| 2H-0HA3-10                    | 2H-XHA3-10                                    | III  | 10"; 25 cm            | Medium |
| 2H-0HA3-12                    | 2H-XHA3-12                                    | III  | 12"; 30 cm            | Long   |
| 2H-0HA3-14                    | 2H-XHA3-14                                    | III  | 14"; 35 cm            | X-Long |
| 2H-0HA4-08                    | 2H-XHA4-08                                    | IV   | 8"; 20 cm             | Short  |
| 2H-0HA4-10                    | 2H-XHA4-10                                    | IV   | 10"; 25 cm            | Medium |
| 2H-0HA4-12                    | 2H-XHA4-12                                    | IV   | 12"; 30 cm            | Long   |
| 2H-0HA4-14                    | 2H-XHA4-14                                    | IV   | 14"; 35 cm            | X-Long |
| 2H-0HA5-08                    | 2H-XHA5-08                                    | V    | 8"; 20 cm             | Short  |
| 2H-0HA5-10                    | 2H-XHA5-10                                    | V    | 10"; 25 cm            | Medium |
| 2H-0HA5-12                    | 2H-XHA5-12                                    | V    | 12"; 30 cm            | Long   |
| 2H-0HA5-14                    | 2H-XHA5-14                                    | V    | 14"; 35 cm            | X-Long |
| 2H-0HA6-08                    | 2H-XHA6-08                                    | VI   | 8"; 20 cm             | Short  |
| 2H-0HA6-10                    | 2H-XHA6-10                                    | VI   | 10"; 25 cm            | Medium |
| 2H-0HA6-12                    | 2H-XHA6-12                                    | VI   | 12"; 30 cm            | Long   |
| 2H-0HA6-14                    | 2H-XHA6-14                                    | VI   | 14"; 35 cm            | X-Long |

\*30–40mmHg A/K  
is not available  
in size 0.

# CUSTOM PROSTHETIC SHRINKERS



When ordering the shrinker, state any additional circumference measurements that would be necessary to fit the special need and the distance of those measurements from the distal end.

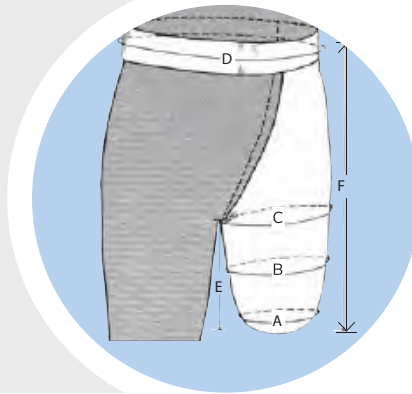

## ACCOMMODATE IRREGULARITIES AND SPECIAL FITTING PROBLEMS

- ✓ Medium compression approximately 30–35mmHg pressure at 50% stretch
- ✓ Aids in the reduction of edema and the shaping of the limb
- ✓ Unique construction knit with grey polyester/Ionic+™/LYCRA® core-spun yarn
- ✓ Soft, comfortable and easy to put on
- ✓ Minimum lead time 24-48 hours

**CAUTION: COMPRESSION PROXIMALLY SHOULD NOT BE GREATER THAN COMPRESSION DISTALLY.**

## MADE TO ORDER A/K SHRINKER

| ORDER FORM — SEND THE FOLLOWING MEASUREMENTS: |   |
|---|---|
| A   | DENOTE MEASUREMENTS IN                  |
|   | IN. CM.                                 |
| B   | CIRCUMFERENCE 1" TO DISTAL END (2.5 CM) |
| C   | CIRCUMFERENCE 3" TO DISTAL END (7.6 CM) |
| D   | CIRCUMFERENCE AT GROIN                  |
| E   | CIRCUMFERENCE AT WAIST                  |
| F   | LENGTH FROM GROIN TO DISTAL END         |
| G.  | LENGTH FROM WAIST TO DISTAL END         |
|   | LEFT RIGHT                              |
| H.  | QUANTITY                                |
| I.  | MEDIUM COMPRESSION                      |

PATIENT I.D.

NOTE: When ordering, please have the above information ready along with your customer number, P.O. number, and any additional instruction you deem important.

## MADE-TO-ORDER SHRINKER

| ORDER FORM — SEND THE FOLLOWING MEASUREMENTS: |   |
|---|---|
| A   | DENOTE MEASUREMENTS IN  |
|   | IN. CM.   |
| B   | CIRCUMFERENCE AT DISTAL END                                     |
| C   | CIRCUMFERENCE AT ANKLE AT AXILLA OR 3" (7.6 IN) FROM DISTAL END |
| D   | CIRCUMFERENCE AT TOP OR SHOULDER WIDTH DESIRED                  |
| E   | LENGTH OF SOCK  |
|   | 1) LENGTH FROM DISTAL TO AXILLA                                 |
|   | 2) LENGTH FROM AXILLA TO SHOULDER                               |
|   | 3) CHEST CIRCUMFERENCE AT ARMPIT                                |
|   | LEFT RIGHT  |
|   | QUANTITY  |

MEDIUM COMPRESSION


PATIENT I.D.

## CARE TIPS


May be machine washed, dried, warm temperature. No bleach.




Unique Construction: LYCRA® /Spandex is in the core of the yarn and in each knit loop, locked in!

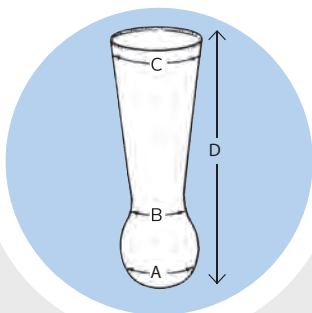

Symes

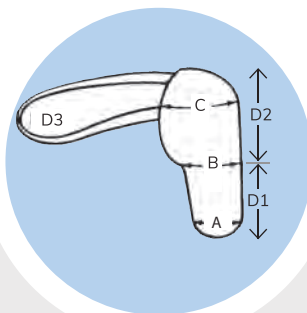

Arm

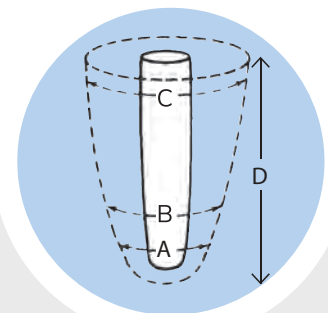

Smaller or larger than standard sizes

# SHRINKER PRIVATE LOGO PROGRAM



**AVAILABLE IN ALL B/K AND A/K VERSIONS OF COMPRESSOGRIP® AND 4-WAY SHRINKERS!**

- ✓ **One-time setup fee: 1 Color for      or 2 Colors for**
- ✓ NO additional charge per SKU or fee after initial one time setup.
- ✓ Lead times: Initial lead time is 4 weeks on receipt of final signed artwork approval. Reorders are approximately 5-7 working days. Working days do not include shipping time.
- ✓ **Graphic Requirements: Art must be provided in editable vector format (AI or EPS files only) with all fonts converted to outlines. JPEGs and PDFs will not be accepted. Knit-Rite can re-draw in the correct format for a one time fee of      if not available in the required format.**
- ✓ Maximum graphic size: 2.5"W x 1.5"H. Knit-Rite will scale art for approval.
- ✓ Consider imprint color(s) that will contrast with the fabric. Grey works well if ordering more than one fabric color.
- ✓ Order minimum: ≥ 20 total units. No minimum per SKU.
- ✓ Completed form via e-mail or fax required. Transfers will be ordered after receiving an order, applicable setup fee(s) invoiced, and art approved.



PROSTHETIC SHRINKERS



## ORDER FORM

|   |              |                              |  |                                |
|---|--------------|------------------------------|--|--------------------------------|
| Bill To Address   |              | Ship To Address              |  |                                |
| Knit-Rite Acct. #   | Order Date   | P.O.#                        | Customer Contact Information   |                                |
| Is order through a Knit-Rite Distributor?<br>Yes <input type="checkbox"/> No <input type="checkbox"/> |              | If Yes, which Distributor?   |  | Practice Distributor Account # |
| Logo Color(s)* Pantone<br>Colors preferred  | Logo Name    | Desired Graphic<br>Dimension | Initial Order:<br>Yes <input type="checkbox"/> No <input type="checkbox"/> |                                |
| <b>KR PRODUCT CODE</b>  | <b>STYLE</b> | <b>PRODUCT DESCRIPTION</b>   | <b>SIZE</b>  | <b>TOTAL UNITS</b>             |
|   |              |                              |  |                                |
|   |              |                              |  |                                |
|   |              |                              |  |                                |
|   |              |                              |  |                                |
|   |              |                              |  |                                |
|   |              |                              |  |                                |
|   |              |                              |  |                                |

\*Consider imprint color(s) that will contrast the fabric. Grey works well if ordering more than one fabric color. Please reference this catalog for product codes and descriptions as some products may not be available in all materials or fabric colors.

For more information please contact Knit-Rite by phone at **800-821-3094**, by fax at **800-462-4707**, or by e-mail **cs@knitrite.com**. You can also contact your Knit-Rite distributor.


# PROSTHETIC SUSPENSION *and* DONNING

## PROSTHETIC SUSPENSION AND DONNING | 36–41

- 37** Power Belt® Auxiliary Suspension Belt
- 38** Cool-Flex™ Auxiliary Suspension Belt
- 38** B/K Suspension Belts
- 39** Power Pull® Donning Sock
- 39** Shrinker Donning Tube
- 40** Sleeve Saver™
- 40** Foot Shell Sock
- 41** Cosmetic Stockings




# POWER BELT® AUXILIARY SUSPENSION BELT

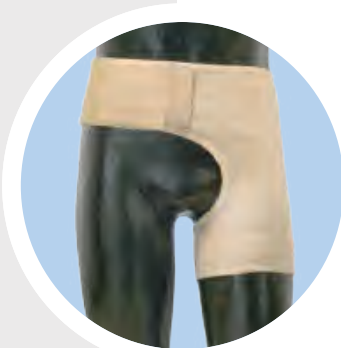

ORIGINAL POWER BELT®  
SUSPENSION BELT  
WITH HAND LOOP



POWER BELT®  
SUSPENSION BELT WITH  
GUNSLINGER WAIST

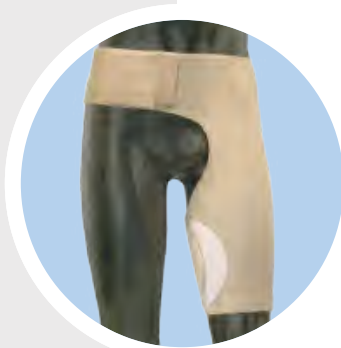

POWER BELT®  
SUSPENSION BELT  
WITH ADJUSTABLE  
WAIST AND LEG

- ✓ High acceptance, greater functional benefit
- ✓ Comfort and confidence for patient
- ✓ Patented dynamic material application is uniquely low profile, and comfortable
- ✓ Designed to allow wear in a more natural transverse (versus oblique) orientation
- ✓ Breathable material is cool and comfortable
- ✓ Latex free materials: 68% nylon, 25% spandex, 3% polyester, 2% silicone, 1% ester, 1% rubber (rubber in leg elastic only)
- ✓ For daily use, at least two belts of each size and style selected should be provided to permit normal wearer hygiene

## ORIGINAL POWER BELT® WITH HAND LOOP

| SIZE     | LEFT       | RIGHT      | FITS WAIST | FITS LEG |
|----------|------------|------------|------------|----------|
| Small    | 3L-928L-SM | 3L-928R-SM | 23"-30"    | 15"-20"+ |
| Medium   | 3L-928L-MD | 3L-928R-MD | 30"-42"    | 17"-23"+ |
| Large    | 3L-928L-LG | 3L-928R-LG | 34"-46"    | 19"-25"+ |
| X-Large  | 3L-928L-XL | 3L-928R-XL | 38"-52"    | 20"-28"+ |
| XX-Large | 3L-928L-2X | 3L-928R-2X | 42"-56"    | 22"-30"+ |

## POWER BELT® WITH GUNSLINGER WAIST

| SIZE    | UNIVERSAL  | FITS WAIST | FITS LEG |
|---------|------------|------------|----------|
| Medium  | 3L-977U-MD | 30"-42"    | 17"-23"+ |
| Large   | 3L-977U-LG | 34"-46"    | 19"-25"+ |
| X-Large | 3L-977U-XL | 38"-52"    | 20"-28"+ |

## POWER BELT® WITH ADJUSTABLE WAIST AND LEG

| SIZE    | LEFT       | RIGHT      | FITS WAIST | FITS LEG |
|---------|------------|------------|------------|----------|
| X-Small | —          | 3L-991R-XS | 21"-28"    | 12"-15"  |
| Small   | 3L-991L-SM | 3L-991R-SM | 23"-30"    | 15"-20"+ |
| Medium  | 3L-991L-MD | 3L-991R-MD | 30"-42"    | 17"-23"+ |
| Large   | 3L-991L-LG | 3L-991R-LG | 34"-46"    | 19"-25"+ |
| X-Large | 3L-991L-XL | 3L-991R-XL | 38"-52"    | 20"-28"+ |

## BELT EXTENDER

| ITEM #     | DESCRIPTION             | SIZE                    |
|------------|-------------------------|-------------------------|
| 3L-914E-12 | Adds up to 12" to waist | For sizes medium and up |



Power Belt®  
Belt Extender

Patent #6,652,596  
Manufacturer's suggested  
L-Code: L5695



# COOL-FLEX™ AUXILIARY SUSPENSION BELT


## COOL-FLEX™ AUXILIARY SUSPENSION BELT

| ITEM #     | SIZE   | PROXIMAL BRIM      | COLOR |
|------------|--------|--------------------|-------|
| 3L-CFAK-SM | Small  | 16"-20" (40-50 cm) | Black |
| 3L-CFAK-MD | Medium | 20"-24" (50-60 cm) | Black |
| 3L-CFAK-LG | Large  | 24"-30" (60-75 cm) | Black |

- ✓ The comfortable alternative to Neoprene
- ✓ "Gun Slinger" Hip Hugger design improves stability, comfort, and ease of application
- ✓ COOL-FLEX™ material is breathable, moisture-wicking, and lets air in to keep the wearer cool and comfortable
- ✓ Universal design for right or left application
- ✓ Lightweight, durable, and conformable
- ✓ Non-Allergenic



## B/K SUSPENSION BELTS

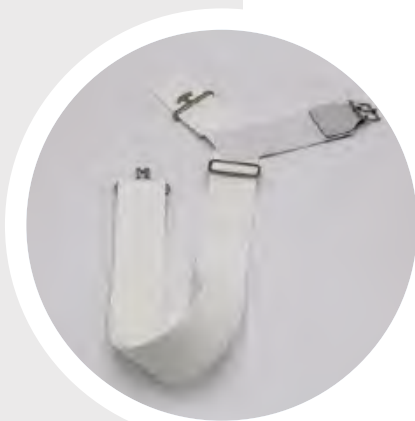

ADJUSTABLE B/K BELT WITH QUICK RELEASE CLASP FASTENER



B/K SUSPENSION BELT

### ADJUSTABLE B/K BELT WITH QUICK RELEASE CLASP FASTENER

| ITEM #     | SIDE  | MEASUREMENTS               |
|------------|-------|----------------------------|
| 3L-BK2-0LF | Left  | Adjustable from 29" to 46" |
| 3L-BK2-0RT | Right | Adjustable from 29" to 46" |


- ✓ Features a full range of length adjustment from 29" to 46"
- ✓ Reduces the necessity to stock a wide assortment of sizes
- ✓ 2" prosthetic webbing with 2" elastic anterior support strap
- ✓ 1" strap buckle

### B/K SUSPENSION BELT

| ITEM #     | SIDE  | MEASUREMENTS  |
|------------|-------|---------------|
| 3M-2KSB-LF | Left  | Up to 44" hip |
| 3M-2KSB-RT | Right | Up to 44" hip |

- ✓ For immediate post-operative prosthetic fittings
- ✓ Lightweight webbing waist belt with elastic suspension strap and cast attachment strap with buckle
- ✓ One size

# POWER PULL® DONNING SOCK

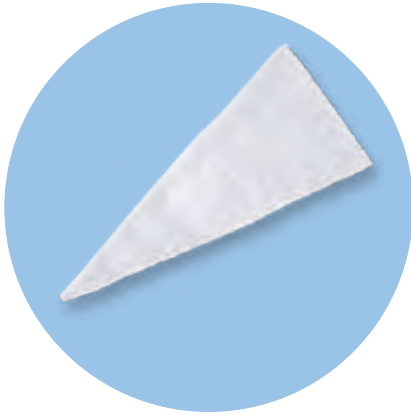

POWER PULL®  
DONNING SOCK

- ✓ Made of a slick, supple fabric for the easiest donning ever
- ✓ Single piece integral design provides a more direct and effective pull transfer, and eliminates stress points for a more durable product
- ✓ 100% Nylon.

## POWER PULL® DONNING SOCK

| ITEM #     | SIZE    | TOP | PROXIMAL CIRCUM. |
|------------|---------|-----|------------------|
| 3M-PP10-SM | Small   | 10" | 16"-20"          |
| 3M-PP12-MD | Medium  | 12" | 20"-24"          |
| 3M-PP14-LG | Large   | 14" | 24"-28"          |
| 3M-PP16-XL | X-Large | 16" | 28"-32"          |

Suggested L-code: L7600



# SHRINKER DONNING TUBE

THE DONNING TUBE IS THE ESSENTIAL TOOL FOR APPLICATION OF SHRINKERS ON THE SENSITIVE LIMB.

- ✓ Great for physical and occupational therapist for training and reapplication of Compressogrip shrinkers
- ✓ By loading the shrinker onto the tube and then sliding onto the limb, one person can apply compression pain free
- ✓ Acrylic material



8" and 6" Donning Tubes

## SHRINKER DONNING TUBE

| ITEM #     | LENGTH |
|------------|--------|
| 6S-TUBE-06 | 6"     |
| 6S-TUBE-08 | 8"     |

## CARE TIPS

Clean with a mild detergent and wipe with a clean, soft cloth. Do not use glass cleaner or paper towels.

## DONNING TUBE INSTRUCTIONS



# SLEEVE SAVER™

## SLEEVE SAVER™

| ITEM #     | STYLE  | LENGTH |
|------------|--|--------|
| 3M-SPB0-10 | Regular                                      | 10"    |
| 3M-SPB0-NA | Narrow                                       | 10"    |
| 3M-SPBE-12 | Regular w/ Elasticized Silicone Gripper Cuff | 12"    |

- ✓ Increases longevity for prosthetic suspension sleeves
- ✓ A simple double layer nylon sleeve worn between the socket and the suspension sleeve
- ✓ Protects the inner surface of the suspension sleeve from abrasive rubbing of the socket brim, because the two nylon layers of the Sleeve Saver™ help absorb friction



Narrow 10"  
Sleeve Saver™



Regular 10"  
Sleeve Saver™



Regular 12" Sleeve  
Saver™ with Elasticized  
Silicone Gripper Cuff

## FOOT SHELL SOCK

### FOOT SHELL SOCK

| ITEM #     | SIZE   | FITS                    |
|------------|--------|-------------------------|
| 1F-MD-2225 | Medium | Fits average 22-25 feet |
| 1F-LG-2630 | Large  | Fits average 26-30 feet |

- ✓ Abrasion resistant
- ✓ Dyneema Core-Spun fibers
- ✓ LYCRA® for a more snug fit
- ✓ Long wear resistance
- ✓ Reinforced heel and toe
- ✓ Shorter top and dark color for better appearance
- ✓ U.S. Patent No. 6,158,253



FOOT SHELL SOCK

# COSMETIC STOCKINGS



- ✓ With elasticized silicone gripper cuff
- ✓ 100% nylon (except elastic)
- ✓ Sold in pairs

When replacing, refer to color-coded thread inside proximal stocking

## SIZING COLOR CODE (INSIDE SEAM)




## A/K COSMETIC STOCKINGS

| SIZE              | 0 NATURAL | I PINK   | II BLUE  | III GREEN | IV YELLOW |
|-------------------|-----------|----------|----------|-----------|-----------|
| Top Circumference | 36 cm     | 36 cm    | 38 cm    | 40 cm     | 40 cm     |
| Leg Length        | 72 cm     | 76 cm    | 78 cm    | 80 cm     | 90 cm     |
| Foot Length       | 20 cm     | 20 cm    | 20.5 cm  | 21 cm     | 22 cm     |
| Caucasian         | 2HB9B140  | 2HB9B141 | 2HB9B142 | 2HB9B143  | 2HB9B144  |
| Brown             | 2HR9B140  | 2HR9B141 | 2HR9B142 | 2HR9B143  | 2HR9B144  |
| Latin             | 2HL9B140  | 2HL9B141 | 2HL9B142 | 2HL9B143  | 2HL9B144  |

## B/K COSMETIC STOCKINGS

| SIZE              | I PINK   | II BLUE  | III GREEN |
|-------------------|----------|----------|-----------|
| Top Circumference | 25 cm    | 25 cm    | 27 cm     |
| Leg Length        | 36 cm    | 41 cm    | 46 cm     |
| Foot Length       | 20 cm    | 21 cm    | 22 cm     |
| Caucasian         | 2HB9B161 | 2HB9B162 | 2HB9B163  |
| Brown             | 2HR9B161 | 2HR9B162 | 2HR9B163  |
| Latin             | 2HL9B161 | 2HL9B162 | 2HL9B163  |

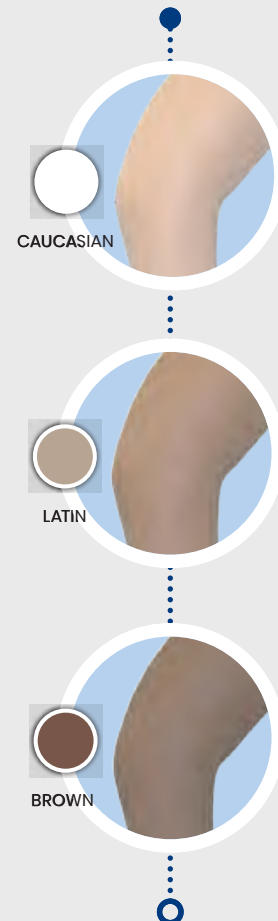

# ORTHOTIC, FRACTURE, *and* COMPRESSION INTERFACES

## ORTHOTIC INTERFACES | 42–51

- 43** Core-Spun by Therafirm® Patterned AFO Socks
- 44** SmartKnit® AFO and KAFO Interface Socks, Child
- 45** SmartKnit® AFO and KAFO Interface Socks, Adult
- 46** SmartKnit® SMO Interface Socks
- 47** SmartKnit® Big Toe Sock
- 48–49** Torso-Interface®
- 50** Cranial Protection
- 51** SmartKnit® Walker Boot Sock
- 51** Hip Orthosis Liner

## FRACTURE INTERFACES | 52–53

- 52** Cast-Rite™ Fracture Cast Sock
- 52** B/K Cast-Rite™ Fracture Brace Sock
- 53** Cast-Rite™ Humeral Sleeve
- 53** Cast-Rite™ Arm Sleeve

## COMPRESSION INTERFACES | 53–55

- 53** 1-Inch Compression Bandage
- 54** SmartKnit® Knee Interface
- 54** Compressogrip® Sleeves
- 55** Compressogrip® Tubular Elastic Bandage



# CORE-SPUN BY THERAFIRM® PATTERNED AFO SOCKS

THESE AFO SOCKS ARE MADE WITH STRETCHY YARNS THAT ENHANCE STRENGTH, DURABILITY, AND STRETCH COMFORT. ADDED THICKNESS PROVIDES ALL-DAY COMFORT AND PROTECTION. THE FORM-FITTING DESIGN HELPS REDUCE IRRITATION FROM WRINKLING AND BUNCHING.

- ✓ Stylish, every day patterns
- ✓ Protective foot and leg interface sock designed for AFO use
- ✓ The non-restrictive, comfort designed top comfortably stays in place all day
- ✓ Superior moisture wicking fibers keep feet cool and dry
- ✓ Smooth toe seam
- ✓ Increases success of orthosis use
- ✓ Physician and orthotist recommended

## CORE-SPUN BY THERAFIRM PATTERNED AFO SOCKS

| ITEM #  | SIZE          | COLOR           | WIDTH / LENGTH       | SOLD   |
|---|---------------|-----------------|----------------------|--------|
| <b>CHILD SIZE AFO INTERFACE SOCKS</b>         |               |                 |                      |        |
| 1SCB-RB-CS                                    | Child Small   | Cheery AFO      | 2.25"x11" (5.5x28cm) | Pair   |
| 1SCB-RB-CR                                    | Child Regular | Cheery AFO      | 2.5"x14" (6.5x36cm)  | Pair   |
| 1SCB-RB-CX                                    | Child X-Long  | Cheery AFO      | 2.75"x20" (7x51cm)   | Pair   |
| 1SCW-GP-CS                                    | Child Small   | Thin Line       | 2.25"x11" (5.5x28cm) | Pair   |
| 1SCW-GP-CR                                    | Child Regular | Thin Line       | 2.5"x14" (6.5x36cm)  | Pair   |
| 1SCW-GP-CX                                    | Child X-Long  | Thin Line       | 2.75"x20" (7x51cm)   | Pair   |
| <b>CHILD SIZE AFO INTERFACE SOCKS, 3-PACK</b> |               |                 |                      |        |
| 1WGB-CS-CS                                    | Child Small   | Halos 3-Pack    | 2.25"x11" (5.5x28cm) | 3-Pair |
| 1WGB-CS-CR                                    | Child Regular | Halos 3-Pack    | 2.5"x14" (6.5x36cm)  | 3-Pair |
| 1WGB-CS-CX                                    | Child X-Long  | Halos 3-Pack    | 2.75"x20" (7x51cm)   | 3-Pair |
| <b>ADULT SIZE AFO INTERFACE SOCKS</b>         |               |                 |                      |        |
| 1SCW-GP-AS                                    | Adult Small   | Classic Diamond | 3"x18" (7.5x46cm)    | Pair   |
| 1SCW-GP-AR                                    | Adult Regular | Classic Diamond | 3.5"x20" (9x51cm)    | Pair   |

Halos 3-Pack, Thin Line, and Cheery AFO available in Child sizes only.  
Classic Diamond available in Adult sizes only.

US Patent No. #7,895,863



Thin Line Core-Spun  
Patterned AFO sock



## CARE TIPS

Machine wash  
warm. Tumble dry  
medium temperature.  
Do not bleach.

Classic Diamond  
AFO sock



Halos  
3-Pack

Cheery  
AFO sock



# SMARTKNIT® SEAMLESS AFO AND KAFO SOCKS, CHILD



EACH AFO AND KAFO SOCK IS MADE THE SAME WAY THAT A CATERPILLAR SPINS ITS COCOON—STARTING AT THE TOE AND CIRCLING OUR WAY TO THE TOP.

The result is a seam-free, wrinkle-free sock that fits like a second layer of skin, provides exceptional softness, and wicks perspiration.

KAFOs longer lengths provide full leg protection.




AFO / KAFO SOCKS



## SMARTKNIT® AFO AND KAFO INTERFACE SOCKS, CHILD

| ITEM #                                   | SIZE          | COLOR    | WIDTH / LENGTH           | SOLD |
|--|---------------|----------|--------------------------|------|
| <b>AFO INTERFACE SOCKS</b>               |               |          |                          |      |
| 1SP1-AF-NS                               | Infant Small  | White    | 1.5" x 7.5" (4 x 19cm)   | Pair |
| 1SP1-AF-NF                               | Infant        | White    | 2" x 9" (5 x 23cm)       | Pair |
| 1SP1-AF-CS                               | Child Small   | White    | 2.25" x 11" (5.5 x 28cm) | Pair |
| 1SP1-AF-CR                               | Child Regular | White    | 2.50" x 14" (6.5 x 36cm) | Pair |
| 1SP1-AF-CX                               | Child X-Long  | White    | 2.75" x 20" (7 x 51cm)   | Pair |
| 1S21-AF-NF                               | Infant        | Pink     | 2" x 9" (5 x 23cm)       | Pair |
| 1S21-AF-CS                               | Child Small   | Pink     | 2.25" x 11" (5.5 x 28cm) | Pair |
| 1S21-AF-CR                               | Child Regular | Pink     | 2.50" x 14" (6.5 x 36cm) | Pair |
| 1S31-AF-NF                               | Infant        | Purple   | 2" x 9" (5 x 23cm)       | Pair |
| 1S31-AF-CS                               | Child Small   | Purple   | 2.25" x 11" (5.5 x 28cm) | Pair |
| 1S31-AF-CR                               | Child Regular | Purple   | 2.50" x 14" (6.5 x 36cm) | Pair |
| 1NX1-AF-NF                               | Infant        | Navy     | 2" x 9" (5 x 23cm)       | Pair |
| 1NX1-AF-CS                               | Child Small   | Navy     | 2.25" x 11" (5.5 x 28cm) | Pair |
| 1NX1-AF-CR                               | Child Regular | Navy     | 2.50" x 14" (6.5 x 36cm) | Pair |
| <b>AFO INTERFACE SOCKS WITH IONIC+™</b>  |               |          |                          |      |
| 1SX1-AF-NF                               | Infant        | Grey     | 2" x 9" (5 x 23cm)       | Pair |
| 1SX1-AF-CS                               | Child Small   | Grey     | 2.25" x 11" (5.5 x 28cm) | Pair |
| 1SX1-AF-CR                               | Child Regular | Grey     | 2.50" x 14" (6.5 x 36cm) | Pair |
| 1SX1-AF-CX                               | Child X-Long  | Grey     | 2.75" x 20" (7 x 51cm)   | Pair |
| 1S51-AF-CS                               | Child Small   | Charcoal | 2.25" x 11" (5.5 x 28cm) | Pair |
| 1S51-AF-CR                               | Child Regular | Charcoal | 2.50" x 14" (6.5 x 36cm) | Pair |
| 1BX1-AF-NF                               | Infant        | Black    | 2" x 9" (5 x 23cm)       | Pair |
| 1BX1-AF-CS                               | Child Small   | Black    | 2.25" x 11" (5.5 x 28cm) | Pair |
| 1BX1-AF-CR                               | Child Regular | Black    | 2.50" x 14" (6.5 x 36cm) | Pair |
| <b>KAFO INTERFACE SOCKS</b>              |               |          |                          |      |
| 1NX1-AF-NF                               | Infant        | Navy     | 2" x 9" (5 x 23cm)       | Pair |
| 1NX1-AF-CS                               | Child Small   | Navy     | 2.25" x 11" (5.5 x 28cm) | Pair |
| 1NX1-AF-CR                               | Child Regular | Navy     | 2.50" x 14" (6.5 x 36cm) | Pair |
| 1SP1-AF-CL                               | Child Long    | White    | 2.50" x 18" (6.5 x 46cm) | Pair |
| 1S21-AF-CL                               | Child Long    | Pink     | 2.50" x 18" (6.5 x 46cm) | Pair |
| 1S31-AF-CL                               | Child Long    | Purple   | 2.50" x 18" (6.5 x 46cm) | Pair |
| 1NX1-AF-CL                               | Child Long    | Navy     | 2.50" x 18" (6.5 x 46cm) | Pair |
| <b>KAFO INTERFACE SOCKS WITH IONIC+™</b> |               |          |                          |      |
| 1SX1-AF-CL                               | Child Long    | Grey     | 2.50" x 18" (6.5 x 46cm) | Pair |
| 1BX1-AF-CL                               | Child Long    | Black    | 2.50" x 18" (6.5 x 46cm) | Pair |

# SMARTKNIT® SEAMLESS AFO AND KAFO SOCKS, ADULT

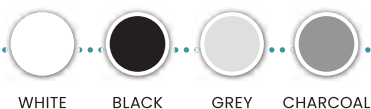

- ✓ No seams means no rubbing on toes
- ✓ Various lengths to allow folding over
- ✓ Wicking fibers keep feet and legs comfortably cool
- ✓ Antimicrobial fibers inhibit odor
- ✓ Heel-less design for a perfect fit
- ✓ Super soft
- ✓ US Patent No. 6,158,253

## CARE TIPS

Machine wash  
warm. Tumble dry  
medium temperature.  
Do not bleach.

## SMARTKNIT® AFO AND KAFO INTERFACE SOCKS, ADULT

| ITEM #                                   | SIZE          | COLOR    | WIDTH / LENGTH                    | SOLD |
|--|---------------|----------|-----------------------------------|------|
| <b>AFO INTERFACE SOCKS — SOLD PAIR</b>   |               |          |                                   |      |
| 1SP1-AF-AS                               | Adult Small   | White    | 3" x 18" (7.5 x 46cm)             | Pair |
| 1SP1-AF-AR                               | Adult Regular | White    | 3.5" x 20" (9 x 51cm)             | Pair |
| <b>AFO INTERFACE SOCKS WITH IONIC+™</b>  |               |          |                                   |      |
| 1SX1-AF-AS                               | Adult Small   | Grey     | 3" x 18" (7.5 x 46cm)             | Pair |
| 1SX1-AF-AR                               | Adult Regular | Grey     | 3.5" x 20" (9 x 51cm)             | Pair |
| 1S51-AF-AS                               | Adult Small   | Charcoal | 3" x 18" (7.5 x 46cm)             | Pair |
| 1S51-AF-AR                               | Adult Regular | Charcoal | 3.5" x 20" (9 x 51cm)             | Pair |
| 1BX1-AF-AS                               | Adult Small   | Black    | 3" x 18" (7.5 x 46cm)             | Pair |
| 1BX1-AF-AR                               | Adult Regular | Black    | 3.5" x 20" (9 x 51cm)             | Pair |
| <b>KAFO INTERFACE SOCKS</b>              |               |          |                                   |      |
| 1SP1-NA-XL                               | Adult Narrow  | White    | 4" x 2.5" x 24" (10 x 6.5 x 61cm) | Each |
| 1SP1-RG-XL                               | Adult Regular | White    | 6" x 3.5" x 28" (15 x 9 x 71cm)   | Each |
| <b>KAFO INTERFACE SOCKS WITH IONIC+™</b> |               |          |                                   |      |
| 1BX1-NA-XL                               | Adult Narrow  | Black    | 4" x 2.5" x 24" (10 x 6.5 x 61cm) | Each |
| 1BX1-RG-XL                               | Adult Regular | Black    | 6" x 3.5" x 28" (15 x 9 x 71cm)   | Each |

### EXTRA ANTIMICROBIAL PROTECTION

**BLACK, CHARCOAL, AND GREY SEAMLESS SOCKS CONTAIN IONIC+™ MINERAL ANTIMICROBIAL FABRIC FOR ITS MOISTURE WICKING PROPERTIES**

Grey AFOs feature  
Ionic+™ Mineral  
Antimicrobial fabric



AFO / KAFO SOCKS .



# SMARTKNIT® SEAMLESS SMO INTERFACE SOCKS



CHARCOAL  
SMARTKNIT® SMO  
INTERFACE SOCK

THESE SUPER SOFT SEAMLESS SOCKS ARE PERFECT FOR CHILDREN WITH LOW PROFILE BRACING (SMOS, UCBLs), BY HELPING TO REDUCE PRESSURE AND INCREASE COMPLIANCE.


- ✓ Made with Core-Spun yarns and LYCRA®
- ✓ No wrinkles or bunching
- ✓ No seams means no rubbing on toes
- ✓ Non-binding top stays in place
- ✓ Various lengths
- ✓ Super soft
- ✓ Moisture-wicking fibers keep feet and legs comfortably cool
- ✓ Heel-less design for a perfect fit and fit various heights
- ✓ Sold as a pair
- ✓ Patent No. 6,158,253.

## SMARTKNIT® SMO INTERFACE SOCKS

| ITEM #                    | COLOR    | WIDTH/LENGTH              |
|---------------------------|----------|---------------------------|
| <b>SMALL (WITH CUFF)</b>  |          |                           |
| 1S11SMSM                  | White    | 2.25" x 5" (5.5 x 13cm)   |
| 1S71SMSM                  | Black    | 2.25" x 5" (5.5 x 13cm)   |
| 1S51SMSM                  | Charcoal | 2.25" x 5" (5.5 x 13cm)   |
| 1S41SMSM                  | Navy     | 2.25" x 5" (5.5 x 13cm)   |
| 1S21SMSM                  | Pink     | 2.25" x 5" (5.5 x 13cm)   |
| 1S31SMSM                  | Purple   | 2.25" x 5" (5.5 x 13cm)   |
| <b>MEDIUM (WITH CUFF)</b> |          |                           |
| 1S11SMMD                  | White    | 2 3/8" x 7" (6 x 18cm)    |
| 1S71SMMD                  | Black    | 2 3/8" x 7" (6 x 18cm)    |
| 1S51SMMD                  | Charcoal | 2 3/8" x 7" (6 x 18cm)    |
| 1S41SMMD                  | Navy     | 2 3/8" x 7" (6 x 18cm)    |
| 1S21SMMD                  | Pink     | 2 3/8" x 7" (6 x 18cm)    |
| 1S31SMMD                  | Purple   | 2 3/8" x 7" (6 x 18cm)    |
| <b>LARGE</b>              |          |                           |
| 1S11SMLG                  | White    | 2.50" x 8.5" (6.5 x 22cm) |
| 1S71SMLG                  | Black    | 2.50" x 8.5" (6.5 x 22cm) |
| 1S51SMLG                  | Charcoal | 2.50" x 8.5" (6.5 x 22cm) |
| 1S41SMLG                  | Navy     | 2.50" x 8.5" (6.5 x 22cm) |
| 1S21SMLG                  | Pink     | 2.50" x 8.5" (6.5 x 22cm) |
| 1S31SMLG                  | Purple   | 2.50" x 8.5" (6.5 x 22cm) |
| <b>X-LARGE</b>            |          |                           |
| 1S11SMXL                  | White    | 2.75" x 10" (7 x 26cm)    |
| 1S71SMXL                  | Black    | 2.75" x 10" (7 x 26cm)    |
| 1S51SMXL                  | Charcoal | 2.75" x 10" (7 x 26cm)    |
| 1S41SMXL                  | Navy     | 2.75" x 10" (7 x 26cm)    |
| 1S21SMXL                  | Pink     | 2.75" x 10" (7 x 26cm)    |
| 1S31SMXL                  | Purple   | 2.75" x 10" (7 x 26cm)    |
| <b>XX-LARGE</b>           |          |                           |
| 1S11SMXX                  | White    | 3" x 12" (7.5 x 31cm)     |
| 1S71SMXX                  | Black    | 3" x 12" (7.5 x 31cm)     |
| 1S51SMXX                  | Charcoal | 3" x 12" (7.5 x 31cm)     |
| 1S41SMXX                  | Navy     | 3" x 12" (7.5 x 31cm)     |
| 1S21SMXX                  | Pink     | 3" x 12" (7.5 x 31cm)     |
| 1S31SMXX                  | Purple   | 3" x 12" (7.5 x 31cm)     |

## CARE TIPS

Machine wash and dry, medium temperature, no bleach.

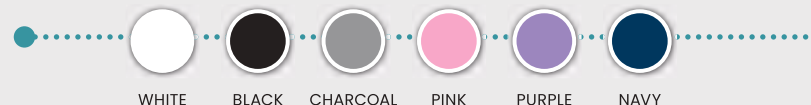

# SMARTKNIT® BIG TOE SOCKS



SMARTKNIT® BIG  
TOE SOCK WITH  
BRACE



SMARTKNIT BIG  
TOE SOCK

BIG TOE AFO SOCKS ARE MADE TO ACCOMMODATE AFO TOE STRAPS, ETC. THESE BIG TOE AFO SOCKS ARE LIKE THEY SHOULD BE—A SEAMLESS, COMFORTABLE, COOL SECOND LAYER OF SKIN BETWEEN YOU AND YOUR AFO.

- ✓ No seams means no rubbing on toes
- ✓ Various lengths to allow folding over
- ✓ Wrinkle-free
- ✓ Super soft
- ✓ Wicking fibers keep feet and legs comfortably cool
- ✓ Antimicrobial fibers inhibit odor
- ✓ Heel-less design for a perfect fit
- ✓ 97.3% Polyester 2.7% LYCRA®
- ✓ Sold as a pair
- ✓ US Patent # 6,158,253

## SMARTKNIT® BIG TOE SOCK

| ITEM #     | SIZE          | COLOR        | WIDTH / LENGTH           |
|------------|---------------|--------------|--------------------------|
| 1BT1-AF-CS | Child Small   | White        | 2.25" x 11" (5.5 x 28cm) |
| 1BT1-AF-CR | Child Regular | White        | 2.50" x 14" (6.5 x 36cm) |
| 1BT1-AF-AS | Adult Small   | White        | 3" x 18" (8 x 46cm)      |
| 1BT1-AF-AR | Adult Regular | White        | 3.50" x 20" (9 x 51cm)   |
| 1BTB-AF-CS | Child Small   | Black/Ionic™ | 2.25" x 11" (5.5 x 28cm) |
| 1BTB-AF-CR | Child Regular | Black/Ionic™ | 2.50" x 14" (6.5 x 36cm) |
| 1BTB-AF-AS | Adult Small   | Black/Ionic™ | 3" x 18" (8 x 46cm)      |
| 1BTB-AF-AR | Adult Regular | Black/Ionic™ | 3.50" x 20" (9 x 51cm)   |



WHITE



BLACK



# TORSO INTERFACE®

Developed by Knit-Rite, the Torso Interface® may be used with all types of spinal orthoses and body jackets. Seamless finishing provides improved skin protection, patient comfort, and aids in preventing skin breakdown.

Slightly heavier than a t-shirt, Torso Interface® provides more cushioning, moisture wicking and anti-microbial properties to keep skin dry, protected, and comfortable.

## CARE TIPS

Hand or machine wash, mild detergent, warm temperature; no bleach. Tumble dry warm or dry flat.




V-Neck Tank



Axilla Flaps



Crew Neck Tank

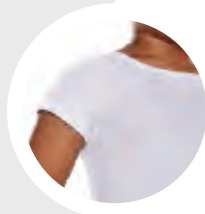

Crew Neck with Sleeves



Strapless Shirt



Turtleneck with Sleeves

IONIC+™ STYLES  
CONTAINS IONIC+™  
MINERAL ANTIMICROBIAL  
FABRIC KNOWN  
FOR ITS MOISTURE  
WICKING PROPERTIES

## ORIGINAL V-NECK TANK

| IONIC+™    | COOLMAX®   | WIDTH/LENGTH            | FIT CIRCUMFERENCE |
|------------|------------|-------------------------|-------------------|
| 1BX1-CV-SM | 1BP1-CV-SM | 6" x 14" (15cm x 35cm)  | 14"-18" (35-45cm) |
| 1BX1-CV-MD | 1BP1-CV-MD | 8" x 16" (20cm x 40cm)  | 18"-24" (45-60cm) |
| —          | 1BP1-CV-ML | 8" x 20" (20cm x 50cm)  | 18"-24" (45-60cm) |
| 1BX1-CV-LG | 1BP1-CV-LG | 10" x 20" (25cm x 50cm) | 22"-30" (55-75cm) |
| 1BX1-CV-LL | 1BP1-CV-LL | 10" x 24" (25cm x 60cm) | 22"-30" (55-75cm) |
| 1BX1-CV-XL | 1BP1-CV-XL | 12" x 24" (30cm x 60cm) | 26"-38" (66-95cm) |

## ORIGINAL V-NECK TANK WITH DOUBLE AXILLA FLAPS

| IONIC+™    | COOLMAX®   | WIDTH/LENGTH            | FIT CIRCUMFERENCE |
|------------|------------|-------------------------|-------------------|
| 1BX1-BD-SM | 1BP1-BD-SM | 6" x 16" (15cm x 40cm)  | 14"-18" (35-45cm) |
| 1BX1-BD-MD | 1BP1-BD-MD | 8" x 18" (20cm x 45cm)  | 18"-24" (45-60cm) |
| 1BX1-BD-ML | 1BP1-BD-ML | 8" x 22" (20cm x 56cm)  | 18"-24" (45-60cm) |
| 1BX1-BD-LG | 1BP1-BD-LG | 10" x 20" (25cm x 50cm) | 22"-30" (55-75cm) |
| 1BX1-BD-LL | 1BP1-BD-LL | 10" x 24" (25cm x 60cm) | 22"-30" (55-75cm) |
| 1BX1-BD-XL | 1BP1-BD-XL | 12" x 24" (30cm x 60cm) | 26"-38" (66-95cm) |

## ORIGINAL V-NECK TANK WITH LEFT AXILLA FLAP

| IONIC+™    | COOLMAX® | WIDTH/LENGTH            | FIT CIRCUMFERENCE |
|------------|----------|-------------------------|-------------------|
| 1BX1-BL-SM | —        | 6" x 14" (15cm x 35cm)  | 14"-18" (35-45cm) |
| 1BX1-BL-MD | —        | 8" x 16" (20cm x 40cm)  | 18"-24" (45-60cm) |
| 1BX1-BL-LG | —        | 10" x 20" (25cm x 50cm) | 22"-30" (55-75cm) |

## ORIGINAL CREW NECK TANK

| IONIC+™    | COOLMAX®   | WIDTH/LENGTH            | FIT CIRCUMFERENCE  |
|------------|------------|-------------------------|--------------------|
| —          | 1BP1-BS-XS | 4" x 12" (10cm x 30cm)  | 10"-12" (25-30cm)  |
| 1BX1-BS-SM | 1BP1-BS-SM | 6" x 14" (15cm x 35cm)  | 14"-18" (35-45cm)  |
| 1BX1-BS-MD | 1BP1-BS-MD | 8" x 16" (20cm x 40cm)  | 18"-24" (45-60cm)  |
| 1BX1-BS-ML | 1BP1-BS-ML | 8" x 22" (20cm x 56cm)  | 18"-24" (45-60cm)  |
| 1BX1-BS-LG | 1BP1-BS-LG | 10" x 20" (25cm x 50cm) | 22"-30" (55-75cm)  |
| 1BX1-BS-LL | 1BP1-BS-LL | 10" x 24" (25cm x 60cm) | 22"-30" (55-75cm)  |
| 1BX1-BS-XL | 1BP1-BS-XL | 12" x 24" (30cm x 60cm) | 26"-38" (66-95cm)  |
| 1BX1-BS-KS | 1BP1-BS-KS | 14" x 24" (35cm x 60cm) | 30"-48" (76-122cm) |
| —          | 1BP1-BS-KX | 14" x 28" (35cm x 71cm) | 30"-48" (76-122cm) |

## ORIGINAL CREW NECK TANK WITH DOUBLE AXILLA FLAPS

| IONIC+™ | COOLMAX®   | WIDTH/LENGTH            | FIT CIRCUMFERENCE |
|---------|------------|-------------------------|-------------------|
| —       | 1BP1-DD-SM | 6" x 14" (15cm x 35cm)  | 14"-18" (35-45cm) |
| —       | 1BP1-DD-MD | 8" x 16" (20cm x 40cm)  | 18"-24" (45-60cm) |
| —       | 1BP1-DD-ML | 8" x 22" (20cm x 56cm)  | 18"-24" (45-60cm) |
| —       | 1BP1-DD-LG | 10" x 20" (25cm x 50cm) | 22"-30" (55-75cm) |
| —       | 1BP1-DD-LL | 10" x 24" (25cm x 60cm) | 22"-30" (55-75cm) |
| —       | 1BP1-DD-XL | 12" x 24" (30cm x 60cm) | 26"-38" (66-95cm) |

## ORIGINAL CREW NECK TANK WITH SLEEVES

| IONIC+™ | COOLMAX®   | WIDTH/LENGTH            | FIT CIRCUMFERENCE |
|---------|------------|-------------------------|-------------------|
| —       | 1BP1-SS-SM | 6" x 16" (15cm x 40cm)  | 14"-18" (35-45cm) |
| —       | 1BP1-SS-MD | 8" x 18" (20cm x 45cm)  | 18"-24" (45-60cm) |
| —       | 1BP1-SS-ML | 8" x 20" (20cm x 50cm)  | 18"-24" (45-60cm) |
| —       | 1BP1-SS-LG | 10" x 20" (25cm x 50cm) | 22"-30" (55-75cm) |
| —       | 1BP1-SS-LL | 10" x 24" (25cm x 60cm) | 22"-30" (55-75cm) |
| —       | 1BP1-SS-XL | 12" x 24" (30cm x 60cm) | 26"-38" (66-95cm) |
| —       | 1BP1-SS-XX | 12" x 28" (30cm x 71cm) | 26"-38" (66-95cm) |



## TORSO INTERFACE®

### ORIGINAL STRAPLESS SHIRT

| IONIC+™ | COOLMAX®   | WIDTH/LENGTH            | FIT CIRCUMFERENCE |
|---------|------------|-------------------------|-------------------|
| —       | 1BP1-B6-14 | 6" x 14" (15cm x 35cm)  | 14"-18" (35-45cm) |
| —       | 1BP1-B8-16 | 8" x 16" (20cm x 40cm)  | 18"-24" (45-60cm) |
| —       | 1BP1-B8-20 | 8" x 20" (20cm x 50cm)  | 18"-24" (45-60cm) |
| —       | 1BP1-10-18 | 10" x 20" (25cm x 50cm) | 22"-30" (55-75cm) |

### ORIGINAL TURTLENECK WITH SLEEVES

| IONIC+™    | COOLMAX® | WIDTH/LENGTH           | FIT CIRCUMFERENCE |
|------------|----------|------------------------|-------------------|
| 1BX1-BT-XS | —        | 6" x 14" (15cm x 35cm) | 10"-12" (25-30cm) |
| 1BX1-BT-SM | —        | 6" x 14" (15cm x 35cm) | 14"-18" (35-45cm) |
| 1BX1-BT-MD | —        | 8" x 16" (20cm x 40cm) | 18"-24" (45-60cm) |

In addition to wearing with all types of spinal orthoses and body jackets, the turtleneck with short sleeves is used in hospitals for a casting procedure on infants and toddlers with scoliosis. It allows parents to be able to take off and put on a clean interface under the cast, with anti-odor Ionic+™ next to the skin.



- ✓ Completely seamless with finished edges
- ✓ High stretch in all directions to accommodate more patients with fewer SKUs
- ✓ Available in white COOLMAX® polyester and grey Ionic+™ styles
- ✓ Contoured, virtually wrinkle-free fit
- ✓ Helps keep skin dry and reduce odors in garment
- ✓ Seamless styling for optimal comfort
- ✓ Moisture wicking fabrics help keep skin dry and odor-free
- ✓ May be used for casting
- ✓ Machine washable

### LIGHTWEIGHT V-NECK TANK

| IONIC+™    | COOLMAX®   | WIDTH/LENGTH            | FIT CIRCUMFERENCE |
|------------|------------|-------------------------|-------------------|
| 1BXL-CV-SM | 1BPL-CV-SM | 6" x 16" (15cm x 40cm)  | 16-27" (41-69cm)  |
| 1BXL-CV-MD | 1BPL-CV-MD | 8" x 18" (20cm x 45cm)  | 22-37" (56-94cm)  |
| 1BXL-CV-ML | 1BPL-CV-ML | 8" x 22" (20cm x 56cm)  | 22-37" (56-94cm)  |
| 1BXL-CV-LG | 1BPL-CV-LG | 10" x 20" (25cm x 50cm) | 30-50" (76-127cm) |
| 1BXL-CV-LL | 1BPL-CV-LL | 10" x 24" (25cm x 60cm) | 30-50" (76-127cm) |
| 1BXL-CV-XL | 1BPL-CV-XL | 12" x 24" (30cm x 60cm) | 36-60" (92-153cm) |

### LIGHTWEIGHT V-NECK TANK WITH LEFT AXILLA FLAP

| IONIC+™    | COOLMAX®   | WIDTH/LENGTH            | FIT CIRCUMFERENCE |
|------------|------------|-------------------------|-------------------|
| 1BXL-BL-SM | 1BPL-BL-SM | 6" x 16" (15cm x 40cm)  | 16-27" (41-69cm)  |
| 1BXL-BL-MD | 1BPL-BL-MD | 8" x 18" (20cm x 45cm)  | 22-37" (56-94cm)  |
| 1BXL-BL-ML | 1BPL-BL-ML | 8" x 22" (20cm x 56cm)  | 22-37" (56-94cm)  |
| 1BXL-BL-LG | 1BPL-BL-LG | 10" x 20" (25cm x 50cm) | 30-50" (76-127cm) |
| 1BXL-BL-LL | 1BPL-BL-LL | 10" x 24" (25cm x 60cm) | 30-50" (76-127cm) |
| 1BXL-BL-XL | 1BPL-BL-XL | 12" x 24" (30cm x 60cm) | 36-60" (92-153cm) |

### LIGHTWEIGHT V-NECK TANK WITH DOUBLE AXILLA FLAPS

| IONIC+™    | COOLMAX®   | WIDTH/LENGTH            | FIT CIRCUMFERENCE |
|------------|------------|-------------------------|-------------------|
| 1BXL-BD-SM | 1BPL-BD-SM | 6" x 16" (15cm x 40cm)  | 16-27" (41-69cm)  |
| 1BXL-BD-MD | 1BPL-BD-MD | 8" x 18" (20cm x 45cm)  | 22-37" (56-94cm)  |
| 1BXL-BD-ML | 1BPL-BD-ML | 8" x 22" (20cm x 56cm)  | 22-37" (56-94cm)  |
| 1BXL-BD-LG | 1BPL-BD-LG | 10" x 20" (25cm x 50cm) | 30-50" (76-127cm) |
| 1BXL-BD-LL | 1BPL-BD-LL | 10" x 24" (25cm x 60cm) | 30-50" (76-127cm) |
| 1BXL-BD-XL | 1BPL-BD-XL | 12" x 24" (30cm x 60cm) | 36-60" (92-153cm) |

## LIGHTWEIGHT TORSO INTERFACE®

DESIGNED WITH THE SAME GREAT BENEFITS OF OUR ORIGINAL TORSO INTERFACE® FOR PATIENTS WHO PREFER A LIGHTER WEIGHT GARMENT.

Lightweight Grey Ionic+™ Torso Interface® with Double Axilla Flaps



# CRANIAL PROTECTION

If skin eruptions have occurred, changing interfaces two or three times a day and washing both head area and garments daily are helpful.



MULLET WITH  
VELCRO CHIN CLOSURE

**DESIGNED FOR CHILDREN WEARING AN ORTHOSIS TO TREAT POSITIONAL PLAGIOCEPHALY AND OTHER HEAD SHAPE DEFORMITIES.**

- ✓ Used to protect the skin and wick moisture away from the skin
- ✓ Stretch allows one size to fit heads of 3 to 18 month old infants

## MULLET CRANIAL INTERFACE WITH VELCRO CHIN CLOSURE

| ITEM #    | SIZE              | MATERIAL | SOLD |
|-----------|-------------------|----------|------|
| 1XPC-MULL | One size fits all | Ionic+™  | Each |

- ✓ Ionic+™ guards against dermatitis
- ✓ Velcro chin overlap closure aids in helmet donning
- ✓ After interface is pulled on, velcro closure fastens under chin to help eliminate wrinkles

**CAUTION:** Leave open while infant is sleeping. velcro closure should be released once helmet is in place for patient comfort.



BALACLAVA

## BALACLAVA CRANIAL INTERFACE FOR CASTING OR SCANNING IMPRESSIONS

| ITEM #    | SIZE              | MATERIAL              | SOLD |
|-----------|-------------------|-----------------------|------|
| 1SP1-BACL | One size fits all | COOLMAX®/Nylon/LYCRA® | Each |

- ✓ Position on head and trim to open up areas as needed
- ✓ Length allows for chin strap if desired

## BALACLAVA SCANNING INTERFACE FOR SCANNING OR FOR USE AS AN INTERFACE

| ITEM #    | SIZE              | MATERIAL                | SOLD |
|-----------|-------------------|-------------------------|------|
| 1CN1-R410 | One size fits all | Microfiber Nylon/LYCRA® | Each |

- ✓ Provides a soft feel in a thinner version

## CRANIAL SCANNING SOCK FOR SCANNING IMPRESSIONS FROM INFANT TO ADULT

| ITEM #     | SIZE     | MATERIAL      | SOLD           |
|------------|----------|---------------|----------------|
| 2C-RANL-03 | 3" x 12" | Nylon/Spandex | 25 per package |

- ✓ Length allows for chin strap if desired
- ✓ Position on head and trim to open areas as needed



CRANIAL  
SCANNING SOCK



# SMARTKNIT® WALKER BOOT SOCK



IONIC+™ SMARTKNIT®  
WALKER BOOT SOCK

THE SMARTKNIT BOOT INTERFACE OFFERS UNIQUE PROTECTION SPECIFICALLY FOR THE WALKING BOOT ENVIRONMENT AND ASSOCIATED LONGER WEAR TIMES.

- ✓ Ionic+™ Mineral Antimicrobial fibers assist in inhibiting odor in the sock
- ✓ The Ionic+™ fibers wick moisture to keep legs drier, healthier and more comfortable
- ✓ Core-Spun yarns provide the unique ability to conform to the anatomy, reducing the risks from pressure-causing wrinkles
- ✓ Features a heel pocket design to ensure a fit with less slipping on the foot
- ✓ Made the same way that a caterpillar spins its cocoon—starting at the toe and circling its way to the top, resulting in a seam-free sock

## CARE TIPS

Machine wash and dry, medium temperature, no bleach.

## SMARTKNIT® WALKER BOOT SOCK

| REGULAR | WIDE  | SIZE    | MEN'S SHOE | WOMEN'S SHOE |
|---------|-------|---------|------------|--------------|
| 73624   | 73664 | Small   | 5–7        | 5–8          |
| 73634   | 73674 | Medium  | 8–10       | 9–11         |
| 73644   | 73684 | Large   | 11–13      | 12–14        |
| 73654   | 73694 | X-Large | 14+        | —            |

Wide size fits circumference 23"–25".

## HIP ORTHOSIS LINER

| ITEM #     | SIZE      | MATERIAL       | SOLD |
|------------|-----------|----------------|------|
| 4HL1-12-34 | 12" x 34" | Acrylic/LYCRA® | Each |

- ✓ Fits most adults
- ✓ Useful with every type of hip orthosis
- ✓ Made of acrylic/LYCRA®
- ✓ Stretch allows the longer leg to be positioned first and the garment to be extended and relaxed carefully over sensitive hip areas
- ✓ The wicking effect of the acrylic fiber moves moisture away from the skin surface
- ✓ 1-ply weight fabric, heavy enough to protect the skin from brace parts
- ✓ Easy wash garment
- ✓ Design provides adequate openings so that no change of garment is required during the day
- ✓ Sold each

## HIP ORTHOSIS LINER




HIP ORTHOSIS  
LINER

## CARE TIPS

Machine wash warm, tumble dry, no bleach



# CAST-RITE™ FRACTURE CAST SOCK

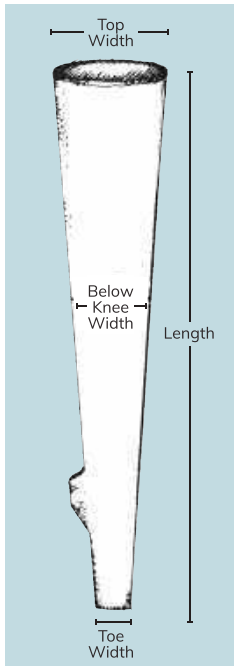

## CAST-RITE™ FRACTURE CAST SOCK FOR FRACTURE MANAGEMENT

| ACRYLIC/LYCRA®    | SIZE    | LENGTH INCLUDING FOOT | TOP WIDTH   | TOE WIDTH      | FITS CALF CIRCUMFERENCE |
|-------------------|---------|-----------------------|-------------|----------------|-------------------------|
| <b>ABOVE KNEE</b> |         |                       |             |                |                         |
| 1FL2-AC-MD        | Medium  | 30" (75cm)            | 7" (17.5cm) | 3" (7.5cm)     | 8"-16" (20-40cm)        |
| 1FL2-AC-LG        | Large   | 34" (85cm)            | 9" (23cm)   | 4" (10cm)      | 11"-22" (28-55cm)       |
| <b>BELOW KNEE</b> |         |                       |             |                |                         |
| 1FL2-BC-XS        | X-Small | 16" (40cm)            | 2" (7.5cm)  | 1.50" (3.75cm) | 2"-4" (5-10cm)          |
| 1FL2-BC-MD        | Medium  | 22" (55cm)            | 5" (12.5cm) | 3" (7.5cm)     | 8"-16" (20-40cm)        |
| 1FL2-BC-LG        | Large   | 26" (65cm)            | 7" (17.5cm) | 4" (10cm)      | 11"-22" (28-55cm)       |

- ✓ Fashion knitted with heel for better fit
- ✓ Core-Spun yarn helps prevent wrinkles
- ✓ Full foot design for procedures where cast encompasses all or part of the foot
- ✓ Toe may be left closed for foot protection, sock may be trimmed distally and/or proximally and folded back into plaster or other casting material, or extra length may be folded down over top of cast
- ✓ Foot may be trimmed off so sock can be used for cylinder cast
- ✓ Can be used with most plastic and plaster casting materials
- ✓ 2-Ply
- ✓ Stretches from 100%
- ✓ May be gas autoclaved — do not steam autoclave
- ✓ Packaged one per envelope
- ✓ 97% acrylic – 3% LYCRA® spandex
- ✓ Fleece inside
- ✓ A/K and B/K

## B/K CAST-RITE™ FRACTURE BRACE SOCK

- ✓ For fracture management
- ✓ Open toe design reduces bulk so that shoe may be worn
- ✓ Gradual taper and heel pocket provide smooth fit
- ✓ Top and toe finishing
- ✓ Over 100% stretchable
- ✓ 97% acrylic – 3% LYCRA®
- ✓ Washable
- ✓ White
- ✓ 2-Ply
- ✓ Sold each

### B/K CAST-RITE™ FRACTURE BRACE SOCK

| ITEM #     | SIZE  | FITS CALF CIRCUMFERENCE | LENGTH INCLUDING FOOT | TOP WIDTH   | TOE WIDTH |
|------------|-------|-------------------------|-----------------------|-------------|-----------|
| 1FL2-BB-LG | Large | 11"-22" (28-55cm)       | 24" (60cm)            | 7" (17.5cm) | 4" (10cm) |



B/K Cast-Rite™  
Fracture Brace Sock

# CAST-RITE™ HUMERAL SLEEVE



Cast-Rite™  
Humeral Sleeve

## CAST-RITE™ HUMERAL SLEEVE

| ITEM #     | LENGTH     | FLAT WIDTH  |
|------------|------------|-------------|
| 1AL2-HM-MD | 12" (30cm) | 3.50" (9cm) |

- ✓ 2-ply cushioning weight
- ✓ 97% acrylic, 3% LYCRA® spandex
- ✓ Washable
- ✓ Sold each

# CAST-RITE™ ARM SLEEVE

## CAST-RITE™ ARM SLEEVE WITH THUMB OR THUMBHOLE

| ITEM #                              | STYLE      | LENGTH     | FLAT WIDTH |
|-------------------------------------|------------|------------|------------|
| ADULT SIZE SLEEVES – SOLD EACH      |            |            |            |
| 1AL1-AT-10                          | With Thumb | 10" (25cm) | 3" (7.5cm) |
| 1AL1-A0-10                          | Thumbhole  | 10" (25cm) | 3" (7.5cm) |
| 1AL1-A0-18                          | Thumbhole  | 18" (45cm) | 3" (7.5cm) |
| CHILD SIZE SLEEVES – SOLD AS A PAIR |            |            |            |
| 1AL1-A0-CH                          | Thumbhole  | 10" (25cm) | 2" (5cm)   |

- ✓ A comfortable, moisture-wicking interface for casts or braces protects skin against irritation and abrasions
- ✓ Acrylic/LYCRA® fabric
- ✓ Washable material stretches and conforms to arm contour
- ✓ Available in 1-ply protective weight



ARM SLEEVE WITH  
THUMB



ARM SLEEVE WITH  
THUMBHOLE

# 1-INCH COMPRESSION BANDAGE



## 1-INCH COMPRESSION BANDAGE

| ITEM #     | WIDTH/LENGTH                | COMPRESSION |
|------------|-----------------------------|-------------|
| 6S-7630-06 | 1" (2.5cm) x 5 yards (4.6M) | Medium      |
| 6S-7630-07 | 1" (2.5cm) x 5 yards (4.6M) | Light       |

- ✓ Washable knitted three dimensional material creates openness and increases healing
- ✓ Use as a single or double bandage for digits
- ✓ 52% nylon; 29% rayon; 19% spandex
- ✓ Latex free materials
- ✓ Sold per 5 yard (4.5m) roll

## SMARTKNIT® KNEE INTERFACE

- ✓ Multidirectional stretch moves with body during activity
- ✓ Can be worn alone or for compressive support
- ✓ Less bulky and better circulation than neoprene garment
- ✓ Latex-free materials
- ✓ Moisture wicking COOLMAX® or Ionic+™ Mineral Antimicrobial helps skin stay cool, dry, and comfortable
- ✓ Ionic+™ Mineral Antimicrobial helps inhibit odor in garment



### SMARTKNIT® KNEE INTERFACE WITH STRAIGHT SLEEVES

| WHITE COOLMAX®/LYCRA  | GREY POLYESTER/IONIC+™ | BLACK POLYESTER/IONIC+™ | FLAT WIDTH | THIGH CIRCUMFERENCE | LENGTH     |
|-----------------------|------------------------|-------------------------|------------|---------------------|------------|
| WITH STRAIGHT SLEEVES |                        |                         |            |                     |            |
| —                     | —                      | 2PCX-04-18              | 4" (10cm)  | 11"-14" (28-36cm)   | 18" (45cm) |
| —                     | —                      | 2PCX-05-18              | 5" (13cm)  | 14"-17" (36-43)     | 18" (45cm) |
| —                     | —                      | 2PCX-06-18              | 6" (15cm)  | 17"-20" (43-50cm)   | 18" (45cm) |
| —                     | —                      | 2PCX-07-18              | 7" (18cm)  | 20"-23" (50-58cm)   | 18" (45cm) |
| —                     | —                      | 2PCX-08-18              | 8" (20cm)  | 23"-29" (58-74cm)   | 20" (50cm) |
| —                     | —                      | 2PCX-09-18              | 9" (23cm)  | 26"-34" (66-86cm)   | 20" (50cm) |
| WITH TAPERED SLEEVES  |                        |                         |            |                     |            |
| —                     | 2BCX-T5-18             | 2PCX-T5-18              | 5" x 4"    | 14"-17" (36-43)     | 18" (45cm) |
| 2PCS-T6-18            | 2BCX-T6-18             | 2PCX-T6-18              | 6" x 4"    | 17"-20" (43-50cm)   | 18" (45cm) |
| 2PCS-T7-18            | 2BCX-T7-18             | 2PCX-T7-18              | 7" x 5"    | 20"-23" (50-58cm)   | 18" (45cm) |
| 2PCS-T8-18            | 2BCX-T8-18             | 2PCX-T8-18              | 8" x 6"    | 23"-29" (58-74cm)   | 18" (45cm) |
| 2PCS-T9-18            | —                      | 2PCX-T9-18              | 9" x 7"    | 26"-34" (66-86cm)   | 18" (45cm) |

## COMPRESSOGRIPT® SLEEVES



### CARE TIPS

Machine wash and dry, warm.  
Do not bleach.

A UNIQUE ELASTIC TUBULAR BANDAGE FOR WRISTS, ELBOWS, ARMS, ANKLES, KNEES, LEGS, AND TORSO THAT APPLIES COMPRESSION FOR COMFORT AND SUPPORT UNDER BRACES.

### COMPRESSOGRIPT® SLEEVES

| ITEM #     | SIZE | FITS CIRCUMFERENCE | LENGTH |
|------------|------|--------------------|--------|
| 2ECA-03-18 | 3    | 9"-12" (22.5-30cm) | 18"    |
| 2ECA-04-18 | 4    | 12"-16" (30-40cm)  | 18"    |
| 2ECA-05-18 | 5    | 15"-20" (38-50cm)  | 18"    |
| 2ECA-06-18 | 6    | 18"-24" (45-60cm)  | 18"    |
| 2ECA-08-18 | 8    | 24"-32" (60-80cm)  | 18"    |
| 2ECA-04-24 | 4    | 12"-16" (30-40cm)  | 24"    |
| 2ECA-05-24 | 5    | 15"-20" (38-50cm)  | 24"    |
| 2ECA-06-24 | 6    | 18"-24" (45-60cm)  | 24"    |
| 2ECA-08-24 | 8    | 24"-32" (60-80cm)  | 24"    |

- ✓ Black
- ✓ Used for controlling and reducing edema, and to minimize hypertrophic scarring
- ✓ Finished edges permit washing and reuse without fraying or raveling
- ✓ Compression is 20-25mmHg at a 50% stretch, and doubled at a 100% stretch
- ✓ 92% Polyester, 8% spandex
- ✓ Latex-free
- ✓ Machine washable
- ✓ Easy removal and reapplication



# COMPRESSOGRIP® TUBULAR ELASTIC BANDAGE

- ✓ More consistent compression than elastic wrap
- ✓ Easy to apply
- ✓ Conforms to contour of limb
- ✓ Proven effective for edema reduction and control, interim control of hypertrophic scarring following burn surgery, dressing and monitoring device retention
- ✓ Smooths irregularities when casting
- ✓ Available in a wide range of sizes and cost effective
- ✓ Moisture-wicking Ionic+™ fibers wick away moisture and an antimicrobial finish helps inhibit odor in the sock
- ✓ Sold in 11-yard (10 meter) rolls
- ✓ Single layer supportive (20–25mmHg) compression; double layer provides therapeutic (35–50mmHg) compression at 50% stretch

## COMPRESSOGRIP® TUBULAR ELASTIC BANDAGE

| RE-ORDER CODES | SIZE | *FITS CIRCUMFERENCE    | POLYESTER WHITE | POLYESTER BEIGE | POLYESTER BLACK | IONIC+™ GREY |
|----------------|------|------------------------|-----------------|-----------------|-----------------|--------------|
| 2" (A)         |      | 5-7" (12.5-17.5cm)     | 6S-7628-02      | 6S-7631-02      | —               | 6S-76XS-02   |
| 2.5" (B)       |      | 7.5-8" (19-20cm)       | 6S-7628-S2      | 6S-7631-S2      | —               | 6S-76XS-S2   |
| 2.75" (C)      |      | 8.5-8.75" (21-22cm)    | 6S-7628-T2      | 6S-7631-T2      | —               | —            |
| 3" (D)         |      | 9-10" (23-25.5cm)      | 6S-7628-03      | 6S-7631-03      | 6S-7633-03      | 6S-76XS-03   |
| 3.5" (E)       |      | 10.5-11.5" (26.5-29cm) | 6S-7628-S3      | 6S-7631-S3      | —               | 6S-76XS-S3   |
| 4" (F)         |      | 12-14" (30.5-35.5cm)   | 6S-7628-04      | 6S-7631-04      | 6S-7633-04      | 6S-76XS-04   |
| 5" (G)         |      | 15-20" (38-51cm)       | 6S-7628-05      | 6S-7631-05      | 6S-7633-05      | 6S-76XS-05   |
| 6" (J)         |      | 18-24" (46-61cm)       | 6S-7628-06      | 6S-7631-06      | 6S-7633-06      | 6S-76XS-06   |
| 8" (K)         |      | 24-32" (61-81cm)       | 6S-7628-08      | 6S-7631-08      | 6S-7633-08      | —            |
| 10" (L)        |      | 30-40" (76-102cm)      | 6S-7628-10      | 6S-7631-10      | —               | —            |
| 12" (M)        |      | 36-48" (91-122cm)      | 6S-6628-12*     | 6S-6631-12*     | —               | —            |
| 14" (N)        |      | 43-57" (110-145cm)     | 6S-6628-14*     | —               | —               | —            |
| 16" (P)        |      | 48"-64" (122-163cm)    | 6S-6628-16*     | —               | —               | —            |

\*Bagged Only

## MATERIALS

White, Beige, Black — 92% polyester, 8% spandex

Polyester/Ionic+™ Grey — 87.4% polyester, 4.6% Ionic+™, 8% spandex

## SIZE CHART

| SIZE | FLAT WIDTH    | FITS CIRCUMFERENCE    |
|------|---------------|-----------------------|
| A    | 2" (5.1cm)    | 5"-7" (12.5-17.5cm)   |
| B    | 2.50" (6.4cm) | 7"-10" (18-25cm)      |
| C    | 2.75" (7.0cm) | 8"-11" (20-28cm)      |
| D    | 3" (7.6cm)    | 9"-12" (23-30cm)      |
| E    | 3.50" (8.9cm) | 10"-13" (25-33cm)     |
| F    | 4" (10.2cm)   | 12"-16" (30.5-40.5cm) |
| G    | 5" (12.7cm)   | 15"-20" (38-51cm)     |
| J    | 6" (15.2cm)   | 18"-24" (46-61cm)     |
| K    | 8" (20.3cm)   | 24"-32" (61-81cm)     |
| L    | 10" (25.4cm)  | 30"-40" (76-102cm)    |
| M    | 12" (30.5cm)  | 36"-48" (91-122cm)    |
| N    | 14" (35.6cm)  | 42"-56" (107-142cm)   |
| P    | 16" (40.6cm)  | 48"-64" (122-163cm)   |



Compressogrip® Tubular Elastic Bandage



## CARE TIPS

Machine wash and dry, warm. Do not bleach. Latex-free. Fabric will relax 20% in length after first washing. Autoclaveable.




# DIABETIC *and* SENSITIVITY SOCKS

## DIABETIC AND SENSITIVITY SOCKS | 56-59


- 57 SmartKnit® Seamless Diabetic Socks Information
- 58 SmartKnit® Seamless Diabetic Socks
- 59 SmartKnit® Seamless Wide Crew Diabetic Socks



# SMARTKNIT® SEAMLESS DIABETIC SOCKS

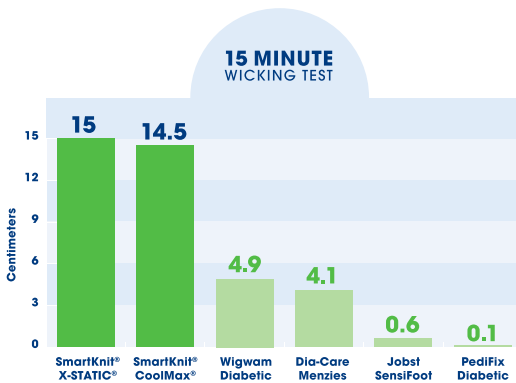

In a 15 minute moisture wicking test, SmartKnit® Seamless Diabetic Socks wick moisture 3–150 times faster than other popular diabetic sock brands. SmartKnit® with Ionic+™ socks wicked moisture to 15 centimeters in only 6 minutes.

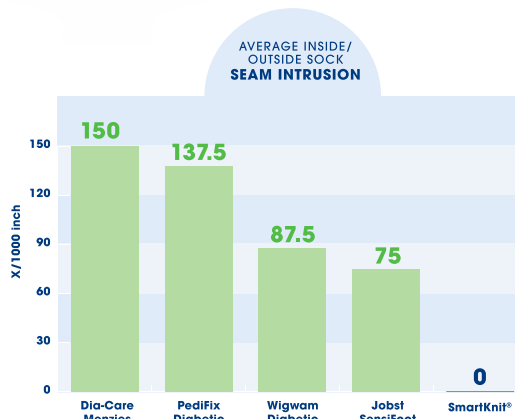

SmartKnit® socks feature 6–150 times less seam protrusion than other popular diabetic sock brands.

Source:  
Precision Testing Labs  
April 7, 2005 Lab  
Report #71507-B

## MOISTURE CONTROL

SmartKnit® Seamless Socks are made with superior moisture wicking yarns. Polyester and Ionic+™ silver fibers wick the moisture away from the skin to the exterior of the sock, moving the moisture up and out of the foot area. The result is a sock that effectively controls moisture and helps reduce blister-causing friction.

## TRULY SEAMLESS

Truly seamless socks—they do exist! SmartKnit® Seamless Socks are made using a knitting process much like how a cocoon is spun, starting at one point and spinning upwards, thus completely eliminating irritating seams and the risks associated with them.



**CORE-SPUN YARNS ARE STRETCHY AND COMFORTABLY HUG YOUR FEET FOR MAXIMUM PROTECTION**

## REDUCES RISKS FROM WRINKLING AND BUNCHING

SmartKnit® Seamless Socks patented knitting process uses high-stretch Core-Spun yarns that “hug” your feet for maximum protection. SmartKnit® socks are not bulky to assure adequate room in the shoe.

## NON-BINDING

SmartKnit® Seamless Socks feature a non-binding Halo Top™ that fits comfortably, minimizing indentations. The Core-Spun yarns stretch to allow a comfortable fit on a range of leg circumferences.

DIABETIC SOCKS



# SMARTKNIT® SEAMLESS DIABETIC SOCKS



White Polyester Crew Socks

THESE SOCKS ARE RECOMMENDED FOR PEOPLE WITH DIABETES, ARTHRITIS, OR ANY CONDITION WHERE PROTECTION FROM IRRITATION IS IMPORTANT.

- ✓ Completely free of seams
- ✓ Uses high-stretch core-spun yarns that hug your feet to help reduce risks from wrinkling and bunching
- ✓ Moisture-wicking fibers wick away moisture from the skin
- ✓ Features a non-binding Halo-Top™ that fits comfortably, minimizing indentations
- ✓ Great casual socks for work or play
- ✓ Perfect for people with diabetes or sensitive feet who are looking for luxuriously soft, non-irritating socks
- ✓ Latex free, 97% Polyester, 3% LYCRA®.

## SIZE CHART

| SIZE   | WOMEN'S SHOE | MEN'S SHOE |
|--------|--------------|------------|
| SMALL  | 5-8          | 5-7        |
| MEDIUM | 8.5-10.5     | 7.5-9.5    |
| LARGE  | 11.5+        | 10+        |

## SMARTKNIT® SEAMLESS DIABETIC SOCKS — SINGLE PAIR

| SIZE    | STYLE         | IONIC+™ BLACK | IONIC+™ GREY | POLYESTER NAVY | POLYESTER WHITE | POLYESTER BLACK | POLYESTER CHARCOAL |
|---------|---------------|---------------|--------------|----------------|-----------------|-----------------|--------------------|
| Small   | Mini-Crew     | —             | 72622        | —              | 72122           | 72722           | —                  |
| Medium  | Mini-Crew     | —             | 72632        | —              | 72132           | 72732           | —                  |
| Large   | Mini-Crew     | —             | 72642        | —              | 72142           | 72742           | —                  |
| X-Large | Mini-Crew     | —             | 72652        | —              | 72152           | 72752           | —                  |
| Small   | Crew          | 71722         | 71622        | 71422          | 71122           | —               | 71522              |
| Medium  | Crew          | 71732         | 71632        | 71432          | 71132           | —               | 71532              |
| Large   | Crew          | 71742         | 71642        | 71442          | 71142           | —               | 71542              |
| X-Large | Crew          | 71752         | 71652        | 71452          | 71152           | —               | 71552              |
| Small   | Over-The-Calf | 73722         | —            | 73422          | 73122           | —               | —                  |
| Medium  | Over-The-Calf | 73732         | —            | 73432          | 73132           | —               | —                  |
| Large   | Over-The-Calf | 73742         | —            | 73442          | 73142           | —               | —                  |
| X-Large | Over-The-Calf | 73752         | —            | 73452          | 73152           | —               | —                  |



Ionic+™ Grey Mini-Crew Socks



Navy Polyester Over-The-Calf Socks



## 3-PACK — MINI CREW

| SIZE    | POLYESTER BLACK | POLYESTER WHITE |
|---------|-----------------|-----------------|
| Small   | 72BBB22         | 72WWW22         |
| Medium  | 72BBB32         | 72WWW32         |
| Large   | 72BBB42         | 72WWW42         |
| X-Large | 72BBB52         | 72WWW52         |

## 3-PACK — CREW

| SIZE    | IONIC+™ BLACK | POLYESTER WHITE |
|---------|---------------|-----------------|
| Small   | 71BBB22       | 71WWW22         |
| Medium  | 71BBB32       | 71WWW32         |
| Large   | 71BBB42       | 71WWW42         |
| X-Large | 71BBB52       | 71WWW52         |

## 3-PACK — OVER-THE-CALF

| SIZE    | IONIC+™ BLACK | POLYESTER WHITE |
|---------|---------------|-----------------|
| Small   | 73BBB22       | 73WWW22         |
| Medium  | 73BBB32       | 73WWW32         |
| Large   | 73BBB42       | 73WWW42         |
| X-Large | 73BBB52       | 73WWW52         |



# SMARTKNIT® SEAMLESS WIDE CREW DIABETIC SOCKS

## SPECIALLY DESIGNED FOR LARGER LEGS AND FEET.

- ✓ Fit like a second skin and stretch up to 25" in circumference
- ✓ Features a non-binding Halo-Top™ that fits comfortably, and won't cut into the legs or inhibit circulation
- ✓ High-stretch core-spun yarns "hug" your feet to help reduce risks from wrinkling and bunching
- ✓ Moisture-wicking fibers wick away moisture and an antimicrobial finish helps inhibit odor in the sock
- ✓ Latex free, 97% Polyester, 3% LYCRA®
- ✓ Great casual socks for work or play
- ✓ Perfect for people with diabetes or sensitive feet who are looking for luxuriously soft, non-irritating socks



WHITE POLYESTER  
WIDE CREW SOCKS



IONIC+™ BLACK WIDE  
CREW SOCKS

## SMARTKNIT® SEAMLESS DIABETIC SOCKS — WIDE CREW

| IONIC+™<br>BLACK | POLYESTER<br>WHITE | STYLE     | SIZE    | WOMEN'S<br>SHOE | MEN'S<br>SHOE |
|------------------|--------------------|-----------|---------|-----------------|---------------|
| 71762            | 71162              | Wide Crew | Small   | 5–8W            | 5–7W          |
| 71772            | 71172              | Wide Crew | Medium  | 9–11W           | 8–10W         |
| 71782            | 71182              | Wide Crew | Large   | 12–14W          | 11–13W        |
| 71792            | 71192              | Wide Crew | X-Large | 14W+            | 14W+          |



IONIC+™ STYLES  
CONTAINS IONIC+™  
MINERAL ANTIMICROBIAL  
FABRIC KNOWN FOR  
ITS MOISTURE WICKING  
PROPERTIES

## SMARTKNIT® SEAMLESS DIABETIC SOCKS WITH MINERAL ANTIMICROBIAL IONIC+™ FIBERS PROVIDE MAXIMUM PROTECTION.

- ✓ Safe, all-natural and lasts the life of the product
- ✓ Antimicrobial to inhibit odor in the sock
- ✓ Moisture-wicking
- ✓ Keeps feet warmer in winter, and cooler in the summer
- ✓ Latex free, 97% Polyester/Ionic+™ blend, 3% LYCRA®

DIABETIC SOCKS




# FABRICATION

*and* CASTING

## FABRICATION | 60-63

- 61** Vacuum Forming Tubing and Hosiery
- 62-63** Stockinettes and Casting

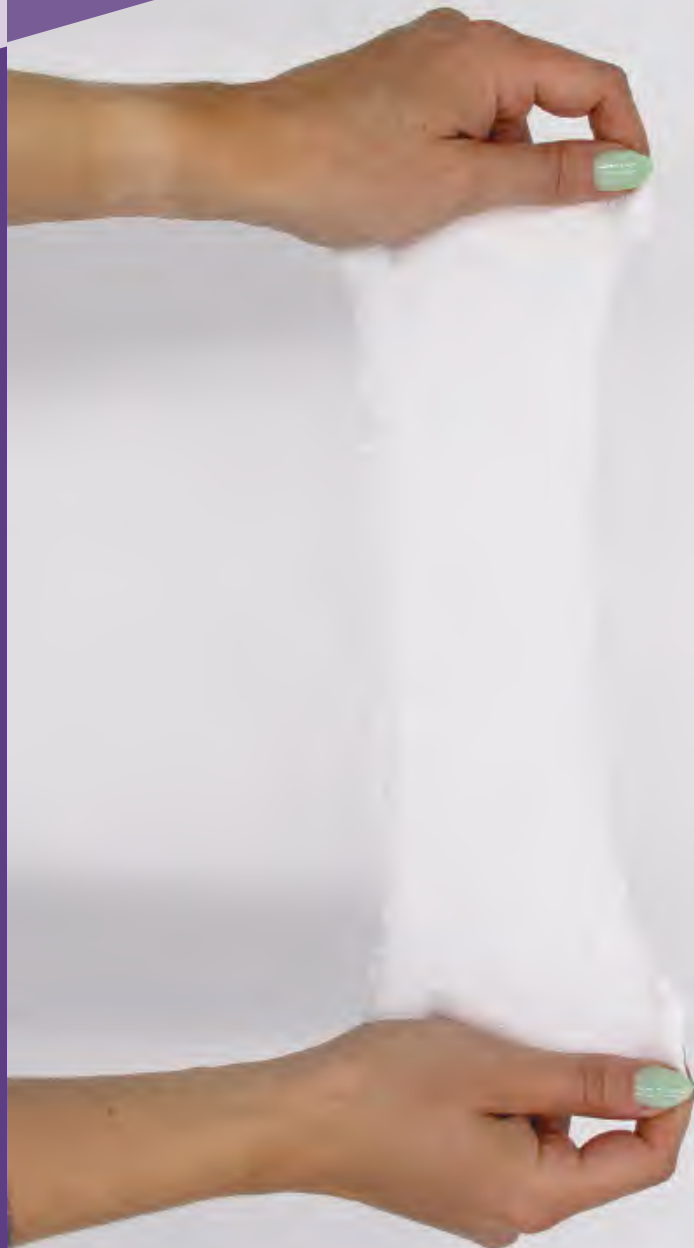

# VACUUM FORMING TUBING AND HOSIERY

A SHEER HOSIERY TUBING USED FOR A VACUUM FORMING INTERFACE ON MODELS AS LARGE AS BODY JACKETS THAT REQUIRE A SMOOTH INNER SURFACE.

- ✓ Stretches up to an incredible 750%
- ✓ Perfect for vacuum forming items with small to a large progression of circumference
- ✓ Sold in boxed 50 yard rolls



## NYLON VACUUM FORMING TUBING


| ITEM #     | STRETCH CIRCUMFERENCE | SIZE                 |
|------------|-----------------------|----------------------|
| 2H-R300-NY | Up to 750%            | 3.50" x 50 yard roll |

## NYSERTS VACUUM FORMING HOSIERY


| ITEM # | TOE    | SIZE                       | STRETCH CIRCUMFERENCE |
|--------|--------|----------------------------|-----------------------|
| 30000  | Closed | 3.75" x 12" (9.5 x 30.5cm) | 8"-20" (20-50cm)      |
| 30001  | Closed | 3.75" x 24" (9.5 x 60cm)   | 8"-20" (20-50cm)      |
| 30002  | Closed | 3.75" x 36" (9.5 x 90cm)   | 8"-20" (20-50cm)      |
| 30003  | Open   | 3.75" x 36" (9.5 x 90cm)   | 8"-20" (20-50cm)      |

MANUFACTURED ESPECIALLY TO PROVIDE VACUUM FACILITATION DURING THERMOFORMING AND LAMINATING PROCEDURES.

### LENGTH COLOR CODE (TOP SEAM)

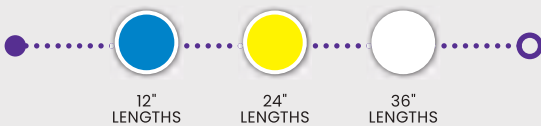

- ✓ 100% stretch nylon available in 12", 24" and 36" lengths
- ✓ 3.50-4" (10cm) unstretched width
- ✓ Sold 200/bag

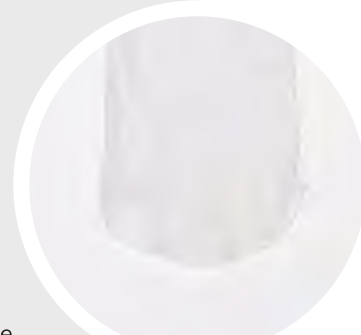


Vacuum Forming Hosiery with sewn toe

## SHEER VACUUM FORMING HOSIERY

| ITEM # | TOE    | SIZE                  | STRETCH CIRCUMFERENCE |
|--------|--------|-----------------------|-----------------------|
| 31000  | Closed | 2.50 x 20" (6 x 50cm) | 6"-20" (15-50cm)      |
| 31001  | Closed | 3.50 x 24" (9 x 60cm) | 8"-24" (20-60cm)      |
| 31002  | Closed | 4 x 30" (10x76cm)     | 9"-30" (24-75cm)      |

### LENGTH COLOR CODE (TOP SEAM)

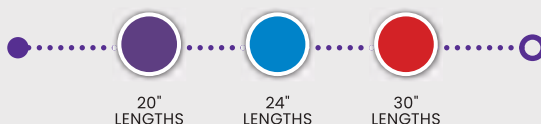

- ✓ 100% stretch nylon available in 20", 24" and 30" lengths
- ✓ Maximum stretch and sheerness
- ✓ Sewn toe
- ✓ Sold 200/bag



12" Vacuum Forming Hosiery with sewn toe



# STOCKINETTES AND CASTING



STANDARD RIB-KNIT STOCKINETTE

## STANDARD RIB-KNIT STOCKINETTE

| NATURAL COTTON | WIDTH        |
|----------------|--------------|
| 6S-7608-01     | 1" (2.5cm)   |
| 6S-7608-02     | 2" (5.1cm)   |
| 6S-7608-S2     | 2.5" (6.4cm) |
| 6S-7608-03     | 3" (7.6cm)   |
| 6S-7608-04     | 4" (10.2cm)  |
| 6S-7608-05     | 5" (12.7cm)  |
| 6S-7608-06     | 6" (15.2cm)  |
| 6S-7608-08     | 8" (20.3cm)  |
| 6S-7608-10     | 10" (25.4cm) |
| 6S-7608-12     | 12" (30.5cm) |

- ✓ 100% cotton
- ✓ Sold boxed 25 yard (23m) roll



VALUE LINE RIB-KNIT STOCKINETTE

## VALUE LINE RIB-KNIT STOCKINETTE

| NATURAL COTTON | WIDTH         |
|----------------|---------------|
| 6S-VL08-02     | 2" (5.1cm)    |
| 6S-VL08-S2     | 2.50" (6.4cm) |
| 6S-VL08-03     | 3" (7.6cm)    |
| 6S-VL08-04     | 4" (10.2cm)   |
| 6S-VL08-05     | 5" (12.7cm)   |
| 6S-VL08-06     | 6" (15.2cm)   |
| 6S-VL08-08     | 8" (20.3cm)   |
| 6S-VL08-10     | 10" (25.4cm)  |
| 6S-VL08-12     | 12" (30.5cm)  |

- ✓ Cost-effective
- ✓ Thinner and stretchier than standard stockinette
- ✓ Sold boxed 25 yard (23m) roll




COTTON RIB-KNIT A/K CASTING PANTS

## A/K CASTING PANTS

| COTTON RIB-KNIT      | WIDTH |
|----------------------|-------|
| COTTON RIB-KNIT      |       |
| 6S-CPCT-03           | 6"    |
| 6S-CPCT-04           | 8"    |
| 6S-CPCT-05           | 10"   |
| 6S-CPCT-06           | 12"   |
| COTTON/COMPRESSOGRIP |       |
| 6S-CPCC-03           | 6"    |
| 6S-CPCC-04           | 8"    |
| 6S-CPCC-05           | 10"   |
| 6S-CPCC-06           | 12"   |

For modesty, have the patient don two cast pants, leave one on patient when the cast is removed.

- ✓ Made from tubular stockinette, split and seamed to form a body section with two open end legs
- ✓ Permits casting to the groin
- ✓ The distal end for the amputated side may be cut and sewn to proper length (optional)
- ✓ Leg length is long enough for knee disarticulation amputations
- ✓ Over 150% stretch
- ✓ Sold each

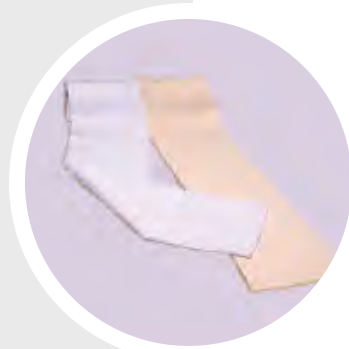

COTTON/COMPRESSOGRIP A/K CASTING PANTS

**Helpful Hint:** First apply a "pull sock" to the limb; then apply the A/K Casting Pants; pull out the "pull sock" to pull down the soft tissue

# STOCKINETTES AND CASTING

## ACCUCAST® STRETCH STOCKINETTE

| ITEM # | WIDTH X LENGTH* | STRETCH RANGE |
|--------|-----------------|---------------|
|--------|-----------------|---------------|

### ACCUCAST® STOCKINETTE ROLLS

|            |                                |                     |
|------------|--------------------------------|---------------------|
| 6S-CLT2-25 | 2.75" (7cm) x 25 yards (23m)   | 5"-22" (12"-56cm)   |
| 6S-CLR6-25 | 6.25" (16cm) x 25 yards (23m)  | 11"-46" (28"-116cm) |
| 6S-CLS8-25 | 8.50" (22cm) x 25 yards (23m)  | 15"-50" (37"-125cm) |
| 6S-CL11-25 | 11.50" (29cm) x 25 yards (23m) | 20"-60" (50"-150cm) |

\*Applying horizontal stretch will shorten the length and vertical stretch will reduce the width circumference.

- ✓ Smooth and 4-way stretchy cotton/LYCRA® casting stockinette
- ✓ Available in rolls and casting socks, it is knit specifically to enable superior limb impressions to be made
- ✓ Smoother surface than traditional "rib-knit" stockinette, providing a model that is free of the ridges
- ✓ Smoother models mean less modification time and greater model accuracy
- ✓ The LYCRA® in AccuCast® allows the casting interface to fit the limb or torso with greater cling and conformity for more accurate impressions
- ✓ Cost is less than one dollar for an average AFO impression



ACCUCAST® STOCKINETTE ROLLS



ACCUCAST® CAST SOCKS

## ACCUCAST® CASTING SOCKS

| ITEM # | WIDTH X LENGTH* |
|--------|-----------------|
|--------|-----------------|

|            |                            |
|------------|----------------------------|
| 6S-CST2-24 | 2.75" (7cm) x 24" (60 cm)  |
| 6S-CSR6-24 | 6.25" (16cm) x 24" (60 cm) |


- ✓ Cotton/LYCRA®
- ✓ Seamed toe

## CAST SOCKS

| ITEM # | STYLE | WIDTH X LENGTH | STRETCH RANGE |
|--------|-------|----------------|---------------|
|--------|-------|----------------|---------------|

|            |                 |             |        |
|------------|-----------------|-------------|--------|
| 6S-7924-LT | B/K Lightweight | 3.50" x 24" | 7-24"  |
| 6S-7936-LT | B/K Lightweight | 3.50" x 36" | 7-24"  |
| 6S-7924-HV | B/K Heavyweight | 3.50" x 24" | 7-20"  |
| 6S-7918-LT | A/K Lightweight | 5.50" x 22" | 11-28" |
| 6S-7918-HV | A/K Heavyweight | 6" x 22"    | 12-30" |

- ✓ Flat, tubular knit, soft, carded cotton yarn
- ✓ Seamed toe
- ✓ Sold per dozen

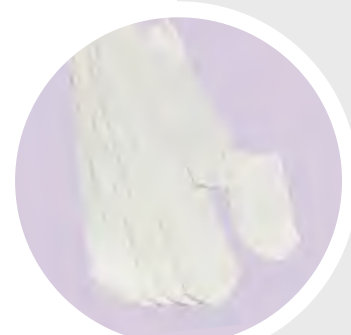

CAST SOCKS



# TERMS AND CONDITIONS

## PAYMENT/REMITTANCE/RETURNS

120 Osage Avenue, Kansas City, KS 66105

### TOLL FREE PHONE

### TOLL FREE FAX

### KNITRITE.COM

## TERMS

### NET 30 FROM DATE OF INVOICE

Forms of payment include check, money order/cashier's check, wire transfer or upon approval, credit card (Visa, MasterCard, and Discover). Orders will be shipped C.O.D or CREDIT CARD until credit terms are approved upon submission of a credit application and tax exemption certificate. Upon credit approval, open terms will be established. Credit approval remains subject to continued review even after credit is approved. We will not DROP SHIP C.O.D. ORDERS. If C.O.D. shipments are refused, all costs incurred plus a charge will be billed. Statements of Account are sent the middle of each month.

Our billing month covers invoices, credit memos and payments for such month. Service charges are based on the amount past due and subject to an interest charge of 1.5% per month, effective annual rate of 18% (or the highest rate permitted by applicable law).

Collection Fee: All collection charges and fees incurred in collecting past due accounts will be billed and paid by customer and must be paid in full before shipment of orders are resumed.

Returned Checks: There will be a service charge on all checks returned unpaid.

## SHIPPING

### FLAT RATE SHIPPING POLICY:

\*\*\*Flat Rate Shipping and Free Freight programs limited to qualified accounts. All shipments will be made using

UPS Ground, UPS Surepost or similar service. All expedited services, including Air services, excluded. Call your Customer Service Specialist for more information.

Drop shipments placed on wholesale.therafirm.com will have no drop ship fee and only a Flat Rate of for shipping. All other drop shipments placed outside the online wholesale portal will have a drop ship fee plus the Flat Rate for shipping.

Additional shipping rates and shipping services will apply to Hawaii, Alaska and Puerto Rico.

## DELIVERY CLAIMS

No claim for shortages will be allowed unless made within 10 working days of receipt of order.

## SALES ONLY IN THE US

Customers shall not directly or indirectly (including, without limitation, via the Internet) market, sell or export the products to customers located outside of the United States or sell, service or promote the products from a location outside of the United States. If customer is a distributor, customers of the distributor may only sell the products to end users in the United States.

## ONLINE RESELLERS

Knit-Rite, LLC /Therafirm customers may not re-sell or distribute to anyone that may re-sell any Knit-Rite/Therafirm products on any online marketplaces, including but not limited to Amazon.com, Wal-Mart.com, eBay.com, without written consent from Knit-Rite. Current online retailers are not automatically granted permission to sell on online marketplaces. All customers must have written consent from Knit-Rite to sell or distribute to customers that sell through these specific channels. Knit-Rite reserves the right to not sell to any customer who violates this policy.

## RETURNED GOODS POLICY

Acceptance of returned merchandise is at the discretion of Knit-Rite/Therafirm. Returned merchandise will not be accepted for credit without a valid return authorization number and the lot number for the product. There is a maximum 45-day time limit in which to return all merchandise, should Therafirm approve the return. Note: The RA# authorizes your company to return the items, however, based on evaluation, credit to be determined.

The following will be required for any return:

\*Full description of the merchandise.

\*Complete reason for the return.

The following applies to all returned merchandise that is not deemed defective or that is not the result of a shipping error:

Please DO NOT return any merchandise that has been worn or tried on. Credit or exchange will be issued at the discretion of Knit-Rite/Therafirm.

\*0% Restock fee if returned within 45 days and no repackaging required.

\*15% Restock fee if repackaging is required.

\*15% Restock fee if product is returned without an RA.

\*15% Restock fee if returned over 45-days with approval.

All returns must be sent to:

120 Osage Avenue, Kansas City, KS 66105

## EXPORT ACCOUNTS

See your Sales Representative for specific Terms and Conditions.

### KNIT-RITE/THERAFIRM MANUFACTURER'S WARRANTY:

All Knit-Rite/Therafirm products are warranted to be free of material or workmanship defects for 45 days from date of purchase. At the discretion of Knit-Rite/Therafirm, any defective merchandise will be replaced or credited at the actual price paid. No refunds.

THIS WARRANTY IS MADE EXPRESSLY IN LIEU OF ANY OTHER WARRANTY AND KNIT-RITE MAKES NO OTHER REPRESENTATION OR WARRANTY OF ANY OTHER KIND, EXPRESS OR IMPLIED, WITH RESPECT TO THE PRODUCTS, WHETHER AS TO MERCHANTABILITY, FITNESS FOR A PARTICULAR PURPOSE INCLUDING ANY PARTICULAR MEDICAL CONDITION, NON-INFRINGEMENT OF THIRD PARTY RIGHTS OR ANY OTHER MATTER. THE WARRANTIES CONTAINED HEREIN SET FORTH CUSTOMER'S SOLE AND EXCLUSIVE REMEDY IN THE EVENT OF A DEFECT IN WORKMANSHIP OR MATERIALS AND SHALL NOT BE DEEMED TO HAVE FAILED ITS ESSENTIAL PURPOSE SO LONG AS KNIT-RITE IS WILLING AND ABLE TO CARRY OUT THE TERMS OF THIS EXCLUSIVE WARRANTY.

Customers should carefully inspect their products as soon as they receive them for manufacturer's defects. Manufacturer's defects are flaws that are existent or apparent when the package is opened.. Note that manufacturer warranties do NOT include: normal wear and tear, pulled threads, runs, wear at the toe, or holes due to fingernails, toenails, or snags; or any problem caused by the wearer in application, improper use, or improper care of the garment. It is recommended that compression products be replaced every 6 months as the compression can break down with normal wear and care. These occurrences will not be considered as defects. If a customer feels that a hole, snag, tear or run is the result of a manufacturer's defect, please contact us. We take quality issues very seriously and we will do our best to resolve the matter.

Knit-Rite/Therafirm products purchased from unauthorized resellers are not covered by our warranty. Knit-Rite/Therafirm products purchased from any marketplace seller will not be covered by our warranty.

## COPYRIGHT NOTICE

AccuCast®, A-Plus™, Cast-Rite™, Compressogrip®, Knit-Rite®, Liner-Liner®, Micro-Cool®, PowerBelt® and PowerPull®, Sleeve Saver™, SmartKnit®, SmartKnitKIDS®, SmartKnitBIGKIDS™, Soft-Sock®, Sterling™, Super Sock®, Super Stretch Sheaths™, Torso-Interface®, Ultra Low Friction Sheath™, and X-Wool® are trademarks and/or registered trademarks of Knit-Rite, LLC.



Ionic™ is a trademark of NOBLE Biomaterials, Inc.  
LYCRA® and COOLMAX® are registered trademarks of The LYCRA Company.



**KNIT-RITE'S MISSION** is to provide medical devices that lift individuals to new levels, and that provide solutions for recovery, self care, health care, and sports care.

**KNIT-RITE'S VISION** is to continue our tradition of producing the highest quality O&P, compression and sensory products that exceed the needs of the patients who use them, distributors who supply them, and physicians who recommend them, and that our employees take pride in producing.

**OUR CORE VALUES** are based on our history of:

- ✓ **Teamwork**  
No one goes it alone
- ✓ **Service to others**  
Our customers and consumers are at the center of our philosophy, operations, and ideas
- ✓ **Service to our communities**  
Giving back to our communities is an integral part of what we are about
- ✓ **Collaboration**  
Together, we learn, improve, and achieve more success for ourselves and our customers and patients
- ✓ **Integrity**  
We seek to always care for and do the right thing for those we serve
- ✓ **Quality**  
Providing competitively priced solutions that work, made by people who care



**Catherine Cassel**  
GENERAL MANAGER

Cell: 083 274 7820  
Office: 012 - 991 2227  
Email: catherine@medroma.co.za  
Orders: orders@medroma.co.za

**MEDROMA**

[www.medroma.co.za](http://www.medroma.co.za)

ORTHOTICS

PROSTHETICS

COMPRESSION GARMENTS

**Marike du Preez**  
HEAD OF PROSTHETICS

Cell: 083 381 0454  
Office: 012 - 991 2227  
Email: marike@medroma.co.za  
Orders: orders@medroma.co.za

**MEDROMA**

[www.medroma.co.za](http://www.medroma.co.za)

ORTHOTICS

PROSTHETICS

COMPRESSION GARMENTS

**KNIT-RITE®**  
THUASNE